

# Behavioral Observation Research Interactive Software (BORIS) user guide



version 2.993

BORIS web site: [www.boris.unito.it](http://www.boris.unito.it)

[BORIS GitHub repository](#)

[BORIS twitter feed](#)

[BORIS discussion group](#)

[BORIS Facebook page](#)

[BORIS Google+](#)

# Contents

<b>Legal</b>	<b>3</b>
<b>Install BORIS</b>	<b>3</b>
<b>Launch BORIS</b>	<b>3</b>
<b>Create a new project with BORIS</b>	<b>4</b>
Set an ethogram	4
Define the subjects	9
Define the Independent variables	10
Observations' tab	10
<b>Open an existing project with BORIS</b>	<b>10</b>
<b>Observations list</b>	<b>10</b>
<b>Delete an observation</b>	<b>11</b>
<b>Create a new observation</b>	<b>11</b>
Live observation	12
Media based observation	13
<b>Coding your media</b>	<b>15</b>
Media controls are available in the toolbar	15
<b>Ethogram</b> widget in the main window	16
<b>Subjects</b> widget in the main window	16
<b>Media player</b> widget in the main window	17
<b>Events</b> widget in the main window	17
Record an event	18
Frame-by-frame mode	19
<b>Exporting events data</b>	<b>19</b>
Export events	20
Export aggregated events	20
Export events as behavioral strings	21
Export events as Praat TextGrid	21
<b>Tools</b>	<b>22</b>
Geometric measurements	22
Extract media file sequences	23
<b>Coding maps</b>	<b>24</b>
Creating a coding map	24
<b>Analysis</b>	<b>25</b>
Time budget	25
Plot events	26
<b>Preferences</b>	<b>28</b>

General preferences	28
FFmpeg framework	29
<b>Various</b>	<b>29</b>
Citing BORIS	29
Bug reports and features request	29
Docking	30
<b>Acknowledgement</b>	<b>30</b>
Valid keys for triggering behavior	30

## Legal

Copyright 2012-2016 Olivier Friard - Marco Gamba

**BORIS** is free software; you can redistribute it and/or modify it under the terms of the GNU General Public License as published by the Free Software Foundation; either version 2, or any later version.

**BORIS** is distributed in the hope that it will be useful, but WITHOUT ANY WARRANTY; without even the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the [GNU General Public License](#) for more details.

## Install BORIS

BORIS can be installed following instructions on the BORIS site at [download page](#)

Unless you use a full bundle version (available only for Microsoft-Windows) BORIS requires a **recent** version (>=2) of the [VLC Media Player](#) from the VideoLAN Organization.

*VLC is a free and open source cross-platform multimedia player and framework that plays most multimedia files as well as DVD, Audio CD, VCD, and various streaming protocols.*

From version 2.8 the open source [FFmpeg multimedia framework](#) is included with BORIS.

If you want to try BORIS without installing it you can download the [VirtualBox virtual appliance](#) from the BORIS site. In this case you must first install [VirtualBox](#)

## Launch BORIS

Start BORIS by clicking on its icon. For launching BORIS from sources see the BORIS web site. When launching BORIS for the first time it may take some time to show up. Please be patient!

### Warning

On Mac OS X 10.9, launching can be stopped according to the **Security & Privacy** settings of your computer. They can be changed opening **System Preferences...** > **Security & Privacy** > **General** and selecting the option **Anywhere** in the frame **Allow apps downloaded from**. Alternatively you can repeat the operation right-clicking on the BORIS icon and then clicking **Open** in the following two dialog boxes.

## **Warning**

On Microsoft-Windows, launching can be stopped by a **Security warning**: *The publisher can not be verified. Are you sure you want to run this software?* Click the **Run** button to launch BORIS.

At the first launch, BORIS will ask you to allow the automatically check for new version. Internet access is required for this feature. If you choose **Yes** BORIS will check for new version about every 15 days. This option can be changed on the **Preferences** window (See [general preferences](#))

## **Create a new project with BORIS**

The BORIS project file is the container for all information related to the project. It contains the ethogram, the independent variables, the subjects and all observations data. The save menu option (File -> Save project) will save all this information. You can also activate the **automatic backup** feature (see **Preferences**). BORIS allows creating an unlimited number of projects but only one project can be opened at a time.

To create a new project, under the menu **File** , select **New project**. You can determine your project name by writing in the **Project name** field in the **Information** tab. Once the project will be saved, the **Project file path** will show the full path to your project file. **Date** will automatically set on the current date and time, but you can alternatively set this info on your media date and time, or whatever you prefer. **Description** can host all the relevant information about your project, can be also left empty. **Time format** can be alternatively set to **seconds** or to **hh:mm:ss.mss**. This choice can be changed at anytime under **File > Preferences** (for MAC users, **BORIS > Preferences**) > **Default project time format**.



## **Set an ethogram**

Switching to the **Ethogram** tab, you can alternatively:

- set your ethogram from scratch;
- import an existing ethogram from another BORIS project;
- import an ethogram from a JWatcher global definition file (.gdf).
- [import an ethogram from a plain text file](#)

edit project

Behavior type	Key	Code	Description	Modifiers	Exclusion	Coding map
1 Point event	J	jump	jumping			
2 State event	G	groom	Animal engages in washing or smoothing	self,others	eat,look for food,sleep,drink,fight,play,locomote	
3 State event	E	eat	Animal consumes food		groom,look for food,sleep,drink,fight,play	
4 State event	O	look for food	Animal searches the environment for food		groom,eat,sleep,drink,fight,play	
5 State event	S	sleep	Animal assumes position for sleep, stays in one place and is not alert to environmental changes		groom,eat,look for food,drink,fight,play,locomote	
6 State event	D	drink	Animal consumes water or others liquids found in its environment		groom,eat,look for food,sleep,fight,play,locomote	
7 State event	F	fight	Animal engages a physical contact with another animal		groom,eat,look for food,sleep,drink,play,locomote	
8 State event	P	play	Animal engages in interactions with others		groom,eat,look for food,sleep,drink,fight,locomote	
9 State event	L	locomote	Animal moves from place to place		groom,sleep,drink,fight,play	

## **Set your ethogram from scratch**

Clicking on the **Add behavior** button you can add a new row in the **Ethogram** table and behavior type will be automatically set to **Point event**. You can switch between **State event** and **Point event** at your convenience from the **Behavior type** column. You can also add a **Coding map** to either a **State event (State event with coding map)** or a **Point event (Point event with coding map)**; see the "Coding map" section for details).

An existing behavior can be duplicated using the **Clone behavior** button. Its code have then to be changed. On a selected behavior, click on the **Remove behavior** button to remove. The **Remove all behaviors** button will clear the **Ethogram** table. Both the above-mentioned operations must be confirmed when prompted.

Behavioral codes (**Code** column) can be sorted alphabetically by checking the **Alphabetical order** checkbox. Alternatively they can be sorted manually by using the **Move up** and **Move down** buttons.

## **Categories of behaviors**

Defining categories of behaviors can be usefull for the analysis of coded events (for example the time budget analysis). Click the **Behavioral categories** button and add a the categories of behaviors. Behaviors can then be included or not in a defined category.



### **Set keys and codes**

For each behavior you have to set a keyboard key (**Key** column) that will be then used to code the behavioral events. You can choose whether you want to set a unique key for each behavior or use the same key for more than one behavior. In the case you set the same key for more than a behavior, BORIS will pause your coding and ask which of the behavior you want to record. The keys are case-insensitive.

#### **Important**

Do not use the / and \* keys! They are reserved for the frame-by-frame mode.

In the **Code** column, you have to add a unique code for each behavior. Duplicated codes are not accepted and BORIS will warn in red about duplicates on the bottom left of the *Ethogram* tab. The code can be an alphanumeric string (which must not include the pipe character | ).

The **Category** column allow you to include the behavior to a predefined category.

The **Description** of your behavior is optional. The **Description** column can be useful to add information about a specific behavior, its characteristics (e.g. to standardise observation between different users) or to refer to external information (e.g. reference to a previous ethogram).

The following three columns (**Modifiers**, **Exclusion**, **Coding map**) cannot be edited from the *Ethogram* table.

### **Set the modifiers**

Modifiers can be used to add attributes to a behavior. A single behavior can have two or more modifiers attached (e.g. "play" may have "solitary" or "social" as modifiers). The use of modifiers can be convenient to significantly reduce the number of keys and simplify the behavioral coding. In BORIS modifiers can also be added in different modifier sets [e.g. "play, social" may have a modifier set (#1) for "brothers" and another (#2) for "sisters"]. In the case of using sets of modifiers, you can attach one modifier for each set.

To add modifiers to a behavior, you need to double-click the **Modifiers** cell corresponding to the behavior you want to add the modifiers to.



When the **Modifiers** window is open you can add and/or remove sets using the buttons **Add set of modifiers** and **Remove set of modifiers**. Within a set of modifiers, you can add and remove modifiers using the **Modifier** field and clicking on the **right-arrow** button (to add) and the **Remove modifier** button (to remove). The selected modifier can be edited using the left-arrow button. The **Key code** box can be used to set a key for the modifier (optional).

The modifier position into the modifiers' set can be manually set using the **Move modifier up** and **Move modifier down** buttons.

The position of a modifiers' set can be customized (using the **Move set left** and **Move set right** buttons)

Modifiers can not contain the following characters: (, ) and | (pipe character).

Click **OK** to save modifiers in the **Ethogram** table.

#### **Set the exclusion matrix**

The occurrence of an event (State or Point) can exclude the occurrence of a state event. This can be set using the **Behaviors exclusion matrix** window, which can be opened clicking on the **Exclusion matrix** button. BORIS will ask for including **Point events** or not and a new **Exclusion matrix** window will open.

Exclusive behavior may be selected by checking on the corresponding cell in the automatically-generated matrix. We suggest to work on the **Exclusion matrix** when all the behaviors have been added to your ethogram.

For example, for the following ethogram:



The **Exclusion matrix** could be:

Behaviors exclusion matrix		
	c	d
a	<input checked="" type="checkbox"/>	<input type="checkbox"/>
b	<input type="checkbox"/>	<input checked="" type="checkbox"/>
c		<input type="checkbox"/>
d	<input checked="" type="checkbox"/>	

- the behavior **a** (Point event) will exclude the behavior **c** (State event) but not **d** (State event)
- the behavior **b** (Point event) will exclude the behavior **d** but not **c**
- the behavior **d** (State event) will exclude the behavior **c**

During the observation, the excluding event will stop all the current excluded state events one millisecond before the occurrence of the event.

### ***Import an ethogram from an existing project***

Behaviors within an ethogram can be imported from an existing BORIS project (.boris) using the **Import behaviours from a BORIS project** button. BORIS will ask to select a BORIS project file and whether imported behaviors should replace or be appended to the **Ethogram** table. Imported behaviors will retain all the previously defined behavior parameters (namely Behavior type, Key, Code, Description, Modifiers and Exclusion information).

## **Import an ethogram from a JWatcher global definition file (.gdf)**

Behaviors can be imported from a JWatcher global definition file (.gdf) using the **Import from JWatcher** button. BORIS will ask to select a JWWatcher file (.gdf) and whether imported behaviors should replace or be appended to the **Ethogram** table. Behavior type and exclusion information for the behaviours imported from JWWatcher have to be redefined.

## **Import an ethogram from a plain text file**

Behaviors can be imported from a plain text file using the **Import from text file** button. The fields must be separated by TAB, comma (,) or semicolon (;). All rows must contain the same number of fields.

The fields will be interpreted as:

- field #1: event type (point or state)
- field #2: key (case insensitive)
- field #3: code (must be unique)
- field #4: behavior category (empty if no category)
- field #5: description (optional)

All fields after the 5th will be ignored.

BORIS will ask to select a plain text file (by default: \*.txt \*.csv \*.tsv) and whether imported behaviors should replace or be appended to the **Ethogram** table. The missing information for the behaviours imported from text file have to be redefined.

## **Define the subjects**



BORIS allows coding behaviors for different subjects within a single observation. The **Subject** table allows specifying subjects using a **Key** (e.g. the "K" on your keyboard), **Subject name** (e.g. "Kanzi"), **Description** (e.g. male, born October 28 - 1980). In this case, pressing "K" will set "Kanzi" as the focal subject of the behavioural coding. Pressing "K" again will deselect "Kanzi" and set to "no focal subject". The definition of one or more subjects is not mandatory. Addition, removal and sorting of the subjects follows the same logic of the **Ethogram** table (see **Set your ethogram from scratch** for info). The subjects can also be imported from an existing BORIS project using the **Import Subjects from a BORIS project**.

## Define the Independent variables



BORIS allows adding information about the observation using **Independent variables**. This can be used to specify factors that may influence the behaviors (e.g. group composition, temperature, weather conditions) but will not change during a single observation within a project. Each *Independent variable* can be defined by **Label** (e.g. weather), **Description** (e.g. weather conditions), **Type** (e.g. **text** for weather condition "sunny", **numeric** for temperatures). A value for each **Independent variable** will be requested when creating a new observation. Addition, removal and sorting of the independent variables follows the same logic of the **Ethogram** table (see **Set your ethogram from scratch** for info). The independent variables can also be imported from an existing BORIS project using the **Import Variables from a BORIS project**.

## Observations' tab

The **Observations** table in BORIS shows information about all the previous observations within a project. A selected "Observation" can be removed using the "Remove observation" button (you will be prompted for confirmation). This operation cannot be undone and deleted observations cannot be recovered once the project is saved. The **Observations** table shows four columns **id Date Description Media**.

## Open an existing project with BORIS

To open an existing BORIS project, under the menu **File**, select **Open project**. A BORIS project file is a container for all information related to a set of observations as the ethogram, the independent variables, and the subjects. BORIS allows creating an unlimited number of projects but only one project can be opened at a time.

## Observations list

The **Observations > Observations list** will show you all the observations contained in the current BORIS project.

Observations list - BORIS								
2 observation(s)								
<input type="button" value="id"/> <input type="button" value="date"/> <input type="button" value="description"/> <input type="button" value="subjects"/> <input type="button" value="media"/> <input type="button" value="meteo"/> <input type="button" value="Group size"/> <input type="button" value="Number of visitors"/>								
1	observation #1	2013-09-06 12:44:42	First observation at Zoom Torino	Bobo	#1:/home/user/Videos/lemurs.mp4	sun	4	125
2	observation #2	2016-05-09 13:50:32	Second observation at Zoom Torino		#1:/home/user/Videos/crop.avi	rain	6	250

The observations list can be filtered selecting a field in the drop-list box.

In the following example observations are filtered: only observations containing the **Bobo** subject are shown:

Observations list - BORIS								
1 observation(s)								
<input type="button" value="subjects"/> <input type="button" value="bobo"/>								
1	observation #1	2013-09-06 12:44:42	First observation at Zoom Torino	Bobo	#1:/home/user/Videos/lemurs.mp4	sun	4	125

## Delete an observation


An observation can be deleted from the project using the following procedure: **File > Edit project > Observation's tab > Select observation to remove > Remove observation** button

Please note that the deletion is irreversible.

## Create a new observation

To create a new observation you must first [Create a new project with BORIS](#) or [Open an existing project with BORIS](#).

Clicking on **Observations > New observation** will show the **New observation** window.

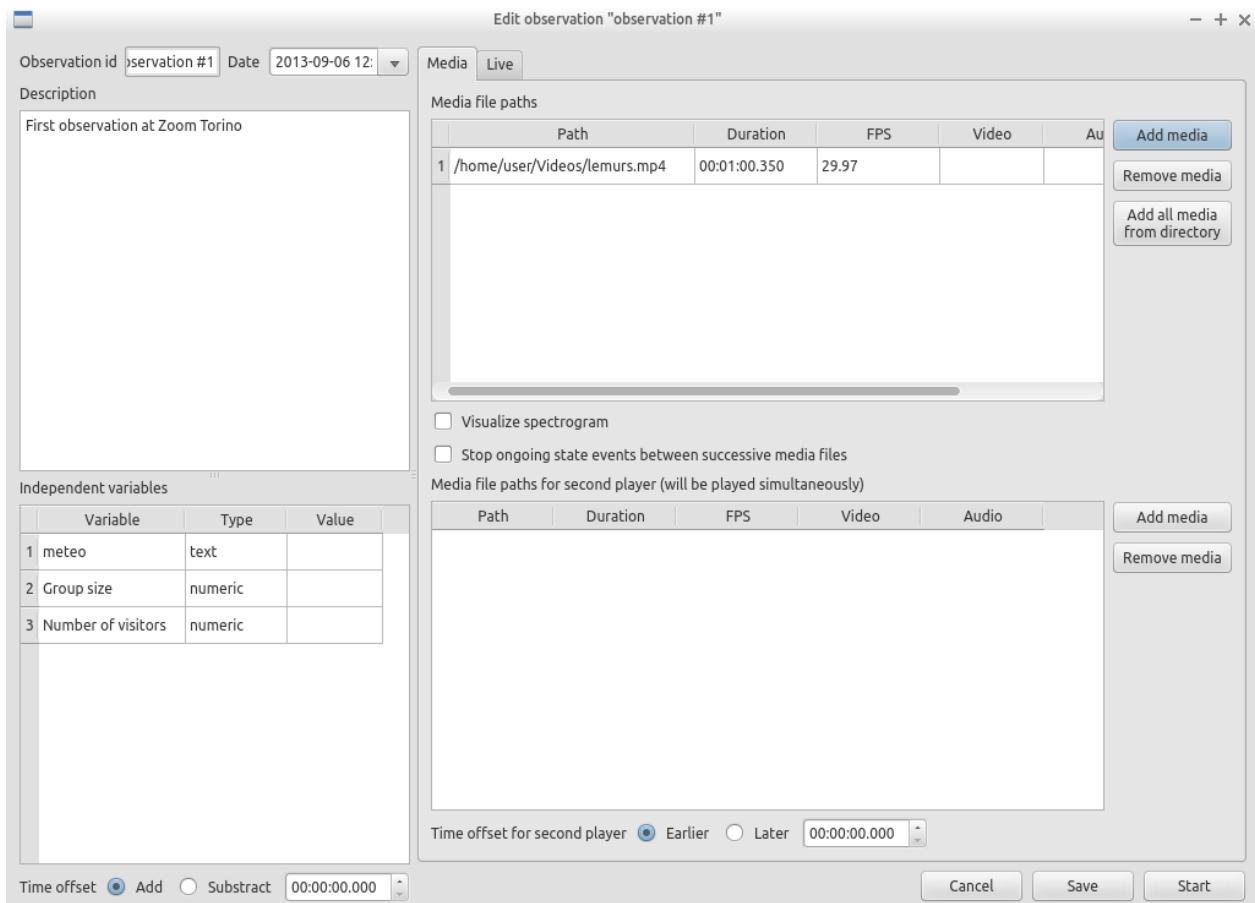

This window allow adding various observation data:

- a mandatory **Observation id** (must be unique across all observations);
- **Date**, which will be automatically set on the current date and time, but you can alternatively set this info on your media date and time, or whatever you prefer.
- **Description**, which can host all the relevant information about your observation, but can be also left empty.
- **Independent variables** (e.g. to specify factors that may influence the behaviors but will not change during the observation within a project). See the [independent variables](#) section for details.
- **Time offset**. BORIS allow specifying a time offset that can be added or subtracted from the media timecode.

You must then indicate if you want to make an observation based on pre-recorded media (audio / video) or a live observation.

## Live observation

During the live observation BORIS will show you a timer that will be used for recording time for coded events.

Choose the **Live tab** to make a live observation.



In the above tab you can select a time for **Scan sampling** observation. In this case the timer will stop at every time offset you indicated and all the coded events will have the same time value.

## Media based observation

Choose the **Media** tab to make an obervation based on media.



In the **Media tab** there are two playlists. In the **Media file paths** playlist you can add one or more media files using the **Add media** button. Information about the selected media file will be extracted and displayed in the media list: media file path, media duration (in seconds), number of frames by second (FPS), video stream, audio stream .

If you have to add various media files, you can use the **Add all media from directory** button, in this case all the media files found in the selected directory will be added in the playlist.

The **Remove media** button can be used to remove the selected media files.

All the media types reported at <http://www.videolan.org/vlc/features.html> can be played in BORIS. The media queued in the *Media file paths* will be played sequentially. This means that an event occurring at time  $t_x$  in the media file queued as second (e.g. second\_video.mp4) in the playlist will be scored as happening at time  $t_1 + t_x$  (where  $t_1$  is the duration of the first media file, e.g. first\_video.mp4).

## Spectrogram visualization

BORIS allow you to visualize the sound spectrogram during the media observation. Activate the *Visualize spectrogram* check box. BORIS will ask you to generate the spectrograms for all media files loaded in the first player.



The spectrogram visualization will be synchronized to the media position during the observation.



## Close current behavior between videos

If this option is selected BORIS will close all ongoing State events between successive media files.

This option can be useful if you have to code various short media files enqueued in the first player (for example the output of a Camera trap).

## Simultaneous play

BORIS also allows simultaneous playback of two media recorded independently (e.g. videos recorded from different points in a room; or a video and an audio recording of the same observation). The videos to be played simultaneously can be loaded in the **Media file paths for second player** playlist. In this case only one video per playlist is accepted. If the two media are not synchronised you can specify a time offset for the second player.

Click OK to start coding. The **Observation** window will be closed and you'll be transferred to the main **BORIS** window.

# Coding your media

When looking at the BORIS main window, the window title bar shows the **Observation id - Project name - BORIS**. The media (the first in the queue) will be loaded in the media player and paused.

## Media controls are available in the toolbar



Key to the symbols:

- **Play**
- **Pause** (the SPACE bar can be used)
- **Rewind** reset your media at the beginning
- **Fast backward** jumps for n seconds backward in your media (See [general preferences](#) to set n)
- **Fast forward** jumps for n seconds forward in your media (See [general preferences](#) to set n)
- **Set the playback at speed 1x**
- **Increase playback speed** (See [general preferences](#) to set the step value)
- **Decrease playback speed** (See [general preferences](#) to set the step value)
- **Jump to the previous media file**
- **Jump to the next media file**
- **Take a snapshot** of current video or frame
- **Switch between VLC and frame-by-frame modes**

In frame-by-frame mode:

- **Move on frame back**
- **Move one frame forward**

The media can be controlled by special keyboard keys:

- **Page Up** key: switch to the next media
- **Page Down** key: switch to the previous media
- **Up arrow** key: jump forward in the current media
- **Down arrow** key: jump backward in the current media
- **ESC**: switch between VLC and frame-by-frame mode

## Ethogram widget in the main window



The **Ethogram** widget provide the user with the list of behaviors defined in the **Ethogram**. It can be used to record an event by double clicking on the corresponding row. The **Key** column indicates the keyboard key assigned to each behavior (if any). Pressing a key will record the corresponding behavior (that will appear in the *Events* widget).

## Subjects widget in the main window

The screenshot shows the BORIS software interface with the title bar "Test 1 - BORIS". Below the menu bar (File, Observations, Playback, Tools, Analyze, Help) is a toolbar with various playback controls (rewind, fast forward, stop, play). The main area is titled "Subjects" and contains a table with the following data:

	Key	Name	Description	Current state(s)
1	1	Kanzi	Male, born October 28 - 1980	
2	2	Bongo	Male, adult	
3	3	Marco	Orangutan	

The **Subjects** widget provide the user with the list of subjects defined in the **Subject** tab in the **Project** window. It can be used to add information about the focal subject on the recorded behaviors by double clicking on the corresponding row. When a subject is selected his/her name appears above the media player. The **Key** column indicates the keyboard key assigned to each subject (if any).

## Media player widget in the main window



The **Media player** widget has two controls: the media position (horizontal slide bar) and the audio volume (vertical slide bar) provide the user with the list of subjects defined in the **Subject** tab in the **Project** window. The horizontal slide bar can be used to navigate the media file.

## Events widget in the main window

Events for observation #4					
	time	subject	code	type	modifier
1	00:00:16.040	Kanzi	eat	START	alone orange
2	00:00:30.320	Kanzi	eat	STOP	alone orange
3	00:00:44.240	Bongo	eat	START	in group banana
4	00:00:53.740	Bongo	eat	STOP	in group banana
5	00:01:11.480	Bongo	face		orb occ
6	00:01:21.920	Bongo	groom	START	self
7	00:01:31.240	Bongo	groom	STOP	self

0:01:36.199 / 00:10:22.760 (paused) Subject: Bongo

The **Events** widget shows all the recorded behaviors with the following parameters (columns):

- **time**, the time at which the event occurred;
- **subject**, the focal subject (if any);
- **code**, the behavior code;
- **type**, in case of a state event indicates whether the time corresponds to the start or to the stop.
- **modifier**, indicates the modifier(s) that was(ere) selected (if any);
- **comment**, is an open field where the user can add notes.

A tracking cursor (red triangle) will show the current event. This cursor can be positioned above the current event, see [tracking cursor position](#) option in Preferences window.

A double-click on a row will reposition the media player to the moment of the corresponding event. See [Time offset for media reposition](#) in Preferences window to customize the time offset for media repositioning.

## Record an event

Once ready to begin your coding, you can start the media player using the **Play** button (or the Space bar). The behaviors can be recorded using the keyboard with the predefined keys (or by double-clicking the corresponding row in the **Ethogram** table).


The screenshot shows the BORIS software interface. At the top is a menu bar with File, Observations, Playback, Tools, Analyze, and Help. Below the menu is a toolbar with playback controls (rewind, play, fast forward, etc.). The main area contains two tables: 'Ethogram' and 'Subjects'.  
**Ethogram Table:**

Key	Code	Type	Description	Modifiers	Excluded
1 T	moves towards	State event			moves backwards
2 B	moves backwards	State event			moves towards
3 J	jumps	Point event	jumps		
4 H	hits	Point event	hits ball or oppo...	ball, opponent	
5 A	goal	Point event	scores goal		
6 A	suffers goal	Point event	suffers goal		

  
**Subjects Table:**

Key	Name	Description	Current state(s)
1 K	korea		
2 G	ghana		

If the pressed key defines a single event, the corresponding event will be recorded in the **Events** table. In the case you have specified the same key for two (or more) events (e.g. key A in the figure below), BORIS will prompt you for the desired behavior.

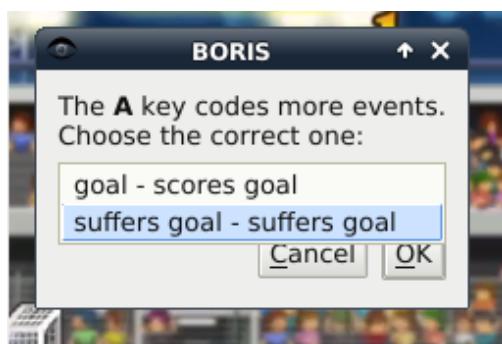

In the case you have specified modifiers (one or more sets), BORIS will prompt you for the desired modifier(s) if any (e.g. **ball** or **opponent** in the figure below). You can select the correct one using the mouse or the keyboard (**b** key or **o** key)



In the case your behavior type is a *Point event with coding map* or a *State event with coding map*, BORIS will show the *Coding map* window and will allow selecting the desired area(s). In case you click a part of the map in which two (or more) areas overlap, the corresponding codes will be recorded.

A recorded event can be edited (once selected) using the *Observations > Edit event* menu option. The resulting *Edit event parameters* allows modifying every parameter (e.g. time, subject, code, modifiers, and comment).

The *Observations > Add event* menu option allows adding a new event by specifying its time and the other parameters.

## Frame-by-frame mode

You can switch between play and frame-by-frame modes using the dedicated button in the toolbar:



In frame-by-frame mode the video will stop playing and the user will visualize the video frame by frame.

You can move between frames by using the arrow keys in the toolbar (on the right) or by using keyboard special keys:

**For the both modes** (VLC and frame-by-frame):

- **Page Up** key: switch to the next media
- **Page Down** key: switch to the previous media
- **Up arrow** key: jump forward in the current media
- **Down arrow** key: jump backward in the current media
- **ESC**: switch between VLC and frame-by-frame mode

**Only for the frame-by-frame mode:**

- **Left arrow** key: go to the previous frame
- **Right arrow** key: go to the next frame

If you have a numeric keypad you can use the following keys in alternative:

- The key **/** will allow you to view the previous frame
- The key **\*** will allow you to view the next frame

To return in the VLC mode press again the frame-by-frame button in the toolbar.

## Exporting events data

The coded events can be exported in various format (**Observations > Export ?**):

## Export events

This function will export the events of selected observations in TSV, XLS or ODS formats. These formats are suitable for further analysis.

A	B	C	D	E	F
1 Observation id	demo #1				
2					
3 Media file(s)					
4					
5 Player #1	/home/olivier/crop.avi				
6 Player #1	/home/olivier/crop2.avi				
7					
8 Observation date	2015-Nov-30 10:39:18				
9					
10 Description					
11					
12 Time offset (s)	0				
13					
14 independent variables					
15 variable	value				
16 group size	7				
17 weather conditions	rain				
18					
19 time	subject	code	Modifier	comment	status
20	4.3 Subject #1	eat	salad	test	START
21	10 Subject #1	eat	salad		STOP
22	26.6 Subject #2	eat	meat		START
23	113.988 Subject #2	eat	meat		STOP
24	116.588 Subject #1	jump			POINT
25	118.988 Subject #1	jump			POINT
26	120.863 Subject #2	jump			POINT
27	122.438 Subject #2	jump			POINT
28					

## Export aggregated events

This function will export the events of the selected observations in the following formats:

- **tabular format** (TSV, CSV, XLS, ODS, HTML)
- **SQL** format for populating a SQL database
- **SDIS** format for analysis with the GSEQ program available at <http://www2.gsu.edu/~psyrob/gseq>

The **State events** are paired and duration is available.

Example of tabular export

A	B	C	D	E	F	G	H	I	J
1 Observation id	Observation date	Subject	Behavior	Modifiers	Behavior type	Start	Stop	Comment start	Comment stop
2 demo #1	2015-Nov-30 10:39:18	Subject #1	jump		POINT	116.588	0		
3 demo #1	2015-Nov-30 10:39:18	Subject #1	jump		POINT	118.988	0		
4 demo #1	2015-Nov-30 10:39:18	Subject #1	eat	salad	STATE	4.3	10.0	test	
5 demo #1	2015-Nov-30 10:39:18	Subject #2	jump		POINT	120.863	0		
6 demo #1	2015-Nov-30 10:39:18	Subject #2	jump		POINT	122.438	0		
7 demo #1	2015-Nov-30 10:39:18	Subject #2	eat	meat	STATE	26.6	113.988		
8									
9									
10									

Example of SQL export:

```
CREATE TABLE events (id INTEGER PRIMARY KEY ASC, observation TEXT,
                     date DATE, subject TEXT, behavior TEXT,
                     modifiers TEXT, event_type TEXT, start FLOAT,
                     stop FLOAT, comment_start TEXT,
                     comment_stop TEXT);

INSERT INTO events (observation, date, subject, behavior, modifiers,
                    event_type, start, stop, comment_start, comment_stop) VALUES
('demo #1', '2015-11-30 10:39:18', 'Subj #1', 'jump', '', 'POINT', 116.588, 0, "", ""),
('demo #1', '2015-11-30 10:39:18', 'Subj #1', 'jump', '', 'POINT', 118.988, 0, "", ""),
('demo #1', '2015-11-30 10:39:18', 'Subj #1', 'eat', 'salad', 'STATE', 4.3, 10.0, 'vvv', ""),
('demo #1', '2015-11-30 10:39:18', 'Subj #2', 'jump', '', 'POINT', 120.863, 0, "", ""),
('demo #1', '2015-11-30 10:39:18', 'Subj #2', 'jump', '', 'POINT', 122.438, 0, "", ""),
('demo #1', '2015-11-30 10:39:18', 'Subj #2', 'eat', 'meat', 'STATE', 26.6, 113.988, "", "");
```

## Export events as behavioral strings

Behavioral string can be used with the BSA service: Behavioral Strings Analysis (BSA)

Example:

```
# observation id: demo#1
# observation description:
# Media file name: crop.avi, crop2.avi

Subject #1:
eat|eat|jump|jump

Subject #2:
eat|eat|jump|jump
```

## Export events as Praat TextGrid

Example:


```
File type = "ooTextFile"
Object class = "TextGrid"

xmin = 4.3
xmax = 113.988
tiers? <exists>
size = 2
item []:
    item [1]:
        class = "IntervalTier"
        name = "Subject #1"
        xmin = 4.3
        xmax = 10.0
        intervals: size = 1
        intervals [1]:
            xmin = 4.3
            xmax = 10.0
            text = "eat"
    item [2]:
        class = "IntervalTier"
        name = "Subject #2"
        xmin = 26.6
        xmax = 113.988
        intervals: size = 1
        intervals [1]:
            xmin = 26.6
            xmax = 113.988
            text = "eat"
```

# Tools

## Geometric measurements

Geometric measurements can only be made in frame-by-frame mode. Distances, areas and angles can be measured. Click on **Tools > Geometric measurements** to activate the measurements. The **Measurements window** will be shown:

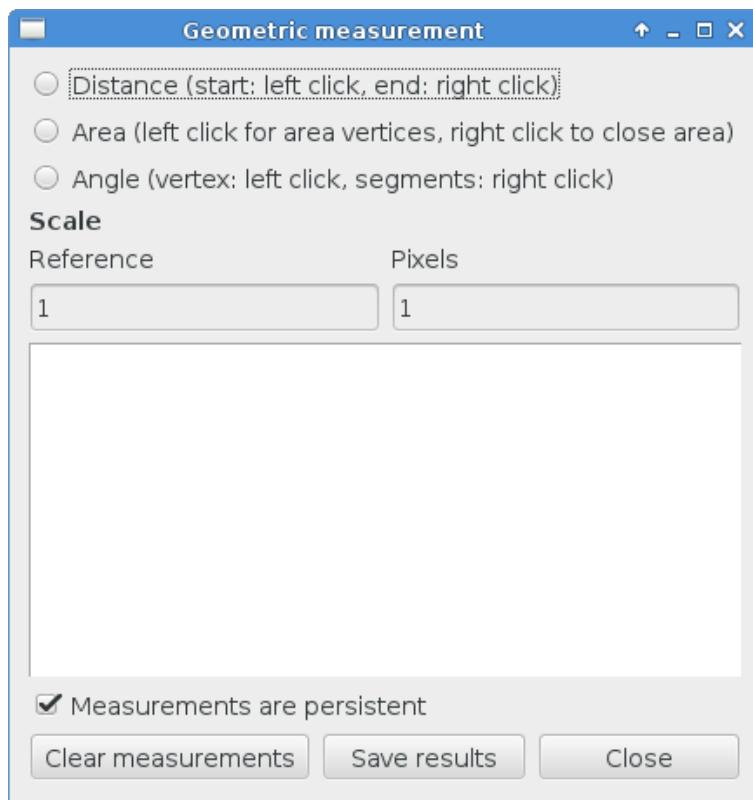


### Setting the scale

For distance and area measurements you can set a scale in order to have results of measurements in a real unit (like centimeters, meters etc). 1) measure a reference object (that have a known size) on the frame (with the distance tool. See next chapter for details) and set the pixel distance in the **Pixel** text box. 2) Set the real size of the reference object in the **Reference** text box (must be a number without unit).

### Distance measurements

Select the **Distance** radio button. Click the left mouse button on the frame bitmap to set the start of the segment that will be measured. A blue circle with a cross will be drawn. Click the right mouse button to set the end. A red circle with a cross will be drawn. The distance between the two selected points will be available in the text area of the **Measurements window**.

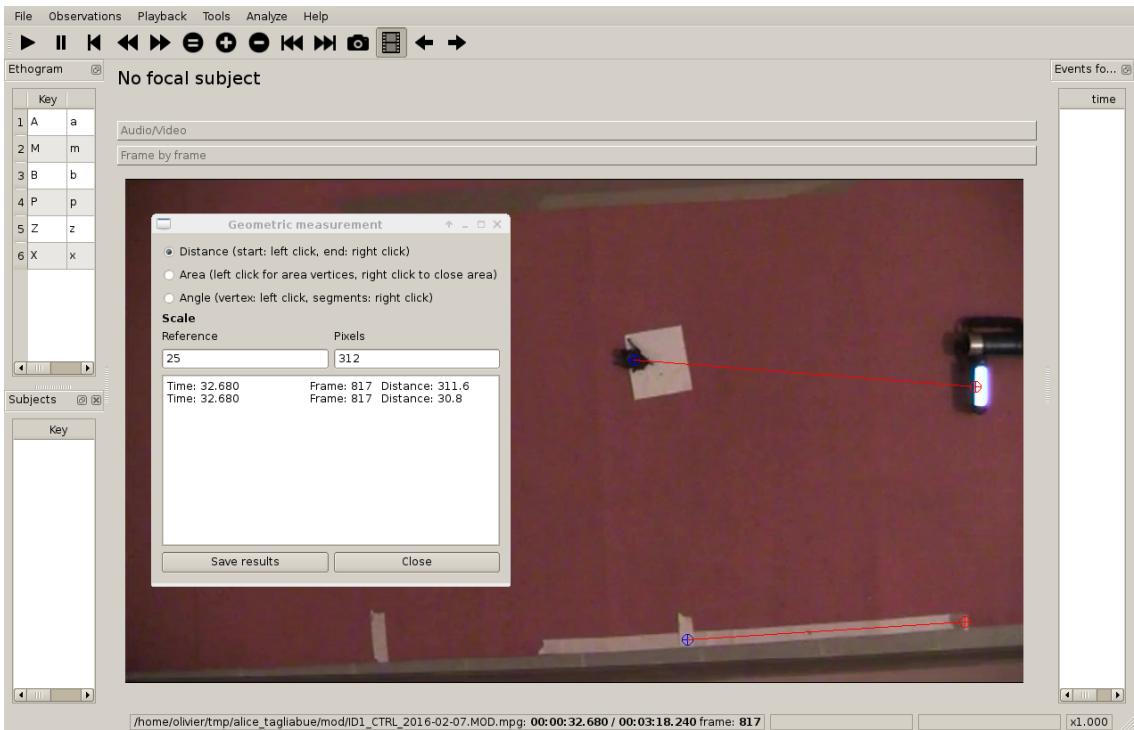

## Area measurements

Select the **Area** radio button. Click the left mouse button on the frame bitmap to set the area vertices. Blue circles with a cross will be drawn. Click the right mouse button to close the area. The area of the drawn polygon will be available in the text area of the **Measurements window**.

## Angle measurements

Select the **Angle** radio button. Click the left mouse button on the frame bitmap to set the angle vertex. A red circle with a cross will be drawn. Click the right mouse button to set the two segments. Blue circles with a cross will be drawn. The angle between the two drawn segments will be available in the text area of the **Measurements window**.

## Persistent measurements

If the **Measurements are persistent** checkbox is checked the measurement schemes will be available on all frames otherwise they will be deleted between frames.

## Extract media file sequences

Sequences of media file corresponding to coded behaviors can be extracted from media file. 1) Click on **Observations > Extract events from media files** option. 2) Choose the observation(s). 3) Select the events to be extracted. 4) Select a destination directory that will contain the extracted media sequences. 5) Select a time offset (in seconds, the default value is 0). The time offset will be subtracted from the starting time of event and added to the stopping time. All the extracted sequences will be saved in the selected directory following the file name format:

{observation id}\_{player}\_{subject}\_{behavior}\_{start time}-{stop time}

# Coding maps

A coding map is a bitmap image with user-defined clickable areas. BORIS allows creating a coding map using the **Map creator** tool (**Tools > Map creator**). Clickable areas may correspond to specific modifiers that can be meaningful for the behavioral coding. Facial expression is the case we thought to when developing this function.

## Creating a coding map

### *Loading a bitmap for a coding map*

To create a new coding map, launch the **Map creator** tool (**Tools > Map creator**). The BORIS main window will be replaced by the **Map creator** window. Click on **Map creator > New Map** and enter a name for the new map in the edit box. You have to load a bitmap image (JPEG or PNG) using the **Load bitmap** button. The loaded image will be displayed.



If the size of your bitmap image is bigger than 512 x 512 pixels BORIS will resize it to 512 x 512 keeping the aspect ratio and store the resized version in the coding map file.

## **Adding areas to a coding map**

To create clickable areas on a coding map, you have to click on the **New area** button and enter an **Area code** in the edit box. The new area can now be defined by clicking on the image. The drawing tool allows defining a irregular polygon (a plane shape with straight sides, which does not have all sides equal and all angles equal) by clicking to determine subsequent vertices. It can be convex or concave. Straight sides must not cross each other. Once selected an area can be deleted using the **Delete area** button. When an area is closed and its name has been defined in the **Area code** field, it can be saved by using the **Save area** button. The areas can partially overlap each other. See the **Using a Coding map** section for more details. Once all areas are added the entire map can be saved using the **Save map** option menu (**Map creator > Save map**). The map is now saved in its own file (.boris\_map) which is NOT part of the BORIS project. A map can be edited at anytime by opening the map file from the **Open map** menu option (**Map creator > Open map**).

## **Adding a coding map to your project**


Creating a Coding map is not automatically adding the map to your project. The Coding map have to be added to your project by selecting the corresponding **Behavior type (Point event with coding map, State event with coding map)**. BORIS will ask to select the file name containing the coding map (.boris\_map) and load the coding map in the project. The coding map name will appear in the **Coding map** column and will be saved in the BORIS project file.

NOTE: If you later modify your coding map you must reload the new version in your BORIS project.

# **Analysis**

## **Time budget**

The **Analyze > Time budget** menu option shows the time budget for the events of one or more observations. You must first select the subjects and behaviors you want to include in the time budget analysis:

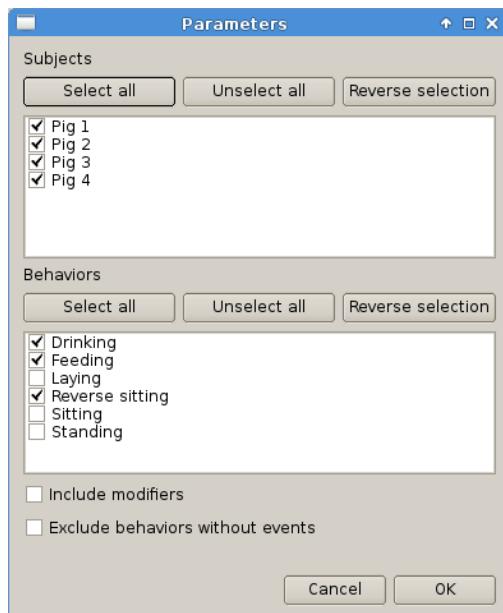

You can choose to select or not the behavior modifiers (if any) and to exclude behaviors without coded events.

Time budget

Selected observations  
DEMO1

Total media length: 00:02:32.240

	Subject	Behavior	Modifiers	Total number	Total duration (s)	Duration mean (s)	Duration std dev	inter-event intervals mean (s)	inter-event intervals std dev	% of total media length
1	Ciro	Huddling	NA	2	77.041	38.521	8.657	45.841	NA	50.6
2	Ciro	Look away	NA	1	-	-	-	NA	NA	-
3	Maurice	Anoint tale	NA	3	-	-	-	39.42	30.971	-
4	Maurice	Flatten ears	NA	1	-	-	-	NA	NA	-
5	Maurice	Flee	NA	1	-	-	-	NA	NA	-
6	Maurice	Looking around	NA	2	7.321	3.661	1.445	29.4	NA	4.8
7	Maurice	Quadrupedal walking	NA	1	5.88	5.88	NA	NA	NA	3.9
8	Maurice	Sitting	NA	1	-	-	-	NA	NA	-
9	Maurice	Stalk	NA	2	-	-	-	32.4	NA	-
10	Maurice	Wave tale	NA	2	-	-	-	35.64	NA	-
11	Totò	Depart	NA	1	-	-	-	NA	NA	-
12	Totò	Hanging	NA	1	27.234	27.234	NA	NA	NA	17.9
13	Totò	Head movement	NA	2	-	-	-	6.96	NA	-
14	Totò	Huddling	NA	2	74.878	37.439	6.449	46.561	NA	49.2
15	Totò	Jumping	NA	2	-	-	-	72.839	NA	-
16	Totò	Look away	NA	2	-	-	-	63.121	NA	-
17	Totò	Passing by	NA	3	-	-	-	44.1	52.524	-
18	Totò	Quadrupedal walking	NA	UNPAIRED	UNPAIRED	UNPAIRED	UNPAIRED	UNPAIRED	UNPAIRED	-

**Save results** **Close**

The percent of total time will be calculated (if available). The total time is intended as the total media length.

The time budget results can be saved in various formats for further analysis: TSV, CSV for plain text file or Microsoft Excel (XLS), Open Document Spreadsheet (ODS).


### Note

If a STATE behavior has an odd number of coded events, BORIS will report "UNPAIRED" instead of results"


## Plot events

The **Analyze > Plot events** menu option plots the events from one observation by subject and behaviours along a time axis.

You must first select the subjects and behaviors you want to include in the plot:

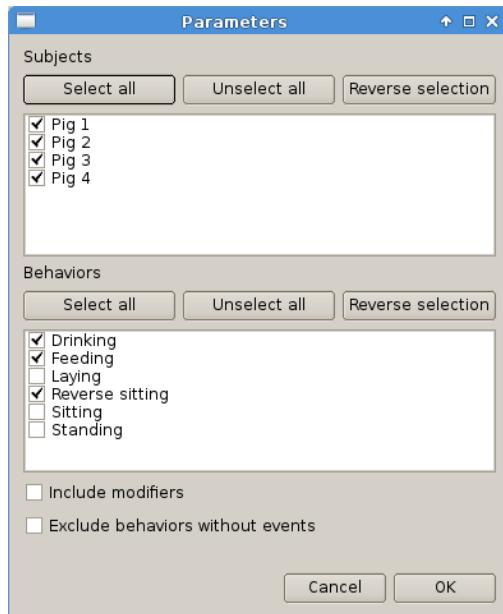

You can choose to select or not the behavior modifiers (if any) and to exclude behaviors without coded events.

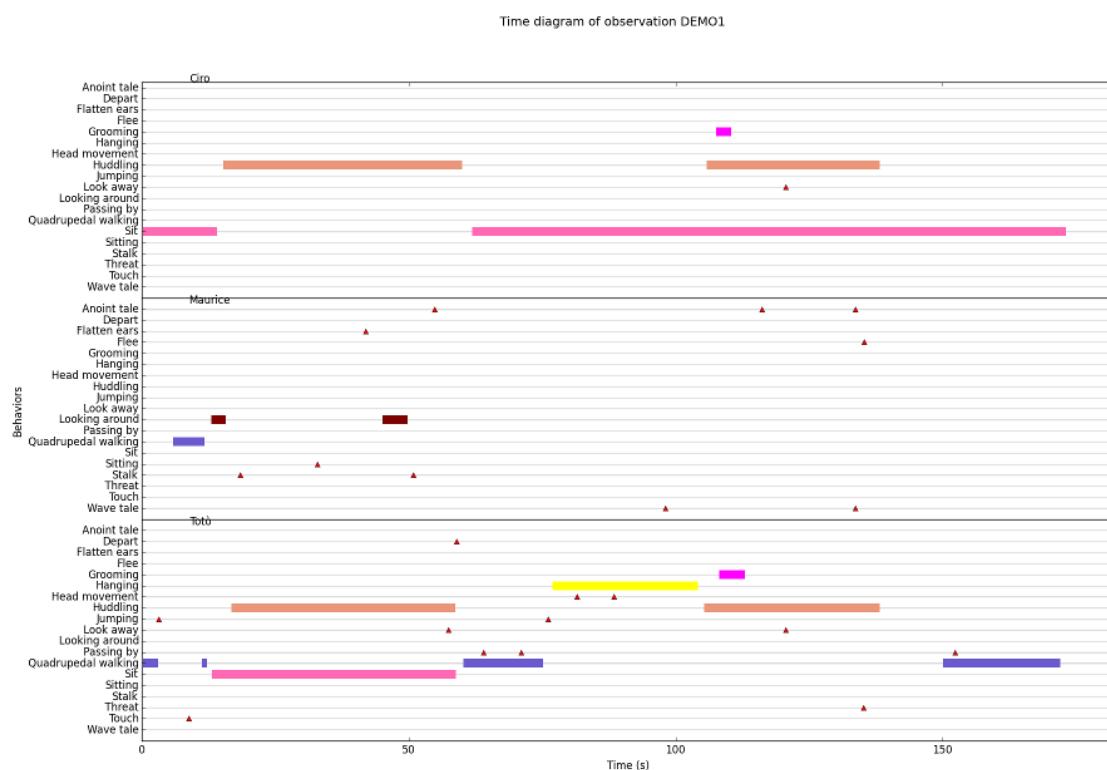

The plot can be exported in various formats like bitmap (PNG, JPG, TIFF) or vectorial graphic (SVG, PDF, EPS, PS). The SVG format can be further edited with the [Inkscape vector graphics editor](#).


### Note

If a STATE behavior has an odd number of coded events, BORIS will give you this error message:  
"The STATE behavior XXX is not paired"

# Preferences

You can customize BORIS using the Preferences window (**File > Preferences**)

## General preferences

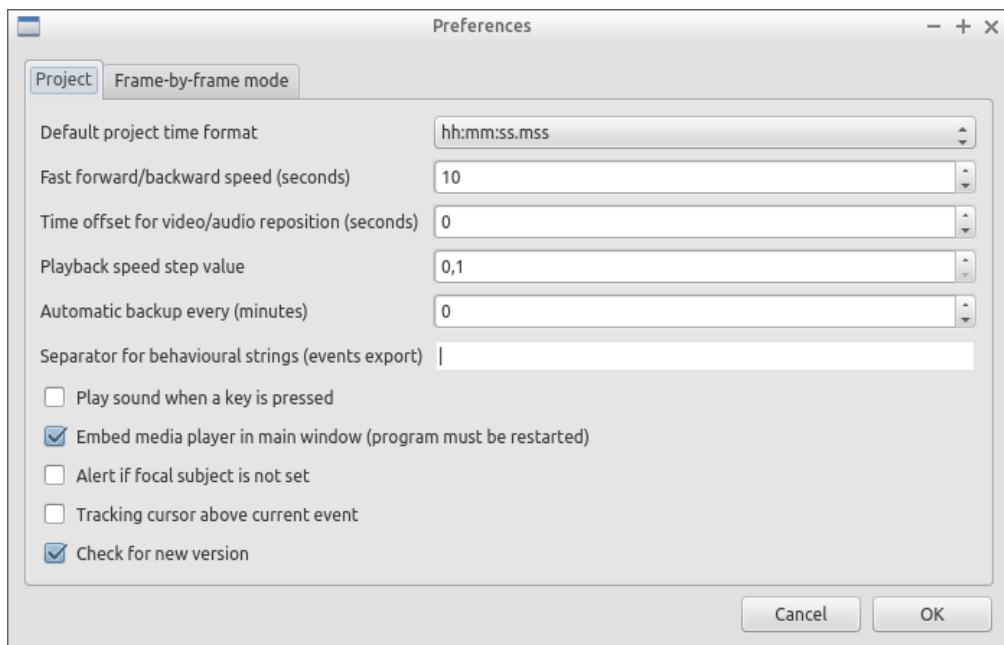

### Default project time format

This option allows the user to choose the format for displaying time in the project. Please note that time is internally always saved in seconds with a precision of 3 decimal digits

### Fast forward/backward speed (seconds)

This option allows the user to customize the amount of time for "jumping" forward or backward in media.

### Time offset for media reposition (seconds)

This value indicates the time offset for repositioning the media after double-click on a row event of the *Events* table. For example -4 seconds indicates that after a double-click the media will be repositioned 4 seconds before the recorded event.'

### Playback speed step value

This value indicate how much the speed will be increased or decreased after pressing the *change playback speed* buttons.

### Automatic backup every (minutes)

If set BORIS will save your project automatically every n minutes. 0 indicate no automatic backup.

### Play sound when a key is pressed

Activate a sound signal after every keypress event

### Embed media player

This option allows the user to detach from the main window or embed the media player in the main window. On Mac OS X the media player can not be detached from main window.

### Alert if focal subject is not set

If this option is activated BORIS will show an alert box if no focal subject is selected


### Tracking cursor above current event

Check this box to position the tracking cursor above the current event in events list table.

## Check for new version

Check for new version on BORIS web site every 15 days (internet access required)

## FFmpeg framework

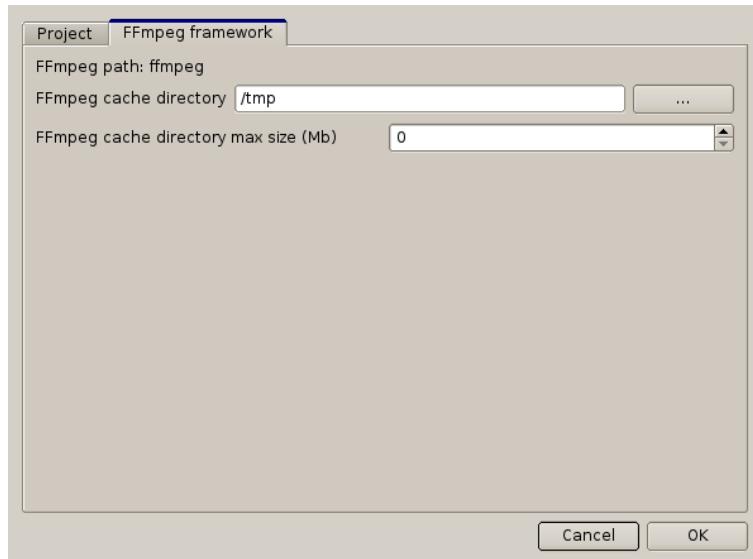

The path for the ffmpeg executable program is displayed. From version 2.8 FFmpeg is included with BORIS.

### FFmpeg cache directory

This indicates the directory that will be used as image cache for frame-by-frame mode and spectrogram visualization. If you do not specify a path, BORIS will use the default temporary directory of your system.

### FFmpeg cache directory max size

Indicate a size limit (in Mb) for the image caching. 0 indicates no limit

## Various

### Citing BORIS

If you have used BORIS for publications, please cite:

Olivier Friard and Marco Gamba  
BORIS: a free, versatile open-source event-logging software for video/audio coding and live observations.  
Methods in Ecology and Evolution. 2016.  
DOI: 10.1111/2041-210X.12584

### Bug reports and features request


Please send bug reports and features request by e-mail (see web site <http://www.boris.unito.it>)

or by using the BORIS GitHub repository (<https://github.com/olivierfriard/BORIS>).

In case of bug report please verify that you are using the last version of BORIS and indicate your operating system, its version and the CPU architecture (32/64 bits). You may also include the BORIS project that gave you an error. Any information you will provide will not be disclosed to any third party.

## Docking

All elements, including the media player can be undocked from the main window and positioned where you prefer (e.g. they can be on the same desktop over one or two screens).

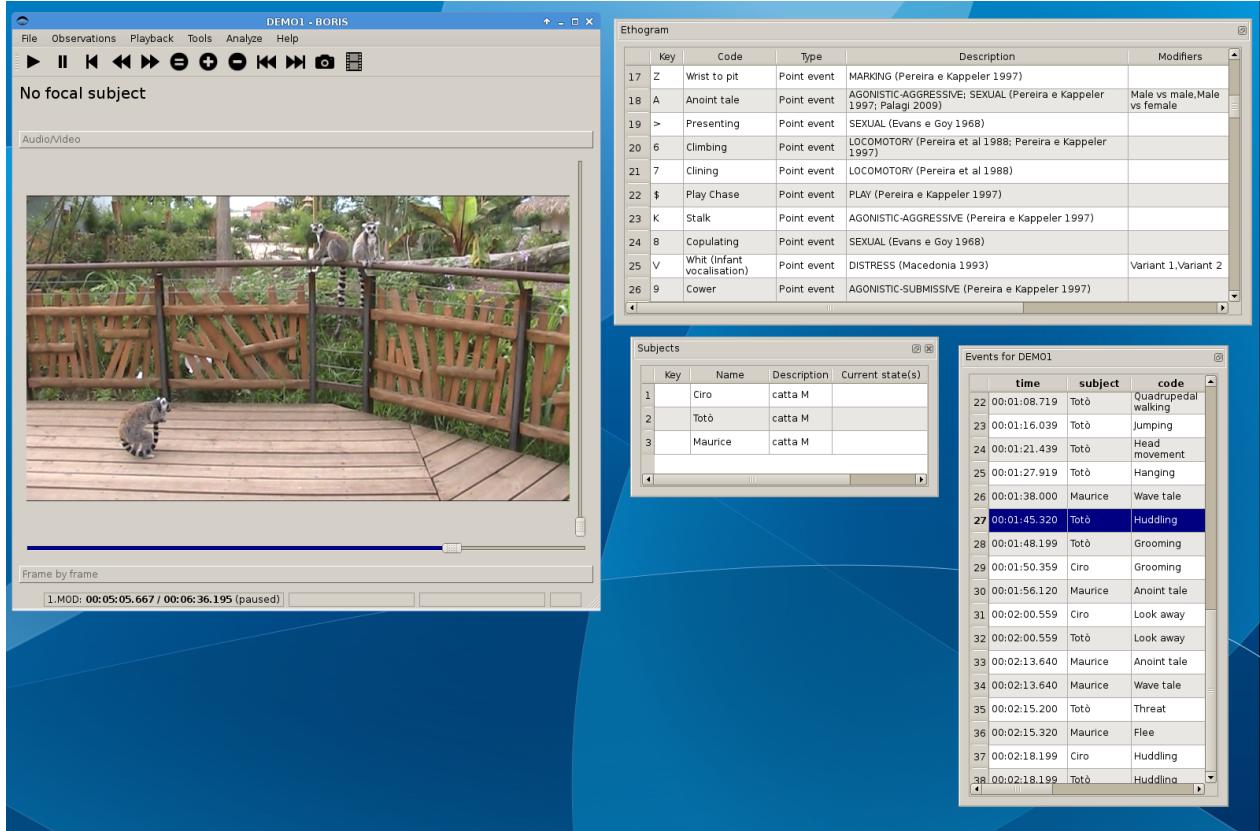

## Acknowledgement

The authors would like to acknowledge Sergio Castellano, Valentina Matteucci and Laura Ozella for their precious help.

## Valid keys for triggering behavior

BORIS will not make difference between lower case and upper case characters

- keys from A to Z
- keys from 0 to 9
- function keys from F1 to F12
- à é è ù ì ç
- ! " £ \$ % & / ( ) = ? ^ [ ] @ | § ° #