

# **Behavioral Observation Research Interactive Software (BORIS) user guide**



BORIS is an easy-to-use event logging software for video/audio coding and live observations.  
BORIS is a free and open-source software available for GNU/Linux, Windows and Mac OS.

**version 7.8**

[www.boris.unito.it](http://www.boris.unito.it)

# Contents

|  |           |
|--|-----------|
| <b>Follow BORIS</b>                        | <b>4</b>  |
| <b>Legal</b>                               | <b>4</b>  |
| <b>Install BORIS</b>                       | <b>4</b>  |
| Linux                                      | 4         |
| Microsoft-Windows                          | 4         |
| Mac OS                                     | 5         |
| <b>Launch BORIS</b>                        | <b>5</b>  |
| <b>Create a new project with BORIS</b>     | <b>6</b>  |
| Set an ethogram                            | 7         |
| Define the subjects                        | 17        |
| Define the Independent variables           | 18        |
| Observations' tab                          | 19        |
| Converters' tab                            | 19        |
| <b>Open an existing project with BORIS</b> | <b>21</b> |
| <b>Import a project</b>                    | <b>22</b> |
| <b>Observations</b>                        | <b>22</b> |
| Create a new observation                   | 22        |
| Observations list                          | 30        |
| Delete observations                        | 33        |
| Import observations from a project file    | 33        |
| <b>Coding your media</b>                   | <b>33</b> |
| Media based coding                         | 33        |
| Live observations                          | 34        |
| <b>Ethogram</b> widget in the main window  | 35        |
| <b>Subjects</b> widget in the main window  | 36        |
| The <b>media player</b> widgets            | 37        |
| The <b>Events</b> widget                   | 38        |
| Record an event                            | 38        |
| Fix unpaired state events                  | 44        |
| Explore project                            | 44        |
| Frame-by-frame mode                        | 44        |
| <b>Exporting events data</b>               | <b>45</b> |
| Export events                              | 45        |
| Export aggregated events                   | 46        |
| Export events for analysis with JWatcher   | 47        |
| Export events as behavioral strings        | 47        |

|  |           |
|--|-----------|
| Export events as Praat TextGrid                                  | 47        |
| Export events as Behavioral Binary Table                         | 48        |
| Extract sequences from media files corresponding to coded events | 49        |
| Extract frames corresponding to coded events                     | 50        |
| Export transitions matrix  | 50        |
| <b>Playback menu</b>   | <b>51</b> |
| Jump   | 51        |
| Zoom   | 51        |
| <b>Tools</b>   | <b>51</b> |
| Focus a video area   | 51        |
| Geometric measurements   | 53        |
| Coding pad   | 54        |
| Subjects pad   | 55        |
| Image overlay on video   | 55        |
| Converters for external data values                              | 56        |
| Transitions flow diagram   | 57        |
| Re-encoding and resizing a video file                            | 58        |
| Rotating a video file  | 59        |
| Project server   | 59        |
| <b>Coding map</b>  | <b>59</b> |
| The Behaviors coding map   | 59        |
| The Modifiers coding map   | 63        |
| <b>Analysis</b>  | <b>64</b> |
| Time budget  | 64        |
| Synthetic time budget  | 67        |
| Plot events  | 67        |
| Inter-rater reliability  | 69        |
| Advanced event filtering   | 69        |
| <b>Preferences</b>   | <b>72</b> |
| General preferences  | 73        |
| FFmpeg framework   | 74        |
| Frame-by-frame mode  | 75        |
| Spectrogram  | 76        |
| Plot colors  | 76        |
| <b>Various</b>   | <b>77</b> |
| Removing path of media files                                     | 77        |
| Citing BORIS   | 77        |
| Docking  | 77        |

|   |           |
|---|-----------|
| Configuration files                       | 79        |
| Bug reports and features request          | 80        |
| <b>Acknowledgement</b>                    | <b>80</b> |
| <b>Valid keys for triggering behavior</b> | <b>80</b> |

## Follow BORIS

BORIS web site: [www.boris.unito.it](http://www.boris.unito.it)

BORIS GitHub repository

BORIS twitter feed

BORIS YouTube channel

## Legal

Copyright 2012-2019 Olivier Friard - Marco Gamba

**BORIS** is free software; you can redistribute it and/or modify it under the terms of the GNU General Public License as published by the Free Software Foundation; either version 2, or any later version.

**BORIS** is distributed in the hope that it will be useful, but WITHOUT ANY WARRANTY; without even the implied warranty of MERCHANTABILITY or FITNESS FOR A PARTICULAR PURPOSE. See the [GNU General Public License](#) for more details.

## Install BORIS

BORIS can be installed following instructions on the BORIS site at [download page](#)

All previous versions of BORIS are available in the **Releases** section of the [GitHub repository](#)

### Linux

The BORIS executable program is available for Linux only for 64-bit systems. If you are using a 32-bit system you have to launch BORIS from sources.

The [VLC media player](#) must be installed on your system.

You can also run BORIS from sources. See the [Linux download page](#) for details about the dependencies to install.

### Microsoft-Windows

BORIS is available for Windows in 32-bit and 64-bit versions. Check your OS to download the right version.

Two versions are available: **BORIS Portable** and **BORIS Setup**

The **BORIS Portable** will NOT install BORIS on your system but will run from a directory (find and launch boris.exe). Use this version if you want to test a new version of BORIS for example. Using the **portable** version more versions can coexist on your system.

The **BORIS Setup** will install BORIS on your system (the default install path is located in your home directory, you do not need administration rights to install BORIS). Installing a newer version of BORIS will replace the currently installed version.

### **BORIS debugging version**

A portable version for debugging is available for Windows 64-bit versions.

Use this version if you want to report bug. In this case you should include screenshots or copy of the terminal output. This version will NOT install BORIS in your system.

You have to extract the archive and launch one of the 3 .bat file:

- the **start\_boris\_no\_log.bat** will display the log file in the command prompt window.
- the **start\_boris\_new\_log.bat** will create a new log file at every launch. The file name will contain a timestamp to avoid overwriting it
- the **start\_boris.bat** will add the log information to the boris.log file

The log files are created in the home directory of the current user:

- /home/USER\_NAME/ for Linux
- C:\Users\USER\_NAME\ for Microsoft Windows
- /Users/USER\_NAME/ for MacOS


## Mac OS

BORIS is available for Mac OS only in 64-bit version (from >= 10.8)

## Launch BORIS

Start BORIS by clicking on its icon. For launching BORIS from sources see the BORIS web site. When launching BORIS for the first time it may take some time to show up. Please be patient!

The main window of BORIS will show up. The toolbar with media player commands is disabled for now.

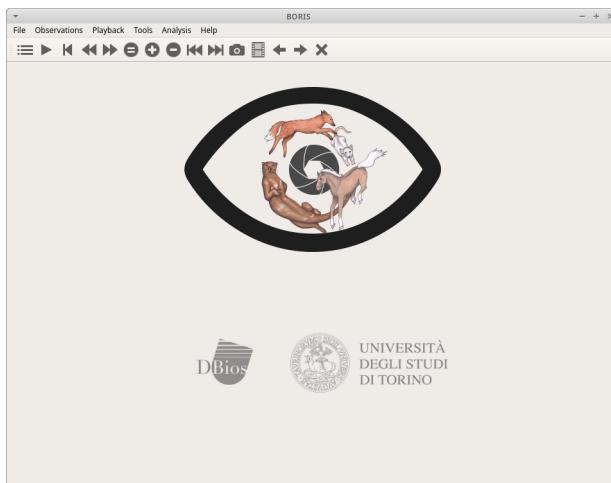

### Warning

On **Mac OS** 10.9 and above, launching can be stopped according to the *Security & Privacy* settings of your computer. They can be changed opening *System Preferences...* > *Security & Privacy* > *General* and selecting the option *Anywhere* in the frame *Allow apps downloaded from*. Alternatively you can repeat the operation right-clicking on the BORIS icon and then clicking *Open* in the following two dialog boxes.

## **Warning**

On **Microsoft-Windows**, launching can be stopped by a *Security warning: The publisher can not be verified. Are you sure you want to run this software?* Click the *Run* button to launch BORIS.

At the first launch, BORIS will ask you to allow the automatically check for new version. Internet access is required for this feature. If you choose **Yes** BORIS will check for new version about every 15 days. This option can be changed on the **Preferences** window (See [general preferences](#))

You can launch BORIS from command line by specifying the project file to open as the first argument.

## **Create a new project with BORIS**

The BORIS project file is the container for all information related to the project. It contains the **ethogram**, the **independent variables**, the **subjects definition**, the **behaviors coding maps**, the **converters** and all **observations** data. The save menu option (**File > Save project** or **Save project as ...**) will save the project in a path in your local file system. You can also activate the [automatic backup](#) feature (see [Preferences](#)).

## **Note**

It is EXTREMELY important to do regular backups of your project files to prevent the lost of data. Software can be reinstalled but your data could quite possibly be lost for ever.

BORIS allows creating an unlimited number of projects but only one project can be opened at a time.

A video tutorial about creating a project is available at <https://www.youtube.com/watch?v=l97Dny5hFOE>

To create a new project, under the menu **File** , select **New project**. You can determine your project name by writing in the **Project name** field in the **Information** tab. Once the project will be saved, the **Project file path** will show the full path to your project file. **Date** will automatically set on the current date and time, but you can alternatively set this info on your media date and time, or whatever you prefer. **Description** can host all the relevant information about your project, can be also left empty. **Time format** can be alternatively set to **seconds** or to **hh:mm:ss.mss**. This choice can be changed at anytime under **File > Preferences** (for MAC users, **BORIS > Preferences**) > **Default project time format**.

**new project**

|   |  |          |                       |              |                      |            |
|---|--|----------|-----------------------|--------------|----------------------|------------|
| Information   | Ethogram   | Subjects | Independent variables | Observations | Behaviors coding map | Converters |
| Project name  | <input type="text"/>   |          |                       |              |                      |            |
| Date  | <input type="text" value="2018-09-07 14:22:49"/>                               |          |                       |              |                      |            |
| Description   | <input type="text"/>   |          |                       |              |                      |            |
| Time format   | <input checked="" type="radio"/> seconds<br><input type="radio"/> hh:mm:ss.mss |          |                       |              |                      |            |
| <input type="button" value="Cancel"/> <input type="button" value="OK"/> |  |          |                       |              |                      |            |

## Set an ethogram

See <https://en.wikipedia.org/wiki/Ethogram> for ethogram definition.

Switching to the **Ethogram** tab, you can alternatively:

- set your ethogram from scratch;
- import an existing ethogram from another BORIS project;
- import an ethogram from a JWatcher global definition file (.gdf).
- [import an ethogram from a plain text file](#)

**edit project**

| Information  | Ethogram | Subjects                 | Independent variables   | Observations | Behaviors coding map | Converters    |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
|--|----------|--------------------------|---|--------------|----------------------|---------------|-----|------|-------------|----------|---------------|---|-------|---|--|---------------|---|-----------|---|--|---------------|---|-------|--|--------------|---------------|---|---------------|---|--|---------------|---|-------|----------------------------|--|---------------|---|----------|---|------------|---------------|---|-----|--|--|---------------|---|-------|---|------------|---------------|---|-----|---------------------|------------|----------------|---|--------------------------|---------------------------------|--|----------------|---|------------|---------------------------------|--|--|--|
| <table border="1"> <thead> <tr> <th>Behavior type</th> <th>Key</th> <th>Code</th> <th>Description</th> <th>Category</th> </tr> </thead> <tbody> <tr><td>1 State event</td><td>A</td><td>Alert</td><td>Otter is stationary and directs its attention towards something or someone (Hasenjager, 2011)</td><td></td></tr> <tr><td>2 State event</td><td>Q</td><td>Allogroom</td><td>Otter licks or scratches with forepaws or hind-paws another river otter's fur (Ethol, 2015)</td><td></td></tr> <tr><td>3 State event</td><td>B</td><td>Breed</td><td>Otters are engaged in copulatory activities (Hasenjager, 2011)</td><td>reproduction</td></tr> <tr><td>4 State event</td><td>C</td><td>Carry objects</td><td>Otter carries objects or food by holding them against the cheek with one front paw while hob...</td><td></td></tr> <tr><td>5 State event</td><td>C</td><td>Chase</td><td>Otter chases other animals</td><td></td></tr> <tr><td>6 State event</td><td>D</td><td>Defecate</td><td>Otter eliminates fecal matter (Ethol, 2015)</td><td>physiology</td></tr> <tr><td>7 State event</td><td>D</td><td>Dig</td><td>Otter uses front legs to move sand, stones on the ground (Ethol, 2015)</td><td></td></tr> <tr><td>8 State event</td><td>D</td><td>Drink</td><td>Otter consumes water found in its environment</td><td>physiology</td></tr> <tr><td>9 State event</td><td>E</td><td>Eat</td><td>Otter consumes food</td><td>physiology</td></tr> <tr><td>10 State event</td><td>I</td><td>Interact with enrichment</td><td>Otter interacts with enrichment</td><td></td></tr> <tr><td>11 State event</td><td>L</td><td>Locomotion</td><td>Otter moves from place to place</td><td></td></tr> </tbody> </table> |          |                          |   |              |                      | Behavior type | Key | Code | Description | Category | 1 State event | A | Alert | Otter is stationary and directs its attention towards something or someone (Hasenjager, 2011) |  | 2 State event | Q | Allogroom | Otter licks or scratches with forepaws or hind-paws another river otter's fur (Ethol, 2015) |  | 3 State event | B | Breed | Otters are engaged in copulatory activities (Hasenjager, 2011) | reproduction | 4 State event | C | Carry objects | Otter carries objects or food by holding them against the cheek with one front paw while hob... |  | 5 State event | C | Chase | Otter chases other animals |  | 6 State event | D | Defecate | Otter eliminates fecal matter (Ethol, 2015) | physiology | 7 State event | D | Dig | Otter uses front legs to move sand, stones on the ground (Ethol, 2015) |  | 8 State event | D | Drink | Otter consumes water found in its environment | physiology | 9 State event | E | Eat | Otter consumes food | physiology | 10 State event | I | Interact with enrichment | Otter interacts with enrichment |  | 11 State event | L | Locomotion | Otter moves from place to place |  | <input type="button" value="Add behavior"/><br><input type="button" value="Clone behavior"/><br><input type="button" value="Remove behavior"/><br><input type="button" value="Remove all behaviors"/><br><input type="button" value="Behavioral categories"/><br><input type="button" value="Convert keys to lower case"/><br><br><input type="button" value="Exclusion matrix"/><br><input type="button" value="Import behaviors from a BORIS project"/><br><input type="button" value="Import from JWatcher"/><br><input type="button" value="Import from text file"/><br><input type="button" value="Export Ethogram"/> |  |
| Behavior type  | Key      | Code                     | Description   | Category     |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
| 1 State event  | A        | Alert                    | Otter is stationary and directs its attention towards something or someone (Hasenjager, 2011)   |              |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
| 2 State event  | Q        | Allogroom                | Otter licks or scratches with forepaws or hind-paws another river otter's fur (Ethol, 2015)     |              |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
| 3 State event  | B        | Breed                    | Otters are engaged in copulatory activities (Hasenjager, 2011)                                  | reproduction |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
| 4 State event  | C        | Carry objects            | Otter carries objects or food by holding them against the cheek with one front paw while hob... |              |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
| 5 State event  | C        | Chase                    | Otter chases other animals  |              |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
| 6 State event  | D        | Defecate                 | Otter eliminates fecal matter (Ethol, 2015)   | physiology   |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
| 7 State event  | D        | Dig                      | Otter uses front legs to move sand, stones on the ground (Ethol, 2015)                          |              |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
| 8 State event  | D        | Drink                    | Otter consumes water found in its environment   | physiology   |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
| 9 State event  | E        | Eat                      | Otter consumes food   | physiology   |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
| 10 State event   | I        | Interact with enrichment | Otter interacts with enrichment   |              |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
| 11 State event   | L        | Locomotion               | Otter moves from place to place   |              |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |
| <input type="button" value="Cancel"/> <input type="button" value="OK"/>  |          |                          |   |              |                      |               |     |      |             |          |               |   |       |   |  |               |   |           |   |  |               |   |       |  |              |               |   |               |   |  |               |   |       |                            |  |               |   |          |   |            |               |   |     |  |  |               |   |       |   |            |               |   |     |                     |            |                |   |                          |                                 |  |                |   |            |                                 |  |  |  |

## **Set your ethogram from scratch**

Clicking on the **Add behavior** button you can add a new row in the **Ethogram** table and behavior type will be automatically set to **Point event**.

The cells with gray background can not be directly edited. You must double-click on them and then select a value.

### **Behavior types**

2 types of behaviors can be defined. Double-click on the cell and select the type of behavior:

- **Point event** behavior when the behavior has **no duration**. The behavior will be coded by pressing the defined keyboard key (see below) or by clicking to the corresponding row in the Ethogram table.
- **State event** behavior when the behavior has a **duration**. The behavior start and stop will be coded by pressing the defined keyboard key (see below) or by clicking to the corresponding row in the Ethogram table. These behaviors **must** have a start event and a stop event.


You can switch between **State event** and **Point event** at your convenience with a double-click on the **Behavior type** cell. You can also add a **Coding map** to either a **State event (State event with coding map)** or a **Point event (Point event with coding map)**; see the "Coding map" section for details).

An existing behavior can be duplicated using the **Clone behavior** button. Its code have then to be changed. On a selected behavior, click on the **Remove behavior** button to remove. The **Remove all behaviors** button will clear the **Ethogram** table. Both the above-mentioned operations must be confirmed when prompted.

Behavioral codes (**Code** column) can be sorted alphabetically by checking the **Alphabetical order** checkbox. Alternatively they can be sorted manually by using the **Move up** and **Move down** buttons.

### **Categories of behaviors**

Defining categories of behaviors can be usefull for the analysis of coded events (for example the time budget analysis). Click the **Behavioral categories** button and add a the categories of behaviors. Behaviors can then be included or not in a defined category.

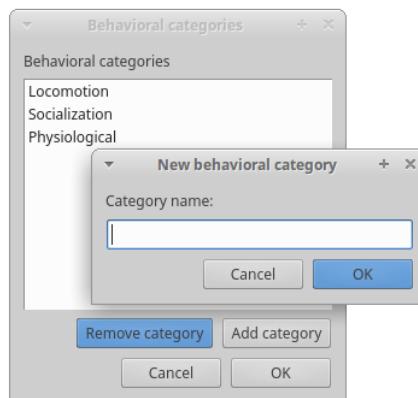

### **Set keys and codes**

For each behavior you have to set a keyboard key (**Key** column) that will be then used to code the behavioral events. You can choose whether you want to set a unique key for each behavior or use the same key for more than one behavior. In the case you set the same key for more than a behavior, BORIS will pause your coding and ask which of the behavior you want to record. From version 7 the keys are **case-sensitive**.

If your project was created with a previous version of BORIS (< v.7) you can use the **Convert keys to lower case** to convert all keys to lower case otherwise you will have to code your observation using upper case key.

### **Important**

If you open a project file created with a version older than v.7 BORIS will ask you to convert the upper case behavior and subject keys to lower case.

### **Important**

**Do not use the / and \* keys! They are reserved for the frame-by-frame mode.**

In the **Code** column, you have to add a unique code for each behavior. Duplicated codes are not accepted and BORIS will warn in red about duplicates on the bottom left of the *Ethogram* tab. The code can be an alphanumeric string (which must not include the pipe character | ).

The **Category** column allow you to include the behavior to a predefined category.

The **Description** of your behavior is optional. The **Description** column can be useful to add information about a specific behavior, its characteristics (e.g. to standardise observation between different users) or to refer to external information (e.g. reference to a previous ethogram).

The following three columns (**Modifiers**, **Exclusion**, **Coding map**) cannot be edited from the *Ethogram* table.

#### **Set the modifiers**


Modifiers can be used to add attributes to a behavior. A single behavior can have two or more modifiers attached (e.g. "play" may have "solitary" or "social" as modifiers). The use of modifiers can be convenient to significantly reduce the number of keys and simplify the behavioral coding.

4 types of modifiers are available: **Single selection**, **Multiple selection**, **Numeric** and **Value from external data file**:


- the **Single selection** type will allow the observer to select only **one** modifier for the current behavior.
- the **Multiple selection** type will allow the observer to select one or more modifiers for the current behavior.
- the **Numeric** type will allow the observer to input a number. For example a distance of interaction.
- the **Value from external data file** type will save the value of a variable from an external data file.

In BORIS modifiers can also be added in different modifier sets [e.g. "play, social" may have a modifier set (#1) for "brothers" and another (#2) for "sisters"]. In the case of using sets of modifiers, you can select one/more modifier for each set.

To add modifiers to a behavior, you need to double-click the **Modifiers** cell corresponding to the behavior you want to add the modifiers to. The following window will show up:

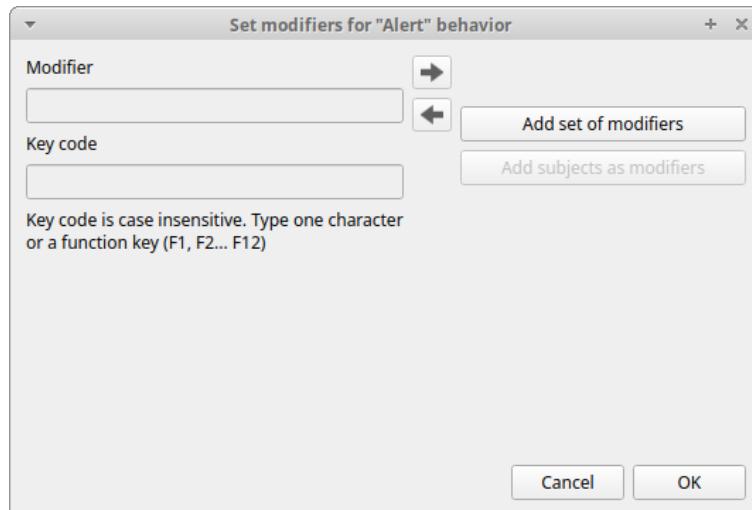

Click the **Add a set of modifiers** button:

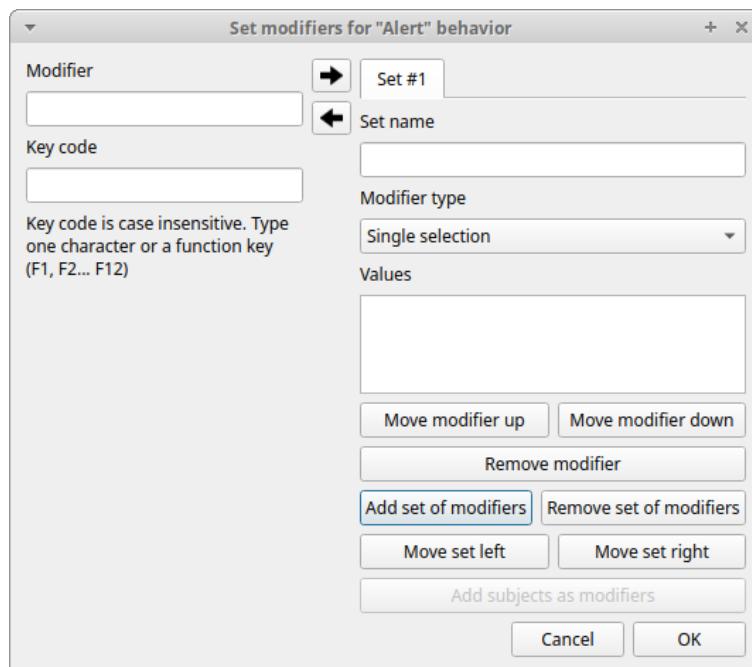


Select the modifier type using the **Modifier type** combo box. You have to choose between **Single selection**, **Multiple selection Numeric** and **Value from external data file**.

#### ***Single selection and Multiple selection modifiers***

Set a name for the new modifiers set by typing it in the **Set name** edit box. Setting a modifiers' set name is not mandatory.


Within a set of modifiers, you can add a modifier by writing the modifier in the **Modifier** edit box. You can choose a shortcut (one character - case sensitive) to this modifier (optional). Then press the **right-arrow** button to add the new modifiers to the set.

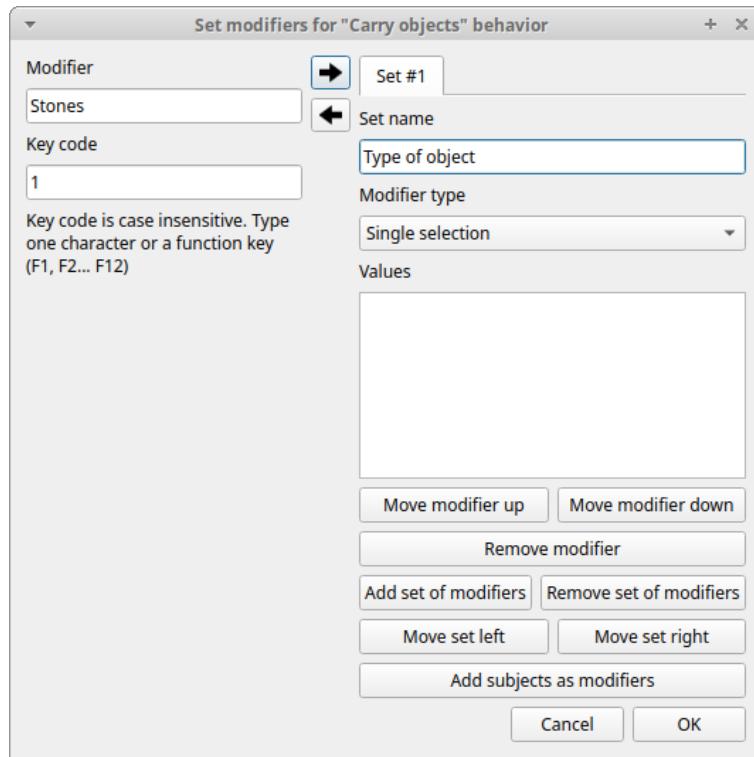

To modify a modifier, select it and press the **left-arrow** button, edit the modifier and press the **right-arrow** button.

A modifier can be removed by pressing the **Remove modifier** button.

After adding all modifiers the window will appear like this:

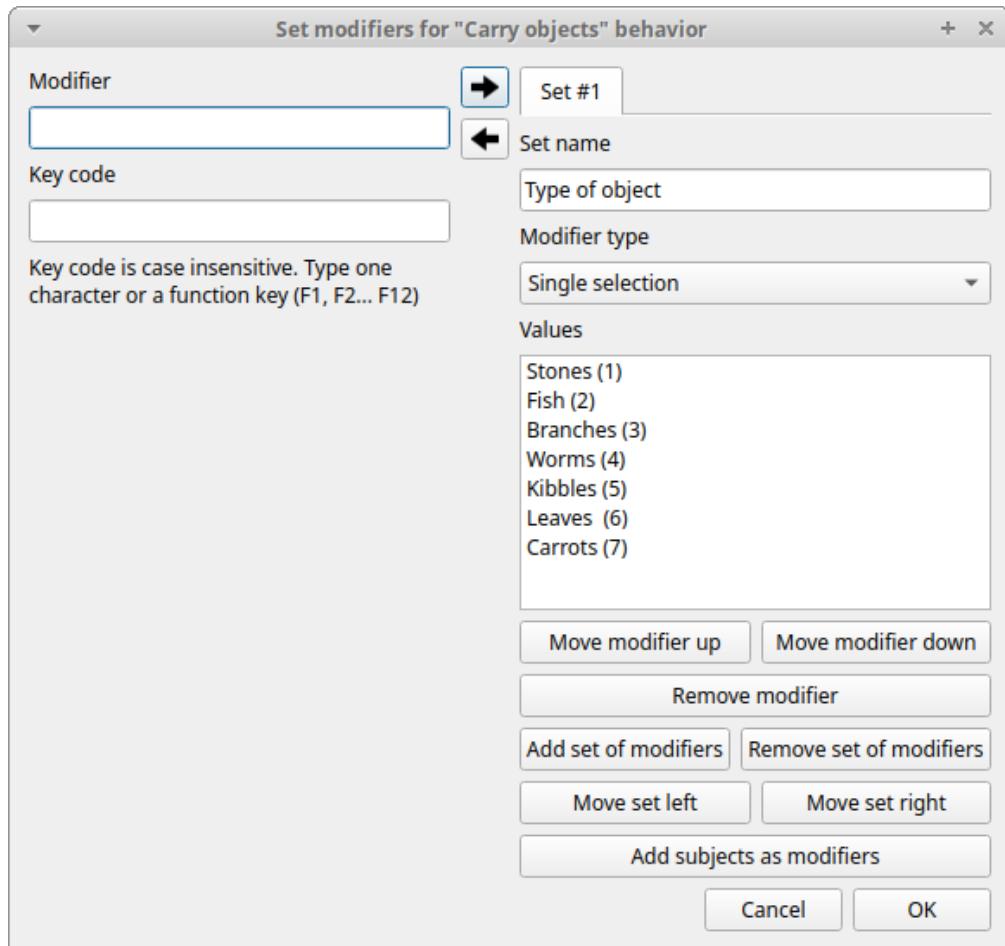

All defined subjects can be added as modifiers using the **Add subjects as modifiers** button. This can help in case of coding the interactions between subjects for example.

The modifiers can be loaded from a plain text file Use the **Load modifiers from file** button.


The modifier position into the modifiers' set can be manually set using the **Move modifier up** and **Move modifier down** buttons. The modifiers can be sorted alphabetically (use the **Sort modifiers** button).

You can add and/or remove sets using the buttons **Add set of modifiers** and **Remove set of modifiers**.

The position of a modifiers' set can be customized (using the **Move set left** and **Move set right** buttons)


Modifiers can not contain the following characters: ((),`~!


Example of a **multiple selection** modifiers set:

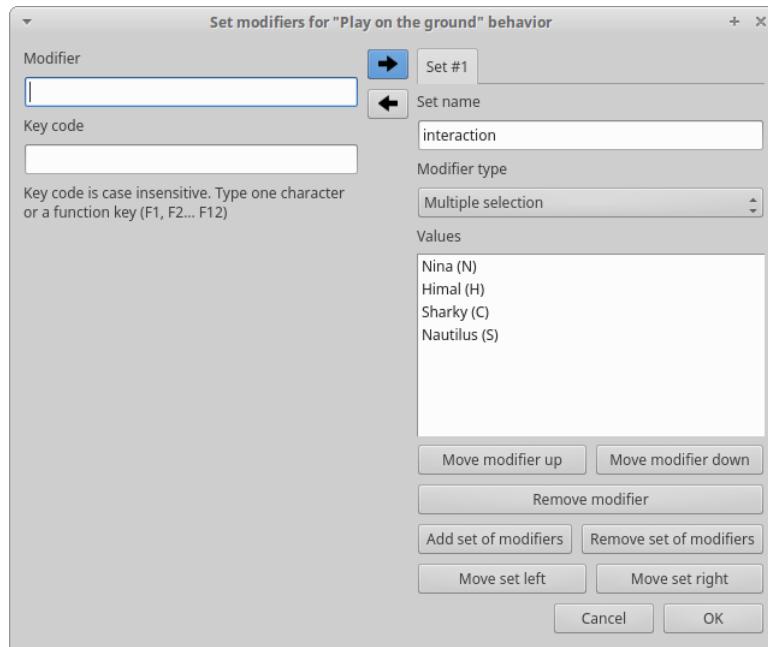

Many values can be selected together.

Example of 2 sets of modifiers:

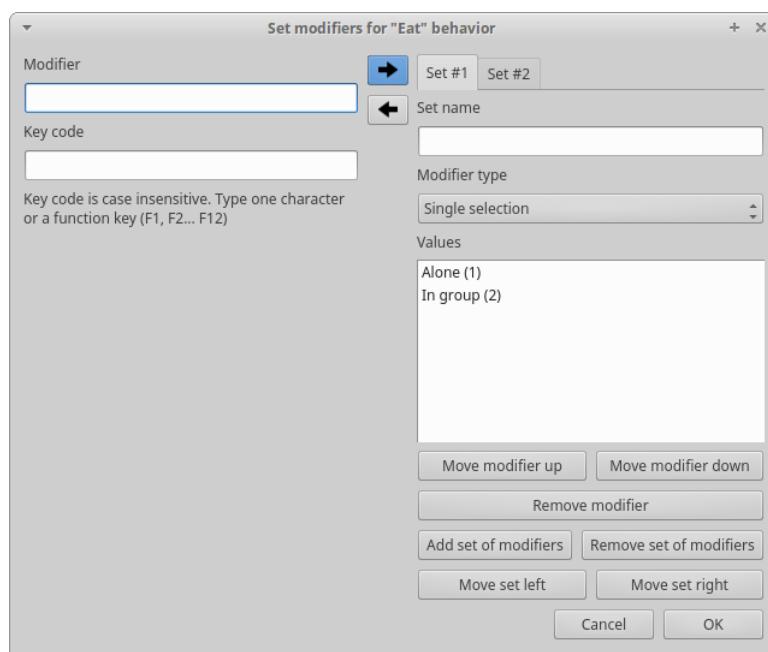

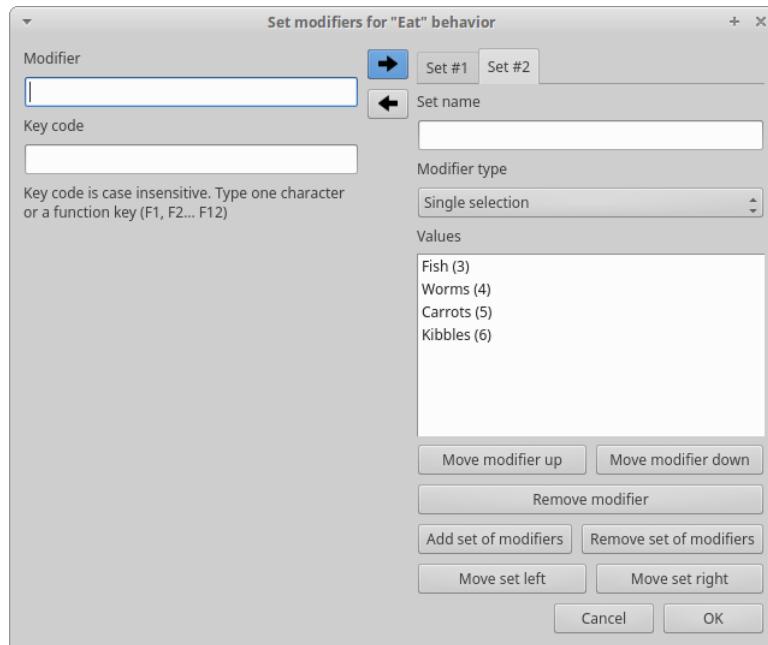

### **Numeric modifier**

Set a name for the new set by typing it in the **Set name** edit box. Setting a modifiers' set name is not mandatory.


When a **Numeric** modifier will trigger, BORIS will ask the observer for a numeric value.

### **Value from external data file modifier**

This modifier can be used to record the value of a variable coming from an external data file (defined during the creation of the observation).

You have to define the variable name in the **Variable name** edit box. This is mandatory and the name of the variable **must** be the same than the variable defined in the observation.

See [External data files](#)

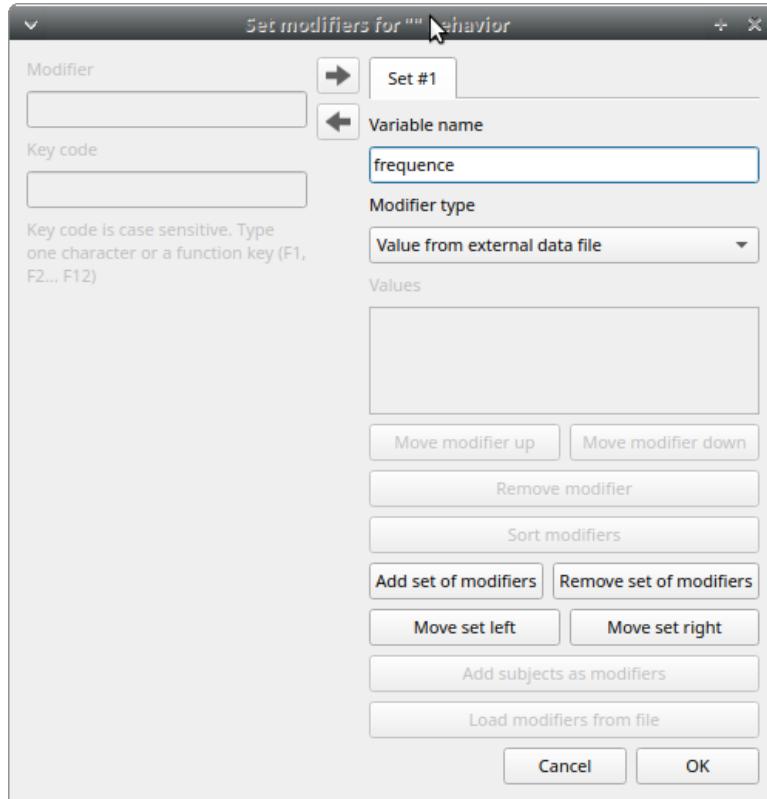

Click **OK** to save modifiers in the **Ethogram** table.

### ***Set the exclusion matrix***

The occurrence of an event (State or Point) can exclude the occurrence of a state event. This can be set using the **Behaviors exclusion matrix** window, which can be opened clicking on the **Exclusion matrix** button. BORIS will ask for including **Point events** or not and a new **Exclusion matrix** window will open.

Exclusive behavior may be selected by checking on the corresponding checkbox in the automatically-generated matrix. We suggest to work on the **Exclusion matrix** when all the behaviors have been added to your ethogram.

All behaviors can be excluded by a particular behavior by selecting the corresponding entire row (click on the row header of the behavior) and by clicking on the **Check selected** button. You can also uncheck all behaviors by selecting the **Uncheck selected** button.

Example of an **exclusion matrix**:

| Behaviors exclusion matrix  |                                     |                                     |                                     |                                     |                                     |                                     |                                     |                                     |                                     |                                     |  |
|-----------------------------|-------------------------------------|-------------------------------------|-------------------------------------|-------------------------------------|-------------------------------------|-------------------------------------|-------------------------------------|-------------------------------------|-------------------------------------|-------------------------------------|--|
| Check behaviors excluded by |                                     |                                     |                                     |                                     |                                     |                                     |                                     |                                     |                                     |                                     |  |
| Alert                       | Allogroom                           | Breed                               | Carry objects                       | Chase                               | Defecate                            | Dig                                 | Drink                               | Eat                                 | Interact with enrichment            | Locomotion                          |  |
| Alert                       | <input checked="" type="checkbox"/> |  |
| Allogroom                   | <input checked="" type="checkbox"/> |                                     | <input checked="" type="checkbox"/> |  |
| Breed                       | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> |                                     | <input checked="" type="checkbox"/> |  |
| Carry objects               | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> |                                     | <input checked="" type="checkbox"/> |  |
| Chase                       | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> |                                     | <input checked="" type="checkbox"/> |  |
| Defecate                    | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> |                                     | <input checked="" type="checkbox"/> |  |
| Dig                         | <input checked="" type="checkbox"/> |                                     | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> |  |
| Drink                       | <input checked="" type="checkbox"/> |                                     | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> |  |
| Eat                         | <input checked="" type="checkbox"/> |                                     | <input checked="" type="checkbox"/> | <input checked="" type="checkbox"/> |  |
| Interact with enrichment    | <input checked="" type="checkbox"/> |                                     | <input checked="" type="checkbox"/> |  |
| Locomotion                  | <input checked="" type="checkbox"/> |                                     |  |
| Look for food               | <input checked="" type="checkbox"/> |  |
| Manipulate                  | <input checked="" type="checkbox"/> |  |
| Play in the water           | <input checked="" type="checkbox"/> |  |
| Play on the ground          | <input checked="" type="checkbox"/> |  |
| Rest                        | <input checked="" type="checkbox"/> |  |

For example in the previous figure, the **Alert** behavior will exclude the following behaviors: **Allogroom, Breed, Carry objects, Chase ...**

During the observation, the excluding event will stop all the current excluded state events one millisecond before the occurrence of the event.

### **Import an ethogram from an existing project**

Behaviors within an ethogram can be imported from an existing BORIS project (.boris) using the **Import behaviours from a BORIS project** button. BORIS will ask to select a BORIS project file and whether imported behaviors should replace or be appended to the **Ethogram** table. Imported behaviors will retain all the previously defined behavior parameters (namely Behavior type, Key, Code, Description, Modifiers and Exclusion information).

### **Import an ethogram from a JWatcher global definition file (.gdf)**

Behaviors can be imported from a JWatcher global definition file (.gdf) using the **Import from JWatcher** button. BORIS will ask to select a JWatcher file (.gdf) and whether imported behaviors should replace or be appended to the **Ethogram** table. Behavior type and exclusion information for the behaviours imported from JWatcher have to be redefined.

### **Import an ethogram from a plain text file**

Behaviors can be imported from a plain text file using the **Import from text file** button. The fields must be separated by TAB, comma (,) or semicolon (;). All rows must contain the same number of fields.

The fields will be interpreted as:

- 1st column: Behavior type **State event** or **Point event** (mandatory)
- 2nd column: Key (one character - case insensitive)
- 3rd column: Behavior code (must be unique)
- 4th column: behavior category (empty if no category)
- 5th column: Description of behavior(optional)

All fields after the 5th will be ignored.

BORIS will ask to select a plain text file (by default: \*.txt \*.csv \*.tsv) and whether imported behaviors should replace or be appended to the **Ethogram** table. The missing information for the behaviours imported from text file have to be redefined.

## **Import an ethogram from a spreadsheet (Google Sheet, Microsoft-Excel or LibreOffice Calc...)**

The ethogram can be imported from a spreadsheet. The spreadsheet must contain one behavior by row and have to be organized as above:

- 1st column: Behavior type: **State event** or **Point event** (mandatory)
- 2nd column: Key (One character - Case sensitive - Optional)
- 3rd column: Behavior code (mandatory - must be unique)
- 4th column: Behavior category (optional)
- 5th column: Description of behavior (optional)

Select all cells of your spreadsheet (CTRL + A), copy to clipboard (CTRL + C). Click the **Import from clipboard** button.

## **Export the ethogram**

The entire ethogram can be exported in various formats (TSV, CSV, XLS, ODS, HTML). See **File > Edit project > Ethogram tab > Export ethogram**

## **Define the subjects**

| edit project |              |  |
|--------------|--------------|--|
| Information  | Ethogram     | Subjects   |
| Key          | Subject name | Description  |
| 1 N          | Nina         | Female, adult, born on 10/03/2013 in Ostrava biopark (Czech Republic), bright white snout          |
| 2 H          | Himal        | Male, adult, born on 04/30/2014 in Amneville biopark (France), larger tale, bigger than the others |
| 3 C          | Sharky       | Male, juvenile, born on 10/30/2015 in Zoom biopark (Italy), bright brown nose and fur              |
| 4 S          | Nautilus     | Male, juvenile, born on 10/30/2015 in Zoom biopark (Italy), dark brown nose and fur                |

BORIS allow coding behaviors for different subjects within a single observation. The **Subject** table allows specifying subjects using a **Key** (e.g. the **k** on your keyboard), **Subject name** (e.g. "Kanzi"), **Description** (e.g. male, born October 28 - 1980).

In this case, pressing **n** will set "Nina" as the focal subject of the behavioural coding. Pressing **n** again will deselect "Nina" and set to "no focal subject".

The definition of one or more subjects is not mandatory. Addition, removal and sorting of the subjects follows the same logic of the *Ethogram* table (see [Set your ethogram from scratch](#) for info).

From version 7 the keys are **case-sensitive**.

If your project was created with a previous version of BORIS (< v.7) you can use the **Convert keys to lower case** to convert all keys to lower case otherwise you will have to code your observation using upper case key.

The subjects can also be imported from an existing BORIS project: use the **Import Subjects from a BORIS project** button.

## Import subject from a spreadsheet

The subjects can be imported from a spreadsheet (Google spreadsheet, Microsoft-Excel, LibreOffice Calc).

The spreadsheet must contain one subject by row and have to be organized as above:

- 1st column: Subject key (One character - Case sensitive - Optional)
- 2nd column: Subject name (mandatory)
- 3rd column: Description of subject (optional)

Select all cells of your spreadsheet (CTRL + A), copy to clipboard (CTRL + C). Click the **Import from clipboard** button.

### Important

If you open a project file created with a version older than v.7 BORIS will ask you to convert the upper case behavior and subject keys to lower case.

## Define the Independent variables

| Label         | Description                            | Type           | Predefined value             | Set of values  |
|---------------|--|----------------|------------------------------|----------------|
| 1 Location    | Location where observations where made | text           | 44° 55' 59" N – 7° 25' 18" E |                |
| 2 Temperature | Average temperature of the day (°C)    | numeric        |                              |                |
| 3 Visitors    | Visitors per day                       | numeric        |                              |                |
| 4 Weather     | Meteorological conditions              | value from set | sun                          | sun,rain,cloud |

Label: Location  
Description: Location where observations where made  
Type: text  
Predefined value: 44° 55' 59" N – 7° 25' 18" E

Add variable  
Remove variable  
Import variables from a BORIS project


BORIS allows adding information about the observation using **Independent variables**. This can be used to specify factors that may influence the behaviors (e.g. group composition, temperature, weather conditions) but will not change during a single observation within a project. Each independent variable can be defined by a **Label** (e.g. weather), a **Description** (e.g. weather conditions), a **Type** (*text*, *numeric*, *value from set* or *timestamp*).

The values of a set are defined in the **Set of values** column separating the available values with a comma (,). Please note that the first value of the set will be selected by default. It should be useful to define a NA value as first value of every set.

The values for the independent variables will be asked when creating a new observation. Addition, removal and sorting of the independent variables follows the same logic of the **Ethogram** table (see **Set**

**your ethogram from scratch** for info). The independent variables can also be imported from an existing BORIS project using the **Import Variables from a BORIS project**.

Example of independent variable defined as "set of values"

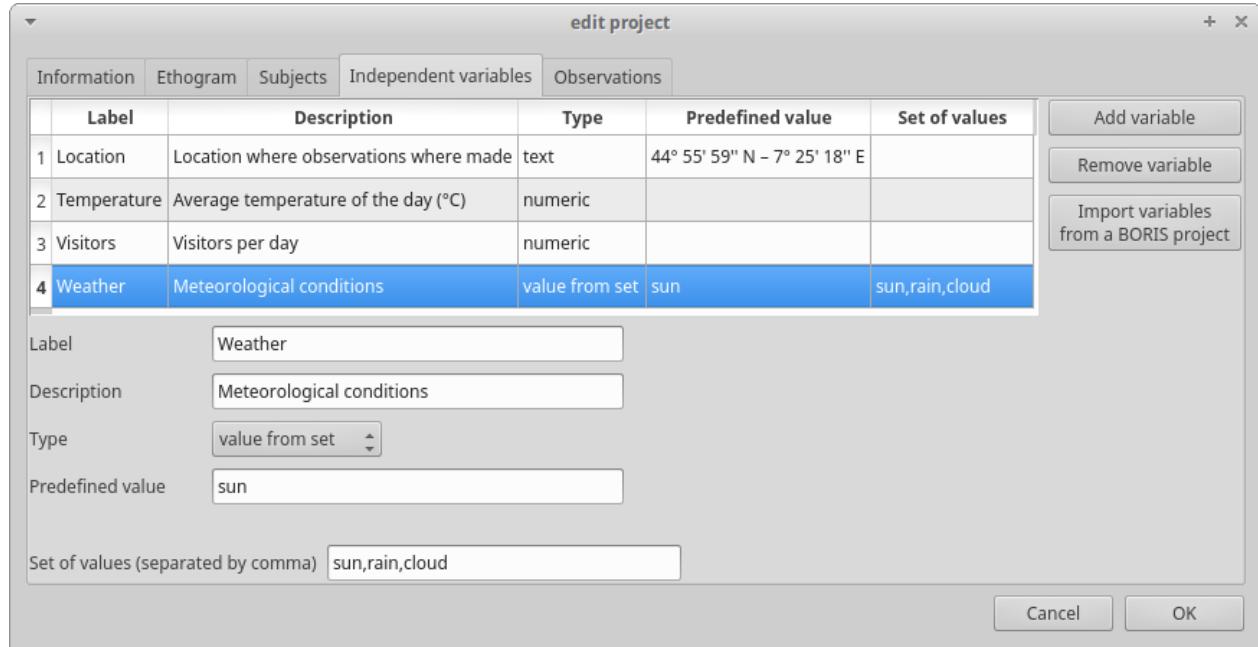


The predefined value must be contained in the set of value.

## Observations' tab

The **Observations** table in BORIS shows information about all the previous observations within a project. A selected "Observation" can be removed using the "Remove observation" button (you will be prompted for confirmation). This operation cannot be undone and deleted observations cannot be recovered once the project is saved. The **Observations** table shows four columns **id Date Description Media**.


## Converters' tab


Converters are used for plotting external data when the timestamp values are not expressed in seconds. Converters can be written by the user, loaded from file or loaded from the repository of the BORIS web site (<http://www.boris.unito.it/static/converters.json>).

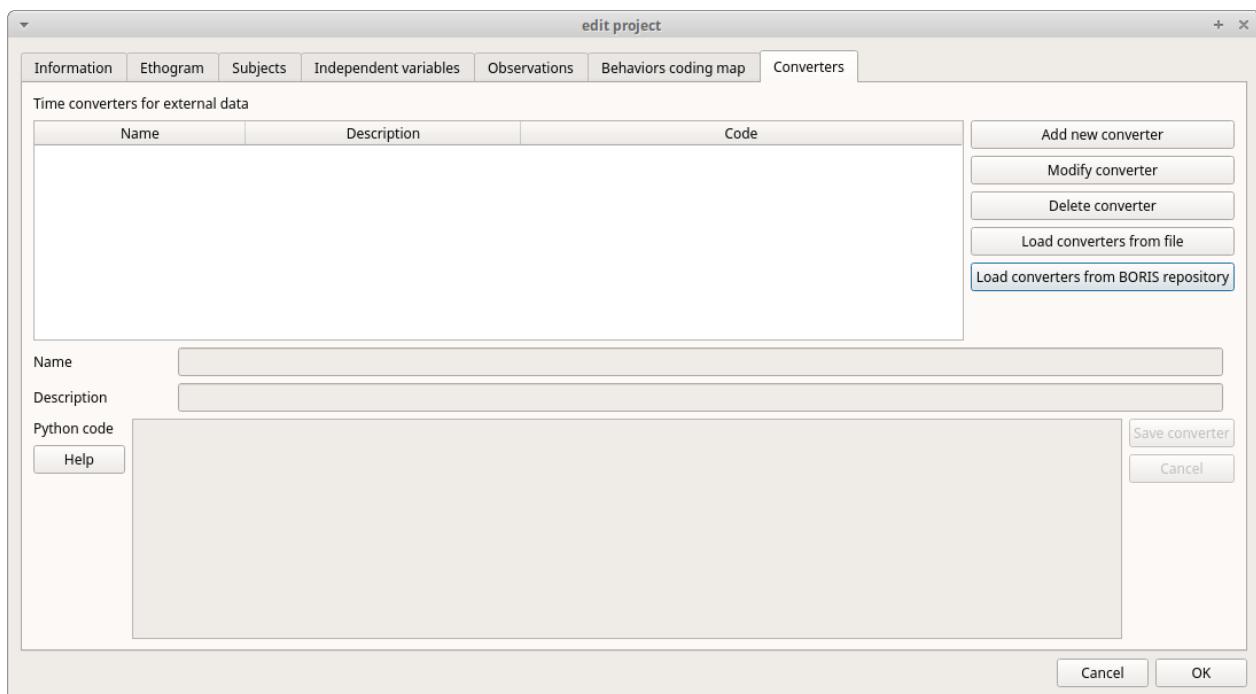

### **Load converters from BORIS web site**

Click **Load converters from BORIS repository** and select the converters to be added to your project.

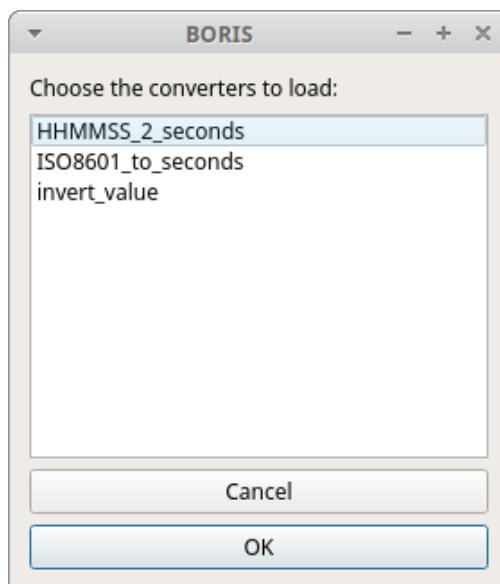

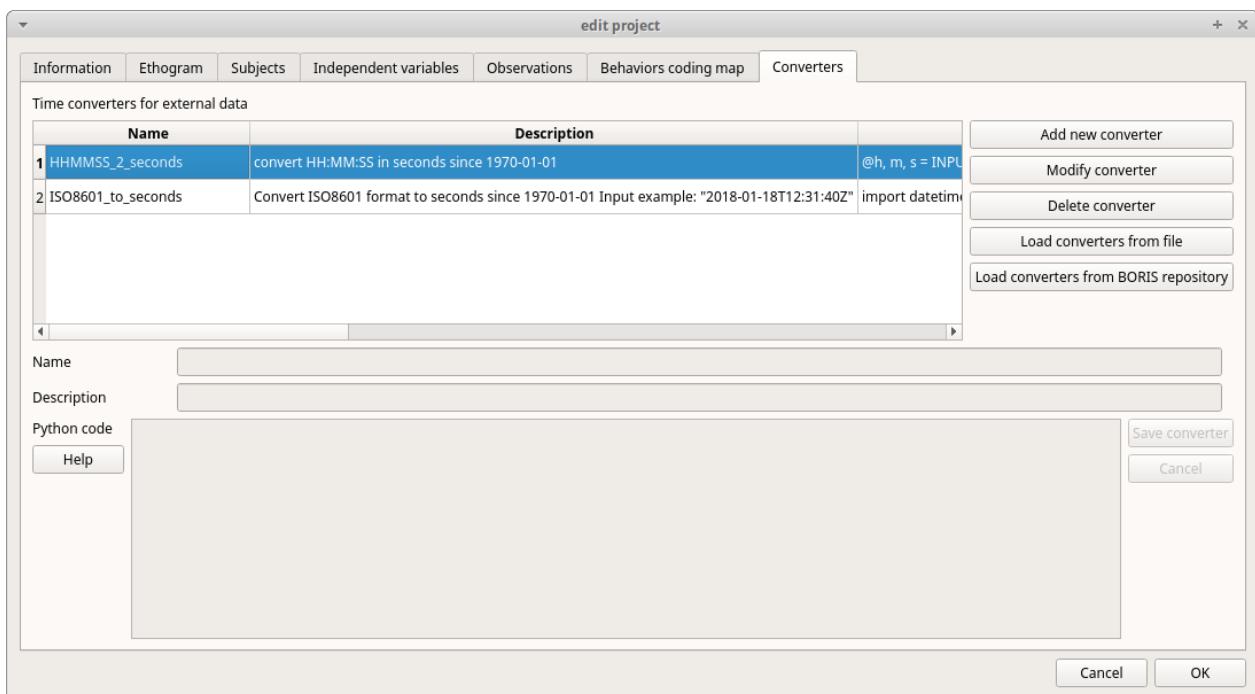

## Writing a converter

See "Converters for external data values"

The converters loaded in your project can be then selected for converting timestamp (or other values) in external data file

See [Converters](#)

## Open an existing project with BORIS

### **Important**

**Due to changes in the project format all the projects created with v. 7+ will not be suitable for previous version of BORIS. The v. 7 will open projects created with previous version and convert them. A copy of your old project will be kept.**

### **Important**

If you open a project file created with a version older than v.7 BORIS will convert your project file and keep an original version of your project file. BORIS will ask you to convert the upper case behavior and subject keys to lower case.

To open an existing BORIS project, select the menu **File > Open project**. A BORIS project file is a container for all information related to a set of observations as the ethogram, the independent variables, and the subjects. BORIS allows creating an unlimited number of projects but only one project can be opened at a time.

# Import a project

You can import a project from a **Noldus The Observer XT Coding Scheme Exchange** file (.otx or .otb extensions). The Observer XT software allows you to save a project as a template (see the reference manual). This template can be imported into BORIS and then saved as a BORIS project (File > Import a project from ... > Noldus The Observer XT template)

The following information will be imported from the template file:

- Ethogram (coding scheme)
- Modifiers
- Behavioral categories
- Subjects
- Independent variables

The Observer is a registered trademark of Noldus Information Technology b.v. See the [Noldus web site for details](#)


# Observations

## Create a new observation

A video tutorial about making an observation is available at <https://www.youtube.com/watch?v=CKeBITCICXc>

To create a new observation you must first [Create a new project with BORIS](#) or [Open an existing project with BORIS](#).

Clicking on **Observations > New observation** will show the **New observation** window.

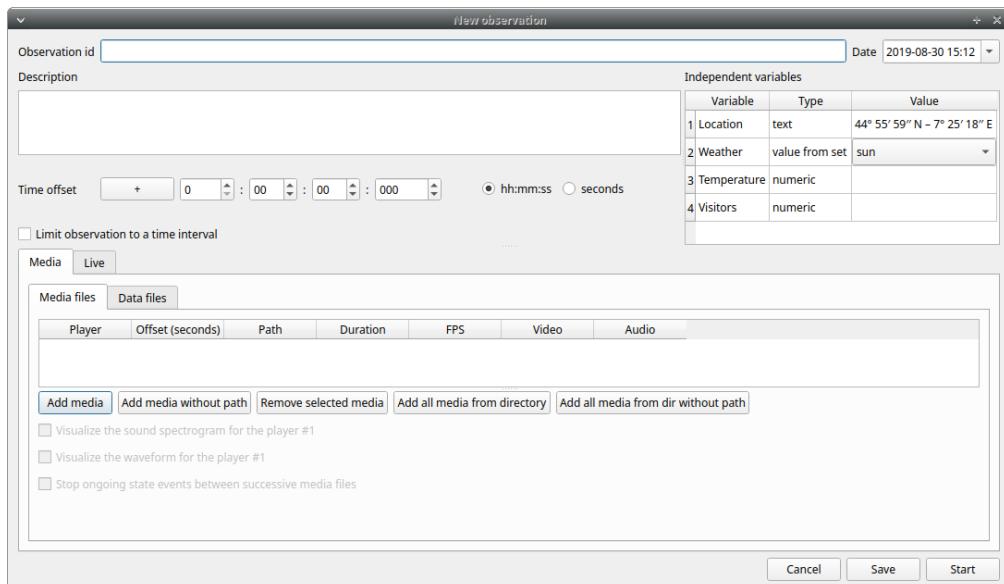

This window allow adding various observation data:

- a mandatory **Observation id** (must be unique across all observations in the open project);
- **Date**, which will be automatically set on the current date and time, but you can alternatively set this info on your media date and time, or whatever you prefer.
- **Description**, which can host all the relevant information about your observation, but can be also left empty.


- **Independent variables** (e.g. to specify factors that may influence the behaviors but will not change during the observation within a project). See the [independent variables](#) section for details.
- **Time offset.** BORIS allow specifying a time offset that can be added or subtracted from the media timecode.
- The **Limit observation to a time interval** option can be used to limit the observation.

You must then indicate if you want to make an observation based on pre-recorded media (audio / video) or a live observation.

## **Live observation**

During the live observation BORIS will show you a timer that will be used for recording time for coded events.

Choose the **Live tab** to make a live observation.

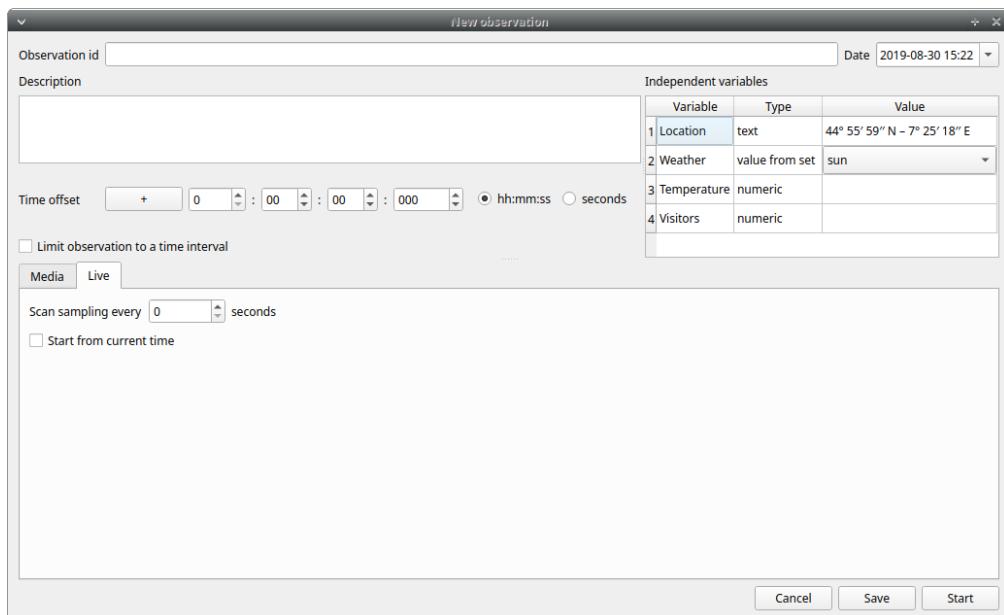

### **Scan sampling**

In the above tab you can select a time for **Scan sampling** observation. In this case the timer will stop at every time offset you indicated and all the coded events will have the same time value.


### **Start from current time**

If you want that the time starts from the current time you can check the **Start from current time** checkbox. The start time will be the computer current time when you will press the **Start** button.

### **Start the observation**

Click the **Start** button to begin the live observation or **Save** to save it in the [Observations list](#).


The main window during a live observation will look like this:

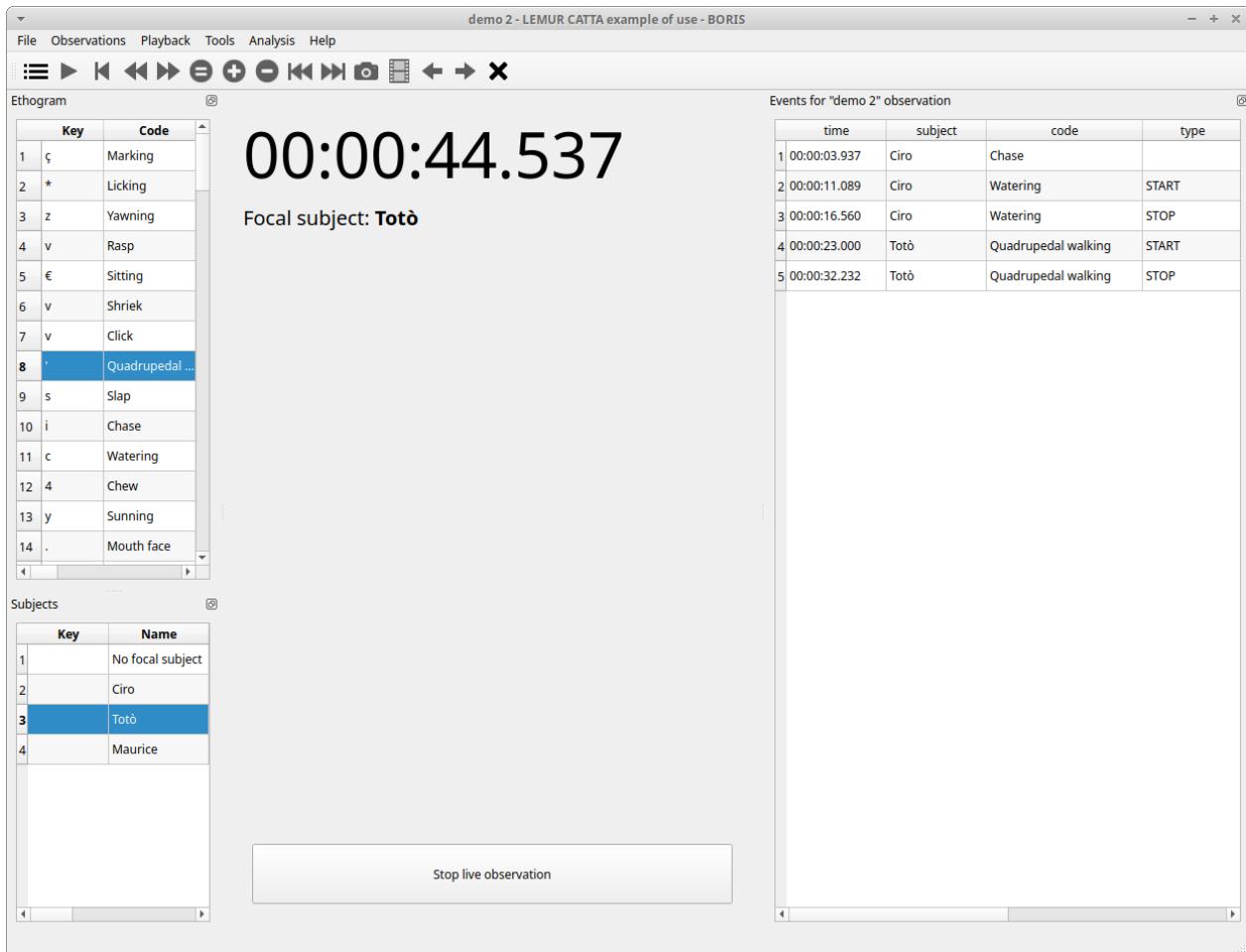

See the [Live observations](#) section to start coding.

## Media based observation

Click on the **Media** tab to make an observation based on one or more media files.

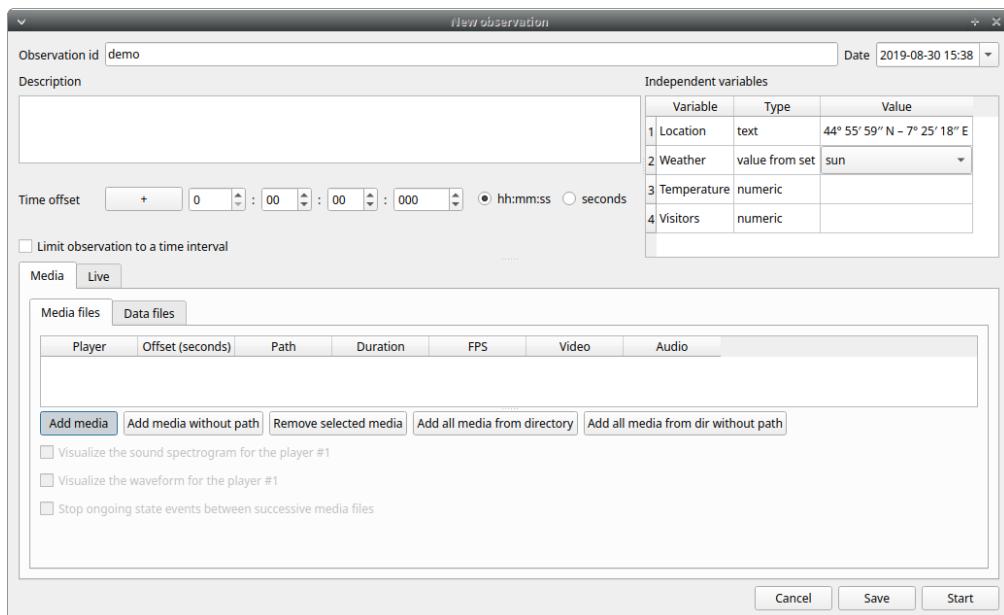

The **Media** tab contains 2 tabs: **Media files** and **Data files**.

Click the **Media files** tab and add one or more media files using the **Add media** or **Add media without path** button. The **Add media without path** button allow you to not store the path of your media file(s) in the BORIS project file. In this case the media file(s) must be in the same directory than the BORIS project file to be found. This option is useful if you are working with various computers on the same project file(s).

If you want to remove the media file paths of an entire project see [Removing path of media files](#).

The **Add all media from directory** and **Add all media from dir without path** allow to add all media files found in a directory to the playlist.

Information about the selected media file will be extracted and displayed in the media list: media file path, media duration, number of frames by second (FPS), the presence of a video stream, the presence of an audio stream .

The screenshot shows the BORIS software interface with the 'Media' tab selected. Below it, the 'Data files' tab is also present. A table displays a single media file entry:

| Player | Offset (seconds) | Path       | Duration     | FPS | Video | Audio |
|--------|------------------|------------|--------------|-----|-------|-------|
| 1   1  | 0                | video1.mp4 | 00:02:59.960 | 25  | True  | False |

Below the table are several buttons: 'Add media', 'Add media without path', 'Remove selected media', 'Add all media from directory', and 'Add all media from dir without path'. There are also two checkboxes at the bottom: 'Visualize spectrogram for the player #1' and 'Stop ongoing state events between successive media files'.

The dropdown list in the first column allow you to choose a player (for a maximum of 8). If you want to observe more media files simultaneously you must use consecutive players (starting from 1). See example below:

The screenshot shows the BORIS software interface with the 'Media' tab selected. Below it, the 'Data files' tab is also present. A table displays three media files assigned to different players:

| Player | Offset (seconds) | Path       | Duration     | FPS | Video | Audio |
|--------|------------------|------------|--------------|-----|-------|-------|
| 1   1  | 0                | video1.mp4 | 00:02:59.960 | 25  | True  | False |
| 2   2  | 0                | video2.mp4 | 00:02:59.960 | 25  | True  | False |
| 3   3  | 0                | video3.mp4 | 00:02:59.960 | 25  | True  | False |

If you have to synchronize 2 (or more) videos you can use the **Offset column** to indicate when the 2nd player should start. For example if the video loaded in the second player starts 15 seconds after the first video you have to input **15** in the **Offset** cell. If the second video starts before the first player you can set a negative value in the **Offset** cell


If you have to play sequentially many videos you have to select the same player (#1) for all video you have loaded. This means that an event occurring at time  $t_x$  in the media file queued as second (e.g. second\_video.mp4) in the playlist will be scored as happening at time  $t_1 + t_x$  (where  $t_1$  is the duration of the first media file, e.g. first\_video.mp4).

The **Remove media** button can be used to remove all the selected media files.


All the media types reported at <http://www.videolan.org/vlc/features.html> can be played in BORIS.

### **Spectrogram visualization**

BORIS allow you to visualize the sound spectrogram during the media observation. Activate the **Visualize spectrogram** check box. BORIS will ask you to generate the spectrograms for all media files loaded in the first player.

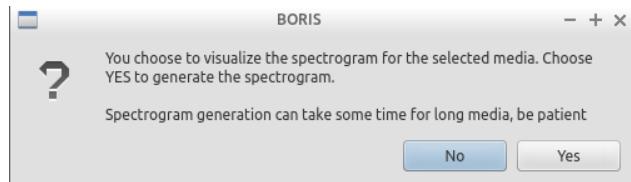

The spectrogram visualization will be synchronized to the media position during the observation.

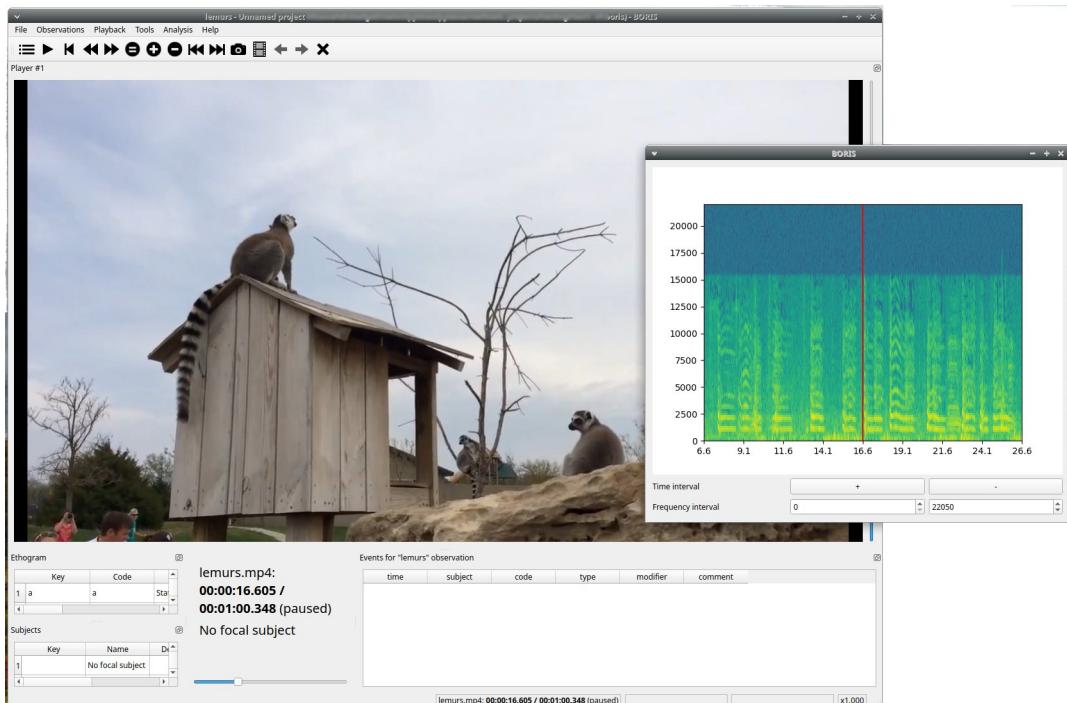

### ***Close current behavior between videos***

This option is disabled for now.

### ***External data files***

**Warning**

**At this time only 2 external data can be plotted with your media file**

You can select one or more external data files to be plotted synchronously with your media. Click the **Data files** tab and use the **Add data file** button to select a data file.

|   |                 |            |               |            |                   |                    |                      |       |
|---|-----------------|------------|---------------|------------|-------------------|--------------------|----------------------|-------|
| Media files   | Data files      |            |               |            |                   |                    |                      |       |
| Data files to plot  |                 |            |               |            |                   |                    |                      |       |
| Path  | Columns to plot | Plot title | Variable name | Converters | Time interval (s) | Start position (s) | Subtract first value | Color |
|   |                 |            |               |            |                   |                    |                      |       |
| <input type="button" value="Add data file"/> <input type="button" value="Add data file without path"/> <input type="button" value="View first rows"/> <input type="button" value="Show plot"/> <input type="button" value="Remove selected data file"/> |                 |            |               |            |                   |                    |                      |       |


The data files must be plain text files with at least **2 columns** separated by comma or TAB. One column must contain a timestamp that will be used to synchronize the plot with the media. The sampling rate can be variable.

Example of a plain text data file with 5 columns separated by comma (,):

```
Display,X Pos,Y Pos,Start Time (secs),Pupil Diameter
1,864,509,549.233,0.00295773451216519
1,863,505,549.25,0.00281810853630304
1,863,503,549.266,0.00287826382555068
1,861,502,549.283,0.0030536837875843
1,858,501,549.3,0.00308083021081984
1,856,499,549.316,0.00306266942061484
1,854,499,549.333,0.00305776367895305
[...]
```

In the above example the 4th column contains the timestamp and the 5th the value to be plotted.

Input the index of the column containing the timestamp and the index of the column containing the value to be plotted. The two indices must be separated by a comma (,). Click **OK** to close the window.

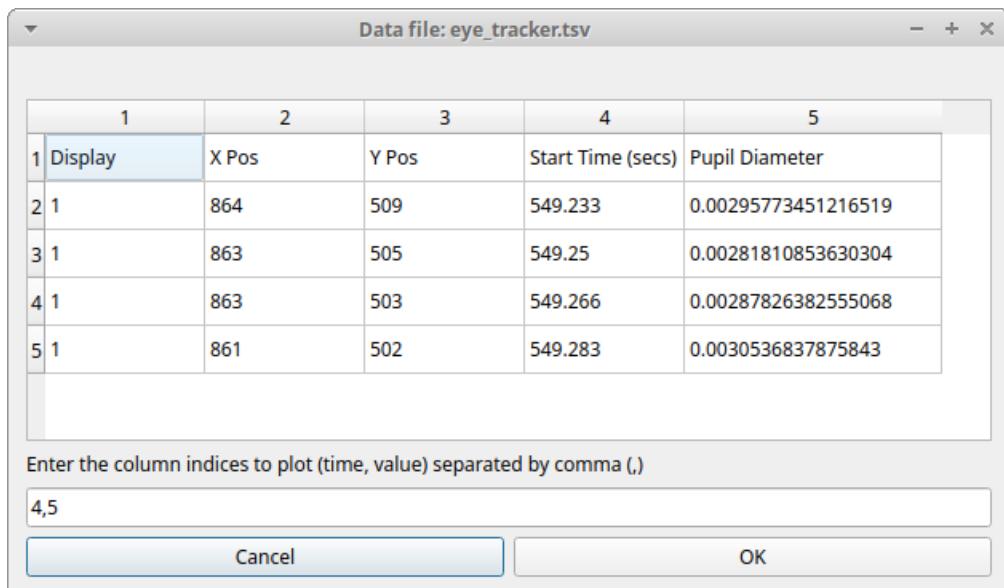

A new row will be added in the data files table.

The screenshot shows the "Data files" tab of a software interface. At the top, there are two tabs: "Media files" and "Data files", with "Data files" being active. Below the tabs, there is a section titled "Data files to plot" with a table. The table has columns for "Path", "Columns to plot", "Plot title", "Variable name", "Converters", "Time interval (s)", "Start position (s)", "Subtract first value", and "Color". There are no rows in the table. At the bottom of this section are four buttons: "Add data file", "Add data file without path", "View first rows", and "Show plot".

You can modify/complete the following parameters by directly typing in the table cells:

- Columns to plot
- Plot title: the title of the plot
- Variable name


- Converters: Used if the timestamp is not expressed in seconds (see below for details)
- Time interval: The time interval that will be plotted (in seconds)
- Start position: the start position of data for synchronisation with the media (in seconds)
- Subtract first value: if the timestamp does not start with a 0 value you can choose to subtract the first value to all timestamp values.
- Color: the color of the color

**NOTE :** if you want to record the value of the plotted variable in a modifier of a behavior (see [Value from external data file modifier](#)) the modifier must have the same **variable name**.

You can check if the data from file can be correctly plotted by using the **Show plot** button. If the data are compatible you will see a plot otherwise you will obtain a message with an explanation.

For now only 2 values can be plotted synchronously with your media file. The values can come from the same file or from two different files.

During the observation the values you have selected in external data files will be plotted synchronously with your media file.

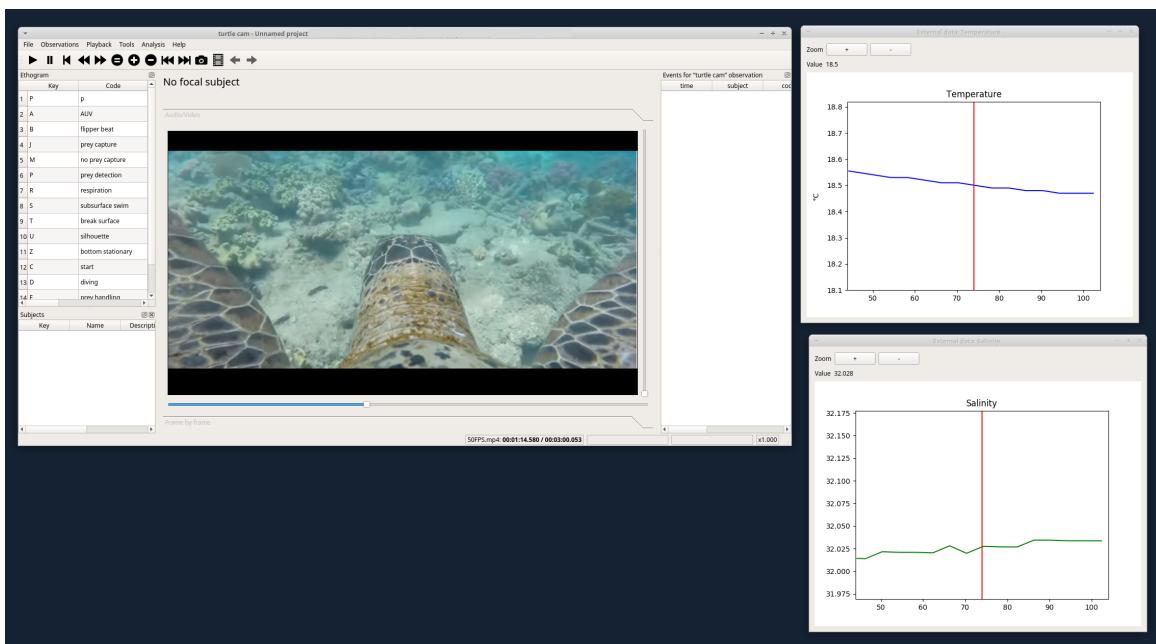


### Converters


If the values in the timestamp column are not expressed in seconds (like 12.45) but in another format (HH:MM:SS, MM:SS, ISO8601 2018-01-18T12:31:40Z ...) you must use a converter that will convert the current format in seconds.

See the Converters' tab in the project configuration: [converters tab](#)

A **double-click** on the converters cell will allow you to select a converter for each column to be plotted

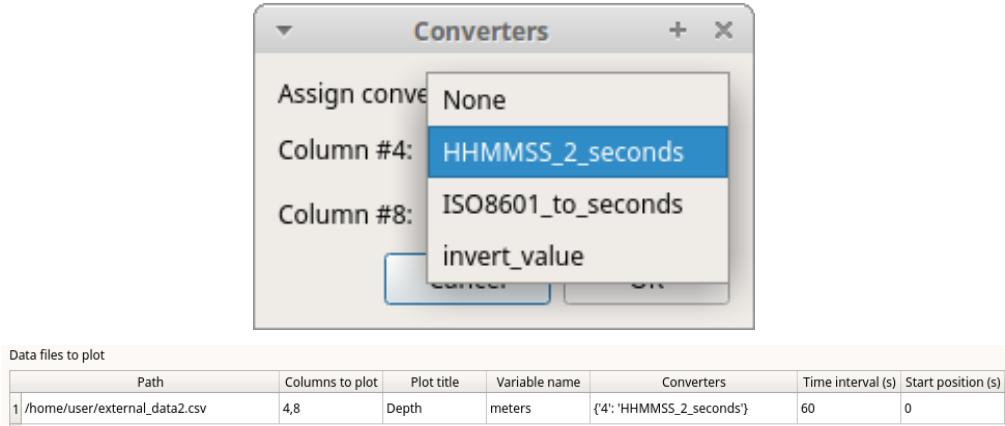

Use the **Show plot** button to verify if your external data can be plotted without problem. The **Close plot** button will close the plot window.

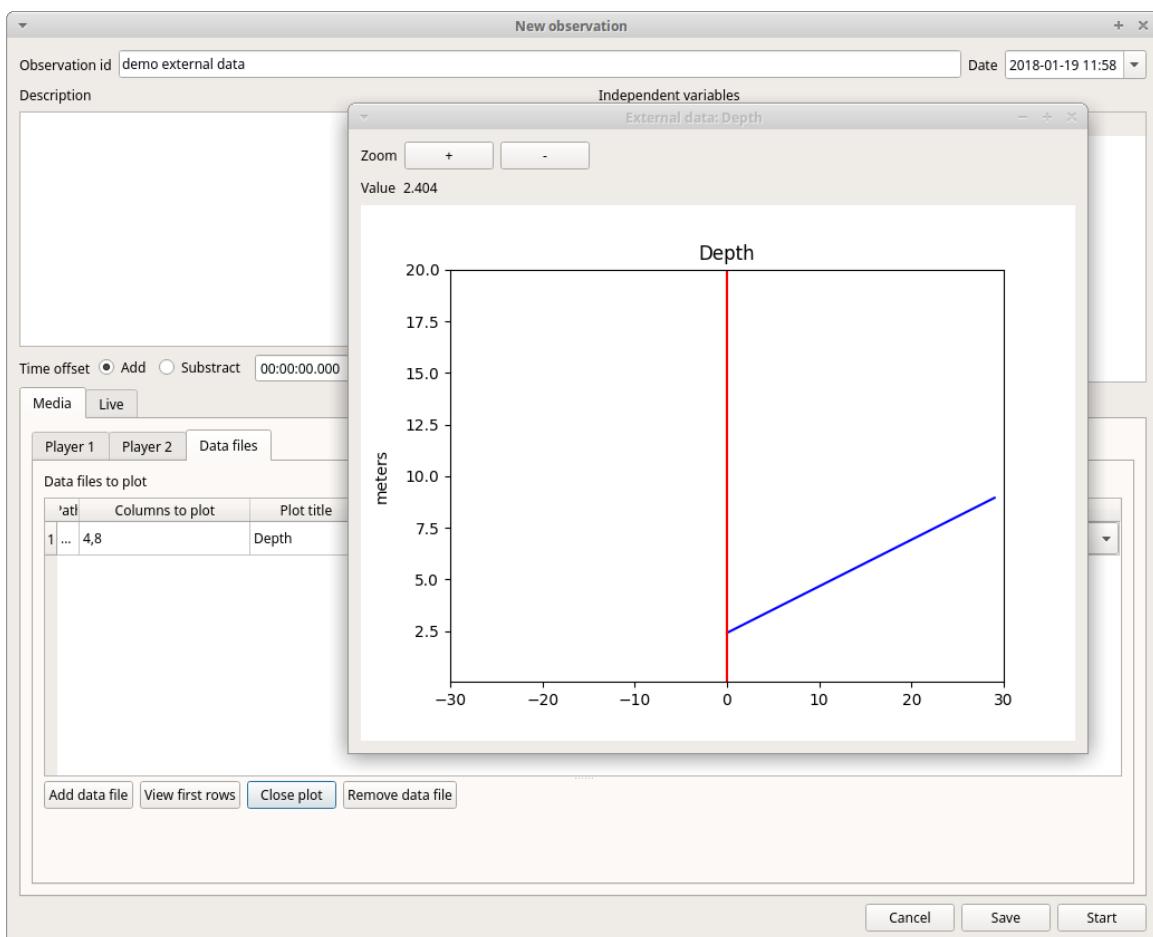

Converters can also be used to convert values that are not time value.


Example of a converter for inverting value:

```
OUTPUT = - float(INPUT)
```

### **Start the observation**

Click the **Start** button to start coding. The **Observation** window will be closed and you'll be transferred to the main **BORIS** window. If you do not want to start the observation click the **Save** button. The observation will be saved in the [observations list](#).


The main window during the observation of a single media file will look like this:

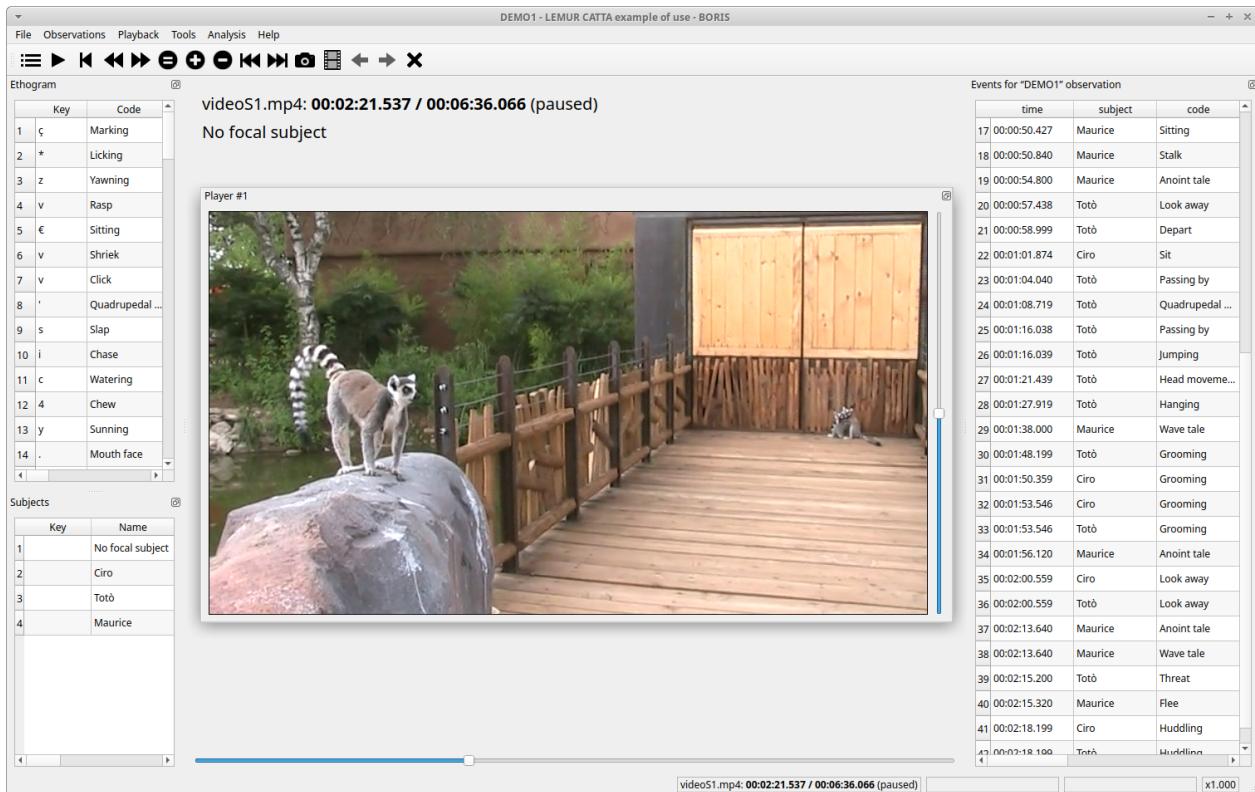

See the [media coding](#) section to start coding.

## ***Limit observation to a time interval***

This option can be used to limit the observation to a time interval for live or media based observations.

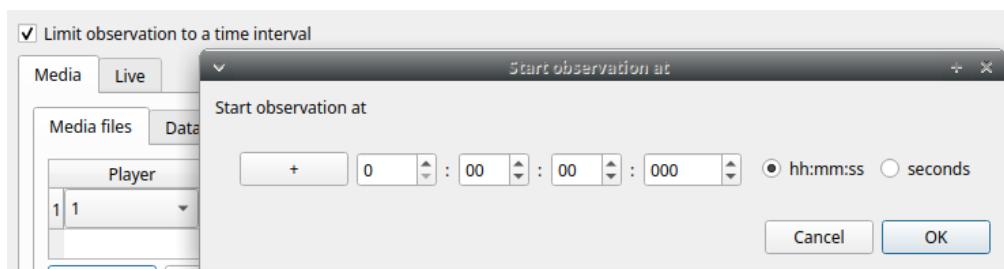

### ***Media based observation***

When the observation will start the media will be automatically positioned to the **Start time** value and the player will stop when the video time will reach the **Stop time** value.

### ***Live observation***

The observation will start at time 0 (or the current time if the **Start from current time** option is enabled). The **Start time** of the time interval will not be applied. The observation will stop when the **Stop time** will be reached.

## **Observations list**

The **Observations > Observations list** will show you all the observations contained in the current BORIS project.

The observations can be sorted by clicking in the desired column header (alphabetic order ascending or descending).

| Observations list - BORIS |                                   |   |                               |                         |                              |         |             |  |
|---------------------------|-----------------------------------|---|-------------------------------|-------------------------|------------------------------|---------|-------------|--|
| 1441 observations         |                                   |   |                               |                         |                              |         |             |  |
|                           | <input type="button" value="id"/> | <input type="button" value="contains"/> |                               |                         |                              |         |             |  |
| id                        | date                              | description                             | subjects                      | media                   | Location                     | Weather | Temperature |  |
| 1 0001_a                  | 2016-05-17 00:00:31               | Vegetation                              | Nautilus, Himal               | #1: 20160517162539.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          |  |
| 2 0001_b                  | 2016-05-17 00:00:31               | Vegetation                              | Nautilus, Himal               | #1: 20160517162539.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 26          |  |
| 3 0002                    | 2016-05-17 00:00:24               | Vegetation                              | Sharky, Nautilus, Himal       | #1: 20160517162540.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 24          |  |
| 4 0003                    | 2016-05-17 00:00:05               | Vegetation                              | Nina, Sharky, Nautilus, Himal | #1: 20160517162641.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          |  |
| 5 0004                    | 2016-05-17 00:00:59               | Central trunks                          | Sharky, Nina, Nautilus, Himal | #1: 20160517162952.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 26          |  |
| 6 0005                    | 2016-05-17 00:00:49               | In the pool                             | Sharky, Nautilus              | #1: 20160517163131.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 24          |  |
| 7 0006                    | 2016-05-17 00:00:42               | In the pool                             | Sharky, Nautilus, Himal       | #1: 20160517163231.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 24          |  |
| 8 0007                    | 2016-05-17 00:00:13               | In the pool                             | Nina, Sharky, Himal           | #1: 20160517163347.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          |  |
| 9 0008                    | 2016-05-17 00:00:17               | In the pool                             | Sharky, Nautilus, Himal       | #1: 20160517163743.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 10 0009                   | 2016-05-17 00:00:10               | In the pool                             | Sharky, Nina, Himal           | #1: 20160517163927.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          |  |
| 11 0010                   | 2016-05-17 00:00:57               | Area near the glass window              | Sharky, Nautilus, Himal       | #1: 20160517164021.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 12 0011                   | 2016-05-17 00:00:50               | Area near the glass window              | Sharky, Nina, Nautilus, Himal | #1: 20160517164106.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 13 0012                   | 2016-05-17 00:00:45               | Area near the glass window              | Sharky, Nina, Nautilus        | #1: 20160517164204.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 26          |  |
| 14 0013                   | 2016-05-17 00:00:25               | Central trunks                          | Nina, Sharky, Nautilus, Himal | #1: 20160517164715.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 15 0014                   | 2016-05-17 00:00:52               | In the pool                             | Nina, Sharky, Nautilus, Himal | #1: 20160517164927.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 16 0015                   | 2016-05-17 00:00:18               | Central trunks                          | Nina, Sharky, Nautilus, Himal | #1: 20160517165631.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 17 0016                   | 2016-05-17 00:00:58               | Central trunks                          | Nina, Sharky, Nautilus, Himal | #1: 20160517165923.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 18 0017                   | 2016-05-17 00:00:37               | Suspended trunks                        | Nina, Sharky, Nautilus, Himal | #1: 20160517170018.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 19 0018                   | 2016-05-17 00:00:18               | Central trunks                          | Nina, Himal                   | #1: 20160517170259.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |

The observations list can be filtered selecting a field and a condition in the drop-list boxes.

In the following example observations are filtered: only observations with **description** containing the **In the pool** subject are shown:

| Observations list - BORIS |  |   |                               |                         |                              |         |             |  |
|---------------------------|--|---|-------------------------------|-------------------------|------------------------------|---------|-------------|--|
| 93 observations           |  |   |                               |                         |                              |         |             |  |
|                           | <input type="button" value="description"/> | <input type="button" value="contains"/> | in the pool                   |                         |                              |         |             |  |
| id                        | date                                       | description                             | subjects                      | media                   | Location                     | Weather | Temperature |  |
| 1 0005                    | 2016-05-17 00:00:49                        | In the pool                             | Sharky, Nautilus              | #1: 20160517163131.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 24          |  |
| 2 0006                    | 2016-05-17 00:00:42                        | In the pool                             | Sharky, Nautilus, Himal       | #1: 20160517163231.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 24          |  |
| 3 0007                    | 2016-05-17 00:00:13                        | In the pool                             | Nina, Sharky, Himal           | #1: 20160517163347.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          |  |
| 4 0008                    | 2016-05-17 00:00:17                        | In the pool                             | Sharky, Nautilus, Himal       | #1: 20160517163743.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 5 0009                    | 2016-05-17 00:00:10                        | In the pool                             | Sharky, Nina, Himal           | #1: 20160517163927.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          |  |
| 6 0014                    | 2016-05-17 00:00:52                        | In the pool                             | Nina, Sharky, Nautilus, Himal | #1: 20160517164927.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 7 0185                    | 2016-05-28 00:00:23                        | In the pool                             | Nina, Sharky, Nautilus        | #1: 20160528124413.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 8 0212                    | 2016-05-31 00:00:53                        | In the pool                             | Sharky                        | #1: 20160531104533.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 9 0217                    | 2016-05-31 00:00:49                        | In the pool                             | Nina, Himal                   | #1: 20160531145324.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 10 0220                   | 2016-05-31 00:00:36                        | In the pool                             | Sharky, Nina, Nautilus, Himal | #1: 20160531153747.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 11 0359                   | 2016-06-06 00:00:54                        | In the pool                             | Nina, Sharky, Nautilus, Himal | #1: 20160606163758.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 12 0411                   | 2016-06-07 00:00:53                        | In the pool                             | Nina, Nautilus                | #1: 20160607103354.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |
| 13 0412                   | 2016-06-07 00:00:21                        | In the pool                             | Nautilus                      | #1: 20160607103630.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |  |

Observations can be filtered with **Independent variables** values.

The following example displays only the observations that do not contain "Sunny" in the **Weather** independent variable :

Observations list - BORIS

1430 observations

Weather does not contain sun

|    | id   | date                | description                | subjects                      | media                   | Location                     | Weather | Temperature |
|----|------|---------------------|----------------------------|-------------------------------|-------------------------|------------------------------|---------|-------------|
| 1  | 0008 | 2016-05-17 00:00:17 | In the pool                | Sharky, Nautilus, Himal       | #1: 20160517163743.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |
| 2  | 0010 | 2016-05-17 00:00:57 | Area near the glass window | Sharky, Nautilus, Himal       | #1: 20160517164021.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |
| 3  | 0011 | 2016-05-17 00:00:50 | Area near the glass window | Sharky, Nina, Nautilus, Himal | #1: 20160517164106.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |
| 4  | 0013 | 2016-05-17 00:00:25 | Central trunks             | Nina, Sharky, Nautilus, Himal | #1: 20160517164715.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |
| 5  | 0014 | 2016-05-17 00:00:52 | In the pool                | Nina, Sharky, Nautilus, Himal | #1: 20160517164927.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |
| 6  | 0015 | 2016-05-17 00:00:18 | Central trunks             | Nina, Sharky, Nautilus, Himal | #1: 20160517165631.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |
| 7  | 0016 | 2016-05-17 00:00:58 | Central trunks             | Nina, Sharky, Nautilus, Himal | #1: 20160517165923.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |
| 8  | 0017 | 2016-05-17 00:00:37 | Suspended trunks           | Nina, Sharky, Nautilus, Himal | #1: 20160517170018.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |
| 9  | 0018 | 2016-05-17 00:00:18 | Central trunks             | Nina, Himal                   | #1: 20160517170259.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |
| 10 | 0019 | 2016-05-17 00:00:23 | Central trunks             | Sharky, Nina, Nautilus, Himal | #1: 20160517170519.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |
| 11 | 0020 | 2016-05-25 00:00:24 | Indoor entrance            | Nina, Himal                   | #1: 20160525145403.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |
| 12 | 0021 | 2016-05-25 00:00:47 | Indoor entrance            | Sharky, Himal                 | #1: 20160525145814.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |
| 13 | 0022 | 2016-05-25 00:00:30 | Indoor entrance            | Nina, Himal                   | #1: 20160525150240.m2ts | 44° 55' 59" N - 7° 25' 18" E |         |             |

Cancel Start View Edit

Observations with a value of **Temperature** independent variable between 18 and 22:

Observations list - BORIS

10 observations

Temperature between (use and to separate terms) 24 and 26

|    | id     | date                | description                | subjects                      | media                   | Location                     | Weather | Temperature |
|----|--------|---------------------|----------------------------|-------------------------------|-------------------------|------------------------------|---------|-------------|
| 1  | 0001_a | 2016-05-17 00:00:31 | Vegetation                 | Nautilus, Himal               | #1: 20160517162539.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          |
| 2  | 0001_b | 2016-05-17 00:00:31 | Vegetation                 | Nautilus, Himal               | #1: 20160517162539.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 26          |
| 3  | 0002   | 2016-05-17 00:00:24 | Vegetation                 | Sharky, Nautilus, Himal       | #1: 20160517162540.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 24          |
| 4  | 0003   | 2016-05-17 00:00:05 | Vegetation                 | Nina, Sharky, Nautilus, Himal | #1: 20160517162641.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          |
| 5  | 0004   | 2016-05-17 00:00:59 | Central trunks             | Sharky, Nina, Nautilus, Himal | #1: 20160517162952.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 26          |
| 6  | 0005   | 2016-05-17 00:00:49 | In the pool                | Sharky, Nautilus              | #1: 20160517163131.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 24          |
| 7  | 0006   | 2016-05-17 00:00:42 | In the pool                | Sharky, Nautilus, Himal       | #1: 20160517163231.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 24          |
| 8  | 0007   | 2016-05-17 00:00:13 | In the pool                | Nina, Sharky, Himal           | #1: 20160517163347.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          |
| 9  | 0009   | 2016-05-17 00:00:10 | In the pool                | Sharky, Nina, Himal           | #1: 20160517163927.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          |
| 10 | 0012   | 2016-05-17 00:00:45 | Area near the glass window | Sharky, Nina, Nautilus        | #1: 20160517164204.m2ts | 44° 55' 59" N - 7° 25' 18" E | sun     | 26          |

Cancel Start View Edit

Observations with a value of **Visitors** independent variable greater than 1000:

| Observations list - BORIS |                     |                            |                               |         |                              |         |             |          |  |
|---------------------------|---------------------|----------------------------|-------------------------------|---------|------------------------------|---------|-------------|----------|--|
| 825 observations          |                     |                            |                               |         |                              |         |             |          |  |
| Visitors                  |                     | >                          |                               | 1000    |                              |         |             |          |  |
| id                        | date                | description                | subjects                      | edi     | Location                     | Weather | Temperature | Visitors |  |
| 4 0003                    | 2016-05-17 00:00:05 | Vegetation                 | Nina, Sharky, Nautilus, Himal | #1: ... | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          | 1046     |  |
| 5 0004                    | 2016-05-17 00:00:59 | Central trunks             | Sharky, Nina, Nautilus, Himal | #1: ... | 44° 55' 59" N - 7° 25' 18" E | sun     | 26          | 1046     |  |
| 6 0005                    | 2016-05-17 00:00:49 | In the pool                | Sharky, Nautilus              | #1: ... | 44° 55' 59" N - 7° 25' 18" E | sun     | 24          | 1046     |  |
| 7 0006                    | 2016-05-17 00:00:42 | In the pool                | Sharky, Nautilus, Himal       | #1: ... | 44° 55' 59" N - 7° 25' 18" E | sun     | 24          | 1046     |  |
| 8 0007                    | 2016-05-17 00:00:13 | In the pool                | Nina, Sharky, Himal           | #1: ... | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          | 1046     |  |
| 9 0008                    | 2016-05-17 00:00:17 | In the pool                | Sharky, Nautilus, Himal       | #1: ... | 44° 55' 59" N - 7° 25' 18" E |         |             | 1046     |  |
| 10 0009                   | 2016-05-17 00:00:10 | In the pool                | Sharky, Nina, Himal           | #1: ... | 44° 55' 59" N - 7° 25' 18" E | sun     | 25          | 1046     |  |
| 11 0010                   | 2016-05-17 00:00:57 | Area near the glass window | Sharky, Nautilus, Himal       | #1: ... | 44° 55' 59" N - 7° 25' 18" E |         |             | 1046     |  |
| 12 0011                   | 2016-05-17 00:00:50 | Area near the glass window | Sharky, Nina, Nautilus, Himal | #1: ... | 44° 55' 59" N - 7° 25' 18" E |         |             | 1046     |  |
| 13 0012                   | 2016-05-17 00:00:45 | Area near the glass window | Sharky, Nina, Nautilus        | #1: ... | 44° 55' 59" N - 7° 25' 18" E | sun     | 26          | 1046     |  |
| 14 0013                   | 2016-05-17 00:00:25 | Central trunks             | Nina, Sharky, Nautilus, Himal | #1: ... | 44° 55' 59" N - 7° 25' 18" E |         |             | 1046     |  |
| 15 0014                   | 2016-05-17 00:00:52 | In the pool                | Nina, Sharky, Nautilus, Himal | #1: ... | 44° 55' 59" N - 7° 25' 18" E |         |             | 1046     |  |
| 16 0015                   | 2016-05-17 00:00:18 | Central trunks             | Nina, Sharky, Nautilus, Himal | #1: ... | 44° 55' 59" N - 7° 25' 18" E |         |             | 1046     |  |
| 17 0016                   | 2016-05-17 00:00:58 | Central trunks             | Nina, Sharky, Nautilus, Himal | #1: ... | 44° 55' 59" N - 7° 25' 18" E |         |             | 1046     |  |

[Cancel](#) [Start](#) [View](#) [Edit](#)

## Delete observations

The observations can be deleted from the project using the following procedure: **File > Edit project > Observation's tab**

Select all observations to remove

Click the **Remove selected observations** button and confirm the deletion.

Please note that the deletion is irreversible.

## Import observations from a project file

The **Observations > Import observations** option allows to import observations from a BORIS project file into the current project. Choose the project file and then the observations to import. BORIS will check if observations with same id are already existing in the current project. BORIS will also check if behaviors and/or subjects used in the imported observations are not defined in the current project.

## Coding your media

When looking at the BORIS main window, the window title bar shows the **Observation id - Project name - BORIS**. The media (the first in the queue) will be loaded in the media player and paused.

## Media based coding

**Media controls are available in the toolbar**



Key to the symbols:

- **List of observations**
- **Play** (become **Pause** when media is played)
- **Rewind** reset your media at the beginning
- **Fast backward** jumps for n seconds backward in your media (See [general preferences](#) to set n)
- **Fast forward** jumps for n seconds forward in your media (See [general preferences](#) to set n)

- **Set the playback speed to 1x**
- **Increase the playback speed** (See [general preferences](#) to set the step value)
- **Decrease the playback speed** (See [general preferences](#) to set the step value)
- **Jump to the previous media file**
- **Jump to the next media file**
- **Take a snapshot** of current video or frame
- **Switch between VLC and frame-by-frame modes**
- **Move on frame back** (enabled in frame-by-frame mode)
- **Move one frame forward** (enabled in frame-by-frame mode)
- **Close current observation**
- **Plot events of current observation**
- **Find in events**

The media can be controlled by special keyboard keys:

- **Page Up** key: switch to the next media
- **Page Down** key: switch to the previous media
- **Up arrow** key: jump forward in the current media
- **Down arrow** key: jump backward in the current media
- **Esc**: switch between VLC and frame-by-frame mode
- **Home** key: Increase the playback speed (See [general preferences](#) to set the step value)
- **End** key: Decrease the playback speed (See [general preferences](#) to set the step value)
- **Backspace**: Set the playback speed to 1x

Only for frame-by-frame mode:

- **Left arrow** key: go to the previous frame
- **Right arrow** key: go to the next frame


## Live observations

During a live observation the media control toolbar is disabled.

Press the **Start live observation** button to start your observation. If some events are already coded BORIS will ask you for deleting them.

A timer will be displayed. The events will be recorded in the events widget.

See the following screenshot.

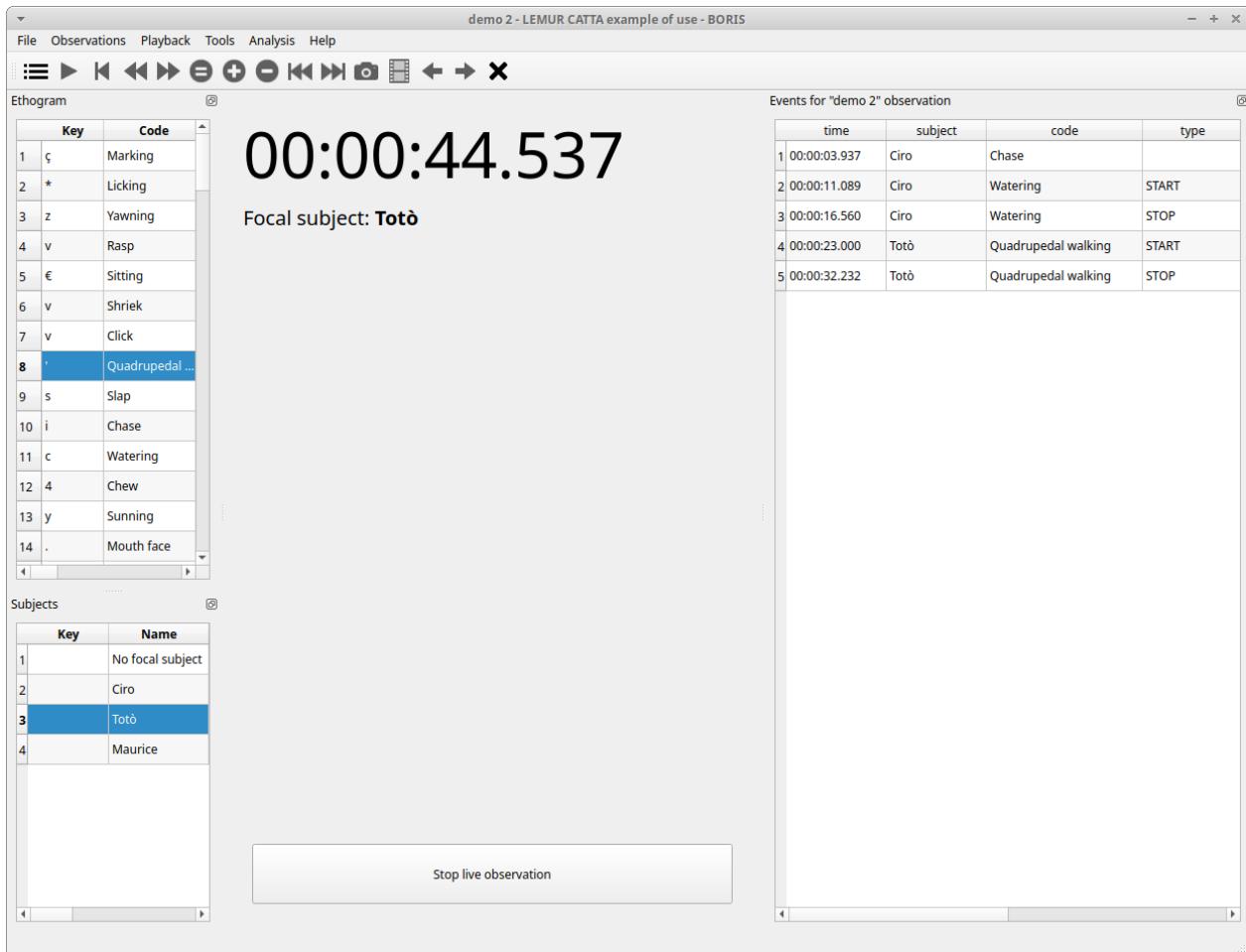

## Ethogram widget in the main window


| Ethogram |      |                    |             |   |                          |                              |
|----------|------|--------------------|-------------|---|--------------------------|------------------------------|
| Key      | Code | Type               | Description | Category  | Modifiers                | Excluded                     |
| 1        | A    | Alert              | State event | Otter is stationary and directs its attention to... | {}                       | Allogroom,Breed,Carry ob...  |
| 2        | Q    | Allogroom          | State event | Otter licks or scratches with forepaws or hind...   | {}                       | Alert,Breed,Carry objects,   |
| 3        | B    | Breed              | State event | Otters are engaged in copulatory activities (H...   | reproduction             | Alert,Allogroom,Carry obj... |
| 4        | C    | Carry objects      | State event | Otter carries objects or food by holding them...    | {'0': {'name': "", ...}} | Alert,Allogroom,Breed,Ch...  |
| 5        | C    | Chase              | State event | Otter chases other animals                          | {'0': {'name': "", ...}} | Alert,Allogroom,Breed,Ca...  |
| 6        | D    | Defecate           | State event | Otter eliminates fecal matter (Ethol, 2015)         | physiology               | Alert,Allogroom,Breed,Ca...  |
| 7        | D    | Dig                | State event | Otter uses front legs to move sand, stones on...    | {}                       | Alert,Allogroom,Breed,Ca...  |
| 8        | D    | Drink              | State event | Otter consumes water found in its environmen...     | physiology               | Alert,Allogroom,Breed,Ca...  |
| 9        | E    | Eat                | State event | Otter consumes food                                 | physiology               | {'0': {'name': "", ...}}     |
| 10       | I    | Interact with e... | State event | Otter interacts with enrichment                     | {'0': {'name': "", ...}} | Alert,Allogroom,Breed,Ca...  |
| 11       | L    | Locomotion         | State event | Otter moves from place to place                     | {'0': {'name': "", ...}} | Alert,Allogroom,Breed,Ca...  |
| 12       | F    | Look for food      | State event | Otter looks for food in the environment             | {}                       | Alert,Allogroom,Breed,Ca...  |

The **Ethogram** widget provide the user with the list of behaviors defined in the **Ethogram**. It can be used to record an event by double clicking on the corresponding row. The **Key** column indicates the keyboard key assigned to each behavior (if any). Pressing a key will record the corresponding behavior (that will appear in the *Events* widget).

The behaviors shown in the ethogram widget can be filtered:

**Right-click** on ethogram widget > **Filter behaviors**

Check/Uncheck single behaviors or double-click on the behavioral category

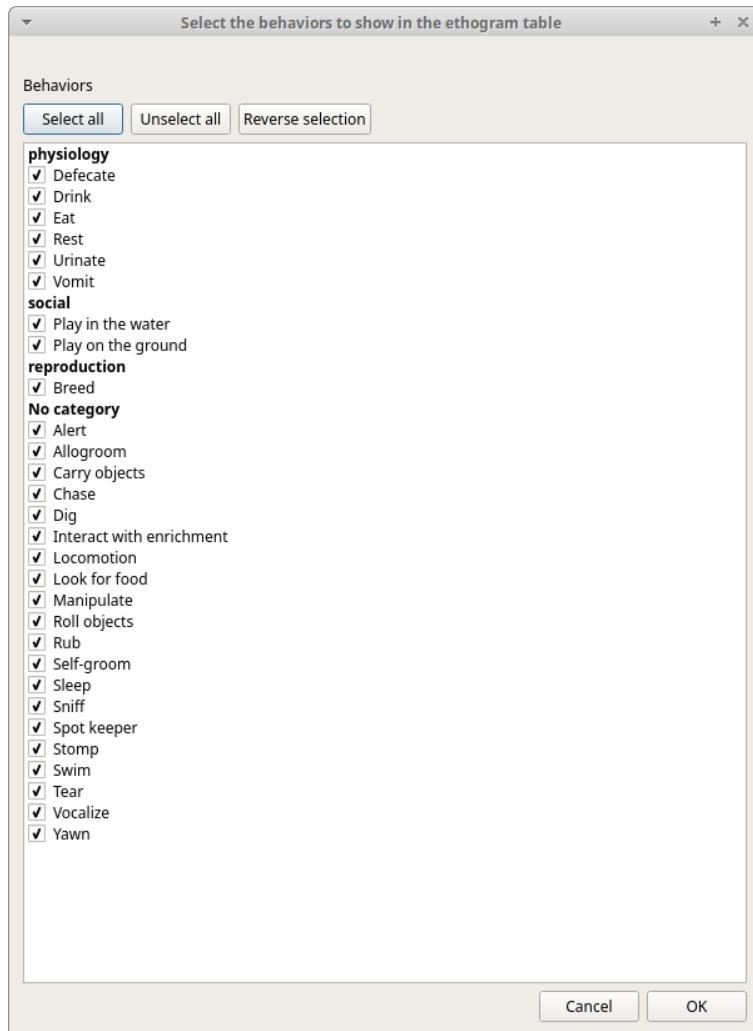

## Subjects widget in the main window


| Subjects |                  |   |                  |
|----------|------------------|---|------------------|
| Key      | Name             | Description                             | Current state(s) |
| 1        | No focal subject |   |                  |
| 2 N      | Nina             | Female, adult, born on 10/03/2013 i...  |                  |
| 3 H      | Himal            | Male, adult, born on 04/30/2014 in ...  |                  |
| 4 C      | Sharky           | Male, juvenile, born on 10/30/2015 i... |                  |
| 5 S      | Nautilus         | Male, juvenile, born on 10/30/2015 i... |                  |

The **Subjects** widget provide the user with the list of subjects defined in the **Subject** tab in the **Project** window. It can be used to add information about the focal subject on the recorded behaviors by double clicking on the corresponding row. When a subject is selected his/her name appears above the media player. The **Key** column indicates the keyboard key assigned to each subject (if any).

The subjects shown in the subjects widget can be filtered:

**Right-click on subjects widget > Filter subjects**


Check/Uncheck single subject

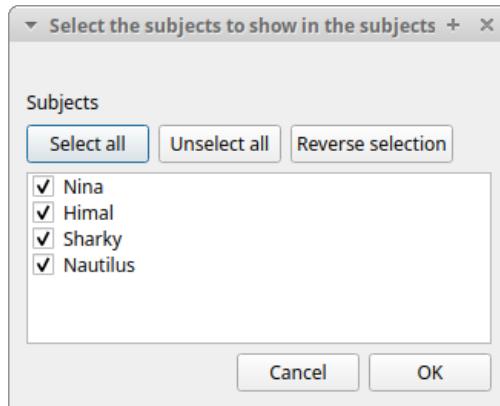

## The media player widgets

The media position can be set using the horizontal slide bar. Each media player has its own audio volume control (the vertical slide bar on the right side of the player).

The organization of the various widget can be customized:

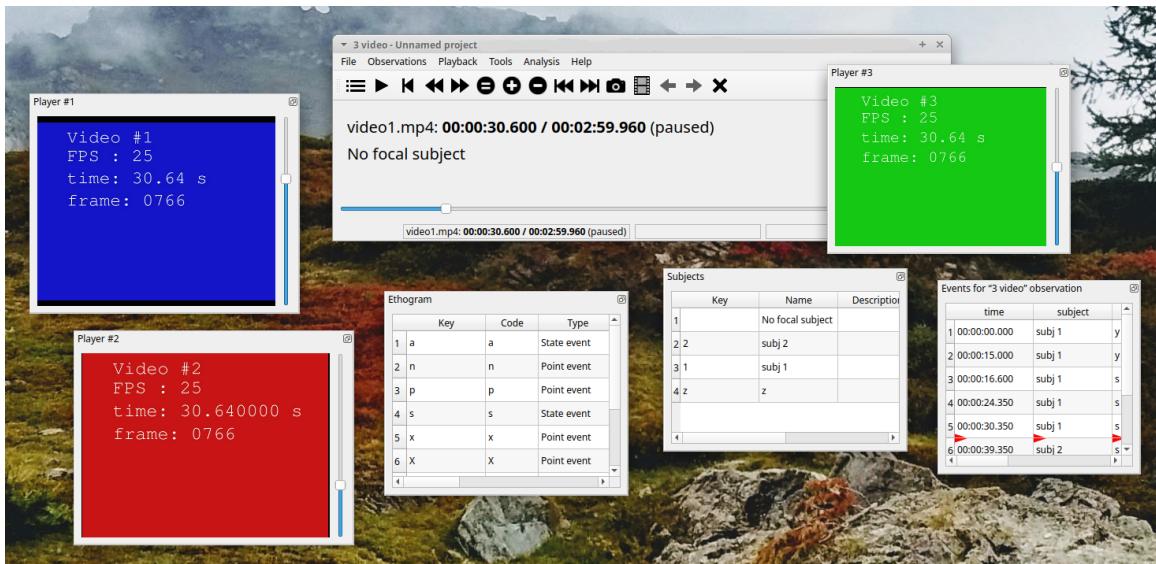

## The Events widget

| Events for "0001_a" observation |          |            |       |          |         |
|---------------------------------|----------|------------|-------|----------|---------|
| time                            | subject  | code       | type  | modifier | comment |
| 1 0.000                         | Himal    | Tear       | START | Branches |         |
| 2 0.000                         | Nautilus | Tear       | START | Branches |         |
| 3 30.199                        | Himal    | Tear       | STOP  | Branches |         |
| 4 30.200                        | Himal    | Locomotion | START | Walk     |         |
| 5 32.400                        | Himal    | Locomotion | STOP  | Walk     |         |
| 6 32.400                        | Nautilus | Tear       | STOP  | Branches |         |

The **Events** widget shows all the recorded behaviors (events) with the following parameters (organized in columns):


- **time**, the time at which the event occurred;
- **subject**, the focal subject (if any);
- **code**, the behavior code;
- **type**, in case of a state event indicates whether the time corresponds to the start or to the stop.
- **modifier**, indicates the modifier(s) that was(ere) selected (if any);
- **comment**, is an open field where the user can add notes.

A tracking cursor (red triangle) will show the current event. This cursor can be positioned above the current event, see [tracking cursor position](#) option in Preferences window.


A double-click on a row will reposition the media player to the moment of the corresponding event. See [Time offset for media reposition](#) in Preferences window to customize the time offset for media repositioning.

## Record an event


Once ready to begin your coding, you can start the media player using the **Play** button (or the **Space bar**). The behaviors can be recorded using the keyboard with the predefined keys, by double-clicking the corresponding row in the **Ethogram** table or by using the **Coding pad** (See [Coding pad](#)).

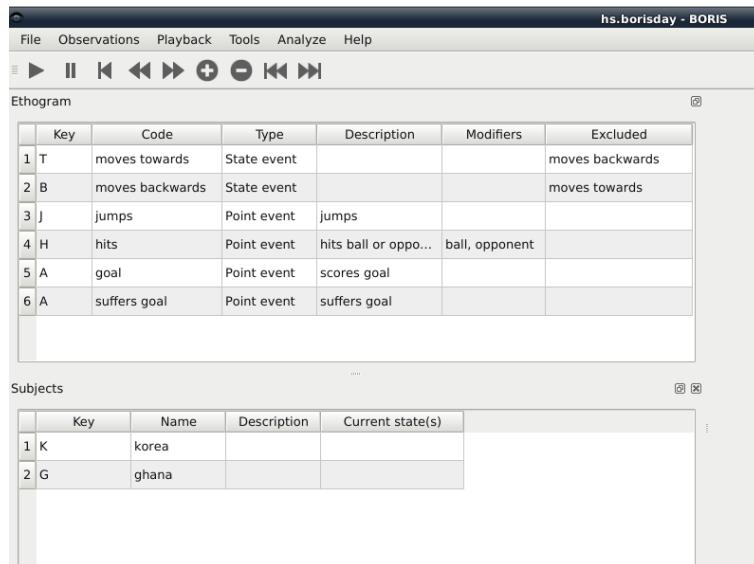

If the pressed key defines a single event, the corresponding event will be recorded in the **Events** table. In the case you have specified the same key for two (or more) events (e.g. key A in the figure below), BORIS will prompt you for the desired behavior.

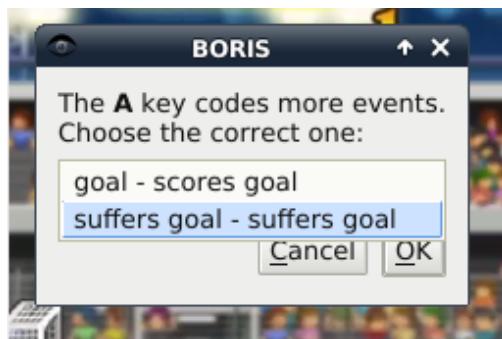

In the case you have specified modifiers (one or more sets), BORIS will prompt you for the desired modifier(s) if any (e.g. **ball** or **opponent** in the figure below). You can select the correct one using the mouse or the keyboard ( **b** key or **o** key)

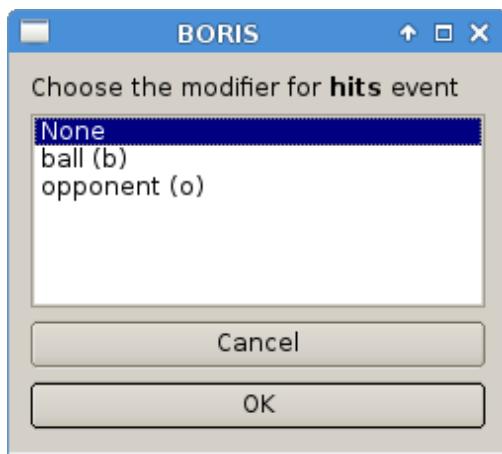

If no keys are defined for the modifier selection, you can type the first character of the modifier and use the **Up arrow** and **Down arrow** keyboard keys to select the correct modifier.


In the case your behavior type is a **Point event with coding map** or a **State event with coding map**, BORIS will show the **Coding map** window and will allow selecting the desired area(s). In case you click a part of the map in which two (or more) areas overlap, the corresponding codes will be recorded.

A recorded event can be edited (once selected) using the **Observations > Edit event** menu option. The resulting *Edit event parameters* allows modifying every parameter (e.g. time, subject, code, modifiers, and comment).

The **Observations > Add event** menu option allows adding a new event by specifying its time and the other parameters.

### **The Events widget context menu**

Some functions are available in the Events widget context menu. Right-click on the Events widget.

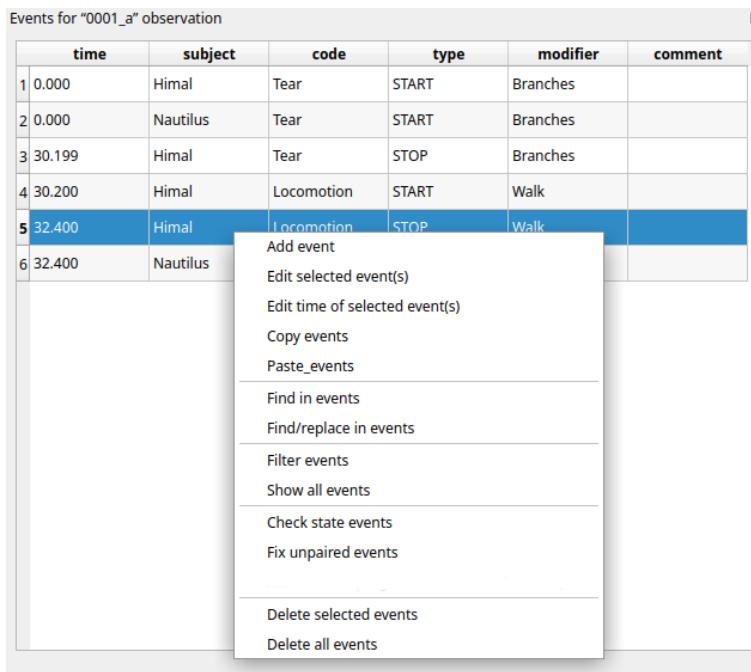

#### **Add event**

This option allows adding a new event by specifying its time and the other parameters.

#### **Edit selected event(s)**

This option allows to edit the selected event(s). When many events are selected you have to choose the field to edit between **Subject**, **Behavior** and **Comment**. In this case the new value will apply to all selected events.

| time     | subject  | code       | type  | modifier | comment |
|----------|----------|------------|-------|----------|---------|
| 1 0.000  | Himal    | Tear       | START | Branches |         |
| 2 0.000  | Nautilus | Tear       | START | Branches |         |
| 3 30.199 | Himal    | Tear       | STOP  | Branches |         |
| 4 30.200 | Himal    | Locomotion | START | Walk     |         |
| 5 32.400 | Himal    | Locomotion | STOP  | Walk     |         |
| 6 32.400 | Nautilus |            |       |          |         |

Edit selected events

Subject  
 Behavior

New value

Nina  
Himal  
Sharky **Nautilus**  
Nautilus

Comment

New comment

### Edit time of selected event(s)

This option allows to add or subtract a time value (in seconds) to all selected events. For subtracting a value use a negative value.

| time     | subject  | code       | type  | modifier | comment |
|----------|----------|------------|-------|----------|---------|
| 1 0.000  | Himal    | Tear       | START | Branches |         |
| 2 0.000  | Nautilus | Tear       | START | Branches |         |
| 3 30.199 | Himal    | Tear       | STOP  | Branches |         |
| 4 30.200 | Himal    | Locomotion | START | Walk     |         |
| 5 32.400 | Himal    | Locomotion | STOP  | Walk     |         |
| 6 32.400 | Nautilus | Tear       | STOP  | Branches |         |

Time value

Value to add or subtract (use negative value):

### Copy events

This option allows to copy the selected events in the clipboard. The clipboard will contain the values of the selected events (except the **type** field) separated by a <TAB> character.

Example of clipboard content:

```
0.0 Himal    Tear    Branches
0.0 Nautilus      Tear    Branches
30.199   Himal    Tear    Branches
30.2     Himal   Locomotion    Walk
32.4     Himal   Locomotion    Walk
32.4     Nautilus      Tear    Branches
```

### Paste events

This option allows to paste the clipboard content into the events widget. The clipboard must respect the format described in the previous section: 5 columns separated by <TAB> character.

### **Find in events**

This option allows to search for a string in the various field of events. Select the fields to be searched. The find/replace operation can be restricted to the selected events.

|   | time   | subject  | code       | type  | modifier | comment |
|---|--------|----------|------------|-------|----------|---------|
| 1 | 0.000  | Himal    | Tear       | START | Branches |         |
| 2 | 0.000  | Nautilus | Tear       | START | Branches |         |
| 3 | 30.199 | Himal    | Tear       | STOP  | Branches |         |
| 4 | 30.200 | Himal    | Locomotion | START |          |         |
| 5 | 32.400 | Himal    | Locomotion | STOP  |          |         |
| 6 | 32.400 | Nautilus | Tear       | STOP  |          |         |

Find in events

Subject  
 Behavior  
 Modifiers  
 Comment

Find

Locom

Find in selected events

Close Find

### **Find/Replace in events**

This option allows to search for a string and replace it by a new value in the various field of events. Select the fields to be searched. The find operation can be restricted to the selected events.

|   | time   | subject  | code       | type  | modifier | comment |
|---|--------|----------|------------|-------|----------|---------|
| 1 | 0.000  | Himal    | Tear       | START | Branches |         |
| 2 | 0.000  | Nautilus | Tear       | START | Branches |         |
| 3 | 30.199 | Himal    | Tear       | STOP  | Branches |         |
| 4 | 30.200 | Himal    | Locomotion | START | Walk     |         |
| 5 | 32.400 | Himal    |            |       |          |         |
| 6 | 32.400 | Nautilus |            |       |          |         |

Find/Replace events

Subject  
 Behavior  
 Modifiers  
 Comment

Find

Himal

Replace


Coco

Find/Replace in selected events

Cancel Find and replace Find and replace all

### **Filter events**

This option allows to filter the events by field value.

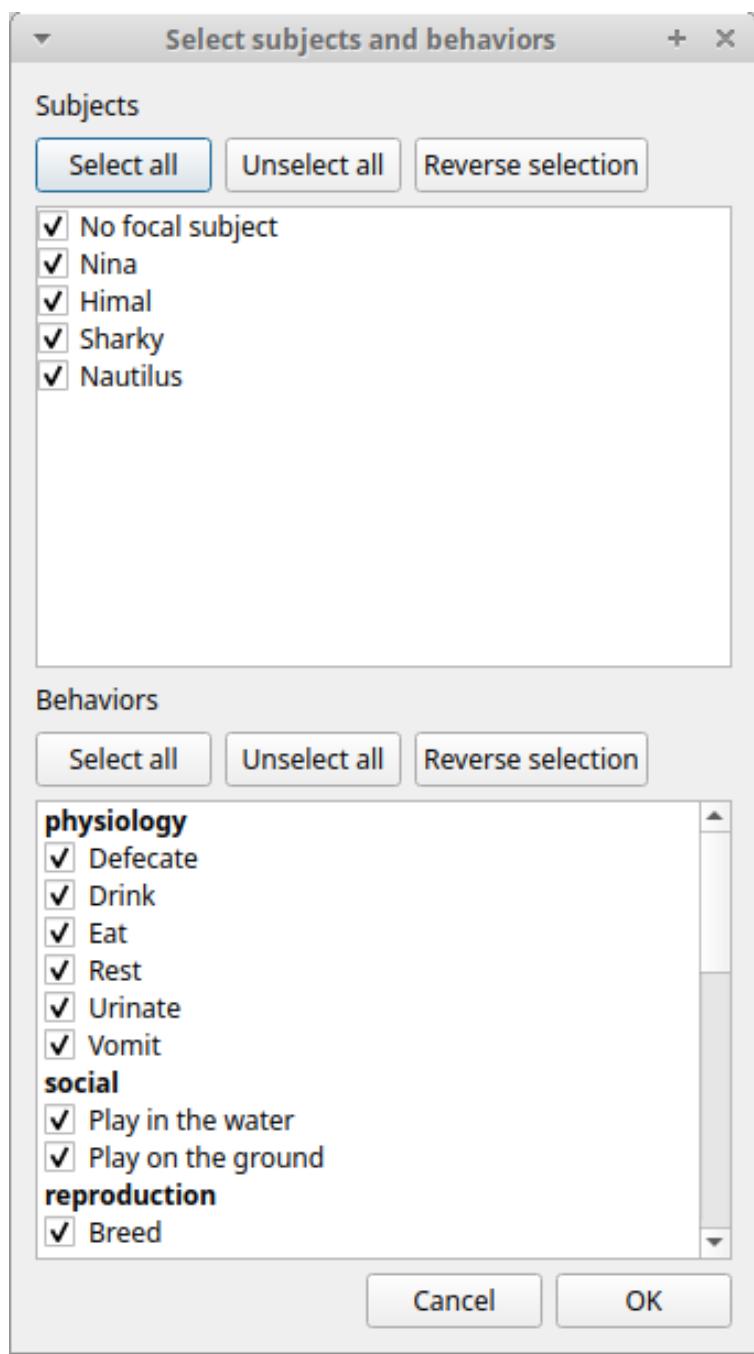

#### **Show all events**

This option reverts the previous one and allows to visualize all coded events

#### **Check state events**

This option allows to check if the **state events** are **PAIRED**, if they have a **START** and a **STOP** occurrences.

#### **Delete selected events**

This option allows to delete the selected events. This operation is irreversible!

#### **Delete all events**

This option allows to delete all then events in the current observation. This operation is irreversible!

## Fix unpaired state events

You can use the "Fix unpaired events" function to fix the **state events** without a STOP event.

**Observations > Fix unpaired events** (keyboard shortcut: CTRL+U)

The program will ask for a time at which insert the STOP events for all unpaired **state events**

This function can be run on a set of selected observations (when no observation is open). In this case the STOP events will be inserted at the end of observation.

## Explore project

You can search information in various fields in all observations in the current project (Observations > Explore project).

The searchable fields are: **subject**, **behavior**, **modifier** and **comment**.

If more than one field is searched a logic AND will apply.



The events that were found are listed in a table. By double-clicking on the row the corresponding observation will be opened and the visualization will be scrolled to the row corresponding to the event.

The screenshot shows the BORIS software interface. On the left, there's an 'Ethogram' table with rows numbered 1 to 15. The 'Subjects' table below it lists three entries: 'No focal subject', 'Nina', and 'Himal'. The main workspace shows a 'Events for "0004" observation' table with columns: time, subject, code, and type. The table contains numerous rows of data, with the 32nd row (Himal, Alert) highlighted in blue. A small 'BORIS' window is also visible in the center of the workspace.

## Frame-by-frame mode

You can switch between the media player and the frame-by-frame mode using the dedicated button in the toolbar:



In frame-by-frame mode the video will stop playing and the user will visualize the video frame by frame.

## **Warning**

Please note that MTS video files should be re-encoded to be used in frame-by-frame mode. Otherwise the extracted frames are not reliable.

You can move between frames by using the arrow keys in the toolbar (on the right) or by using keyboard special keys:

### **For the both modes (media player and frame-by-frame mode)**

- **Page Up** key: switch to the next media
- **Page Down** key: switch to the previous media
- **Up arrow** key: jump forward in the current media
- **Down arrow** key: jump backward in the current media
- **Esc**: switch between VLC and frame-by-frame mode

### **Only for the frame-by-frame mode**

- **Left arrow** key: go to the previous frame
- **Right arrow** key: go to the next frame

If you have a numeric keypad you can use the following keys in alternative:

- The key **/** will allow you to view the previous frame
- The key **\*** will allow you to view the next frame

To return in the media player mode press again the frame-by-frame button in the toolbar.

The frame can be resized before visualization. See the [frame resizing](#) option.

The frame viewer can be detached from the main window (See [File > Preferences](#)). In case on 2 simultaneous players the frame viewers will be automatically detached from the main window.

For every second of the media file the frames are extracted by BORIS with the embedded ffmpeg program and are saved in the [FFmpeg cache directory](#) specified in the [Preferences](#) window.

## **Exporting events data**

The coded events can be exported in various formats:

### **Export events**

#### **Observations > Export events**

This function will export the events of selected observations in TSV, CSV, ODS, XLSX, XLS or HTML formats. If many observations are selected BORIS will ask for a directory to save the various files. For the XLS and XLSX formats the events can be exported on various worksheet in a single workbook. These formats are suitable for further analysis.

| A                        | B                       | C    | D        | E       | F      |
|--------------------------|-------------------------|------|----------|---------|--------|
| 1 Observation id         | demo #1                 |      |          |         |        |
| 2                        |                         |      |          |         |        |
| 3 Media file(s)          |                         |      |          |         |        |
| 4                        |                         |      |          |         |        |
| 5 Player #1              | /home/olivier/crop.avi  |      |          |         |        |
| 6 Player #1              | /home/olivier/crop2.avi |      |          |         |        |
| 7                        |                         |      |          |         |        |
| 8 Observation date       | 2015-Nov-30 10:39:18    |      |          |         |        |
| 9                        |                         |      |          |         |        |
| 10 Description           |                         |      |          |         |        |
| 11                       |                         |      |          |         |        |
| 12 Time offset (s)       | 0                       |      |          |         |        |
| 13                       |                         |      |          |         |        |
| 14 Independent variables |                         |      |          |         |        |
| 15 variable              | value                   |      |          |         |        |
| 16 group size            | 7                       |      |          |         |        |
| 17 weather conditions    | rain                    |      |          |         |        |
| 18                       |                         |      |          |         |        |
| 19 time                  | subject                 | code | Modifier | comment | status |
| 20                       | 4.3 Subject #1          | eat  | salad    | test    | START  |
| 21                       | 10 Subject #1           | eat  | salad    |         | STOP   |
| 22                       | 26.6 Subject #2         | eat  | meat     |         | START  |
| 23                       | 113.988 Subject #2      | eat  | meat     |         | STOP   |
| 24                       | 116.588 Subject #1      | jump |          |         | POINT  |
| 25                       | 118.988 Subject #1      | jump |          |         | POINT  |
| 26                       | 120.863 Subject #2      | jump |          |         | POINT  |
| 27                       | 122.438 Subject #2      | jump |          |         | POINT  |
| 28                       |                         |      |          |         |        |

## Warning

Please note that for some formats (XLS - Excel 97) the name of the sheet will be based on a modified **observation id** in order to not contain forbidden characters (: / ? \* [ or ] ) and shortened to 31 characters.

## Export aggregated events

### Observations > Export aggregated events

This function will export the events of the selected observations in the following formats:

- **tabular format** (TSV, CSV, XLSX, XLS, ODS, HTML)
- **SQL** format for populating a SQL database
- **SDIS** format for analysis with the GSEQ program available at <http://www2.gsu.edu/~psyrab/gseq>

If many observations are selected you can choose to group all results in one file. If you do not want to group results BORIS will ask for a directory to save the various files.

The **State events** are paired and the duration is available.

An arbitrary time interval can be selected (check the **Limit to time interval** option). In this case the ongoing events will be started at start time and stopped at end time in the export file.

Example of tabular export

| A                | B                    | C          | D        | E         | F             | G       | H       | I             | J            |
|------------------|----------------------|------------|----------|-----------|---------------|---------|---------|---------------|--------------|
| 1 Observation id | Observation date     | Subject    | Behavior | Modifiers | Behavior type | Start   | Stop    | Comment start | Comment stop |
| 2 demo #1        | 2015-Nov-30 10:39:18 | Subject #1 | jump     |           | POINT         | 116.588 | 0       |               |              |
| 3 demo #1        | 2015-Nov-30 10:39:18 | Subject #1 | jump     |           | POINT         | 118.988 | 0       |               |              |
| 4 demo #1        | 2015-Nov-30 10:39:18 | Subject #1 | eat      | salad     | STATE         | 4.3     | 10      | test          |              |
| 5 demo #1        | 2015-Nov-30 10:39:18 | Subject #2 | jump     |           | POINT         | 120.863 | 0       |               |              |
| 6 demo #1        | 2015-Nov-30 10:39:18 | Subject #2 | jump     |           | POINT         | 122.438 | 0       |               |              |
| 7 demo #1        | 2015-Nov-30 10:39:18 | Subject #2 | eat      | meat      | STATE         | 26.6    | 113.988 |               |              |
| 8                |                      |            |          |           |               |         |         |               |              |
| 9                |                      |            |          |           |               |         |         |               |              |
| 10               |                      |            |          |           |               |         |         |               |              |

Example of SQL export:

```
CREATE TABLE events (id INTEGER PRIMARY KEY ASC, observation TEXT,
date DATE, subject TEXT, behavior TEXT,
modifiers TEXT, event_type TEXT, start FLOAT,
stop FLOAT, comment_start TEXT,
comment_stop TEXT);
```

```

INSERT INTO events (observation, date, subject, behavior, modifiers,
event_type, start, stop, comment_start, comment_stop ) VALUES
("demo #1", "2015-11-30 10:39:18", "Subj #1", "jump", "", "POINT", 116.588, 0, "", ""),
("demo #1", "2015-11-30 10:39:18", "Subj #1", "jump", "", "POINT", 118.988, 0, "", ""),
("demo #1", "2015-11-30 10:39:18", "Subj #1", "eat", "salad", "STATE", 4.3, 10.0, "vvv", ""),
("demo #1", "2015-11-30 10:39:18", "Subj #2", "jump", "", "POINT", 120.863, 0, "", ""),
("demo #1", "2015-11-30 10:39:18", "Subj #2", "jump", "", "POINT", 122.438, 0, "", ""),
("demo #1", "2015-11-30 10:39:18", "Subj #2", "eat", "meat", "STATE", 26.6, 113.988, "", );

```

### **Warning**

Please note that for some formats (XLS - Excel 97) the name of the sheet will be based a modified **observation id** in order to not contain forbidden characters (: / ? \* [ or ]) and shortened to 31 characters.

## **Export events for analysis with JWatcher**

JWatcher is a powerful tool for the quantitative analysis of behavior.

The events coded with BORIS can be exported to be analyzed with JWatcher.

Click **Observations > Export events > for analysis with JWatcher** to export the coded events.

BORIS will ask for selecting a directory. After this, for each combination of selected observation and selected subject the following files will be created:

- the Focal Data File (.dat)
- the Focal Analysis Master File (.faf)
- the Focal Master File (.fmf)

These files can be used to analyze your observations with JWatcher.

## **Export events as behavioral strings**

Behavioral strings can be used with the **Behatrix** program: [Behatrix \(formerly BSA\)](#)

Example:

```

# observation id: demo#1
# observation description:
# Media file name: video1.mp4, video2.mp4

Subject #1:
eat|jump|eat|jump

Subject #2:
eat|rest|jump|eat|jump

```

## **Export events as Praat TextGrid**

Example:

```

File type = "ooTextFile"
Object class = "TextGrid"

xmin = 4.3
xmax = 113.988
tiers? <exists>
size = 2
item []:
    item [1]:
        class = "IntervalTier"
        name = "Subject #1"
        xmin = 4.3
        xmax = 10.0
        intervals: size = 1
        intervals [1]:
            xmin = 4.3
            xmax = 10.0
            text = "eat"
    item [2]:
        class = "IntervalTier"
        name = "Subject #2"
        xmin = 26.6
        xmax = 113.988
        intervals: size = 1
        intervals [1]:
            xmin = 26.6
            xmax = 113.988
            text = "eat"

```

## Export events as Behavioral Binary Table

A time interval will be asked to the user (in seconds). The observation will be checked every n seconds and the presence of the selected behaviors will be exported in a table for each selected subjects.

Example for a time interval of 1 second:

| time | Alert | Drink | Locomotion | Swim |
|------|-------|-------|------------|------|
| 0.0  | 0     | 1     | 0          | 0    |
| 1.0  | 0     | 1     | 0          | 0    |
| 2.0  | 0     | 1     | 0          | 0    |
| 3.0  | 0     | 1     | 0          | 0    |
| 4.0  | 0     | 1     | 0          | 0    |
| 5.0  | 0     | 1     | 0          | 0    |
| 6.0  | 0     | 1     | 0          | 0    |
| 7.0  | 0     | 1     | 0          | 0    |
| 8.0  | 0     | 1     | 0          | 0    |
| 9.0  | 1     | 0     | 0          | 0    |
| 10.0 | 1     | 0     | 0          | 0    |
| 11.0 | 0     | 0     | 1          | 0    |
| 12.0 | 1     | 0     | 0          | 0    |
| 13.0 | 1     | 0     | 0          | 0    |
| 14.0 | 1     | 0     | 0          | 0    |
| 15.0 | 1     | 0     | 0          | 0    |
| 16.0 | 1     | 0     | 0          | 0    |
| 17.0 | 1     | 0     | 0          | 0    |

|      |   |   |   |   |
|------|---|---|---|---|
| 18.0 | 0 | 0 | 1 | 0 |
| 19.0 | 0 | 0 | 1 | 0 |
| 20.0 | 0 | 0 | 1 | 0 |
| 21.0 | 0 | 0 | 1 | 0 |
| 22.0 | 1 | 0 | 0 | 0 |
| 23.0 | 0 | 0 | 0 | 0 |
| 24.0 | 0 | 0 | 0 | 0 |
| 25.0 | 0 | 0 | 0 | 0 |
| 26.0 | 0 | 0 | 0 | 0 |
| 27.0 | 0 | 0 | 0 | 0 |
| 28.0 | 0 | 0 | 0 | 0 |
| 29.0 | 0 | 0 | 0 | 0 |
| 30.0 | 0 | 0 | 0 | 0 |
| 31.0 | 0 | 0 | 0 | 0 |
| 32.0 | 0 | 0 | 0 | 0 |
| 33.0 | 0 | 0 | 0 | 1 |
| 34.0 | 0 | 0 | 0 | 1 |
| 35.0 | 0 | 0 | 0 | 1 |
| 36.0 | 0 | 0 | 0 | 1 |
| 37.0 | 0 | 0 | 0 | 1 |
| 38.0 | 0 | 0 | 0 | 1 |
| 39.0 | 0 | 0 | 0 | 1 |
| 40.0 | 0 | 0 | 0 | 1 |
| 41.0 | 0 | 0 | 0 | 1 |
| 42.0 | 0 | 0 | 0 | 1 |
| 43.0 | 0 | 0 | 0 | 1 |
| 44.0 | 0 | 0 | 0 | 1 |
| 45.0 | 0 | 0 | 0 | 1 |
| 46.0 | 0 | 0 | 0 | 1 |
| 47.0 | 0 | 0 | 0 | 1 |
| 48.0 | 0 | 0 | 0 | 1 |
| 49.0 | 0 | 0 | 0 | 1 |
| 50.0 | 0 | 0 | 0 | 1 |
| 51.0 | 0 | 0 | 0 | 1 |
| 52.0 | 0 | 0 | 0 | 1 |
| 53.0 | 0 | 0 | 0 | 0 |
| 54.0 | 0 | 0 | 0 | 0 |
| 55.0 | 0 | 0 | 0 | 0 |
| 56.0 | 0 | 0 | 0 | 0 |
| 57.0 | 0 | 0 | 0 | 1 |
| 58.0 | 0 | 0 | 0 | 1 |
| 59.0 | 0 | 0 | 0 | 1 |
| 60.0 | 0 | 0 | 0 | 1 |
| 61.0 | 0 | 0 | 0 | 1 |
| 62.0 | 0 | 0 | 0 | 1 |
| 63.0 | 0 | 0 | 0 | 1 |

## Extract sequences from media files corresponding to coded events

Sequences of media file corresponding to coded events can be extracted from media files:

1. Click on **Observations > Extract events from media files** option.
2. Choose the observation(s).
3. Select the events to be extracted.

4. Select a destination directory that will contain the extracted sequences.

5. Select a time offset (in seconds, the default value is 0).

The time offset will be subtracted from the starting time of event and added to the stopping time. All the extracted sequences will be saved in the selected directory followind the file name format:

{observation id}\_{player}\_{subject}\_{behavior}\_{start time}-{stop time}

## Extract frames corresponding to coded events

The frames corresponding to coded events can be extracted and saved as images.

1. Click on **Observations > Extract frames from media files** option.

2. Choose the observation(s).

3. Select the events to be extracted.

4. Select a destination directory that will contain the extracted sequences.

5. Select a time offset (in seconds, the default value is 0).

## Export transitions matrix

3 transitions matrix outputs are available: The matrix of frequencies of transitions, the matrix of frequencies of transition after each behavior and the matrix of number of transitions.

### ***Matrix of frequencies of transitions***

This matrix contains the frequencies of total transitions. The sum of all frequencies must be 1.

Example of frequencies of transitions matrix:

|       | eat   | sleep | walk  |
|-------|-------|-------|-------|
| eat   | 0.0   | 0.286 | 0.143 |
| sleep | 0.143 | 0.0   | 0.143 |
| walk  | 0.286 | 0.0   | 0.0   |

In this matrix you can see that the **eat** behavior precedes the **sleep** behavior with a frequency of **0.286** of the total number of transitions.

### ***Matrix of frequencies of transitions after behavior***

This matrix contains the frequencies of transitions after each behavior. The sum of each row must be 1.

Example:

|       | eat | sleep | walk  |
|-------|-----|-------|-------|
| eat   | 0.0 | 0.667 | 0.333 |
| sleep | 0.5 | 0.0   | 0.5   |
| walk  | 1.0 | 0.0   | 0.0   |

In this example you can see that **sleep** follows **eat** with a frequency of **0.667** and **walk** follows with a frequency of **0.333**.

### ***Matrix of number of transitions***

This matrix contains the number of transitions after each behavior.

Example:

|       | eat | sleep | walk |
|-------|-----|-------|------|
| eat   | 0   | 2     | 1    |
| sleep | 1   | 0     | 1    |
| walk  | 2   | 0     | 0    |

## Playback menu

### Jump

#### **Jump forward**

Allow to jump forward in the current media file. See **File > Preferences** for setting the jump value.

#### **Jump backward**

Allow to jump backward in the current media file. See **File > Preferences** for setting the jump value.

#### **Jump to specific time**

Allow to go to a specific time in the current media file.

### Zoom

Allow to zoom into the current video file. The available zoom values are: 1:2, 1:4, 2:1 and 1:1. **Fit to window** adapts the zoom value to the current window size. In case of simultaneous playing you can set different values of zoom in player #1 and player #2.

The Zoom option is not available in frame-by-frame mode.

#### **Warning**

**The Zoom function is not working on Mac**

## Tools


### Focus a video area

It is possible to focus an area of the video. Click on the center of the area you want to focus. The video will be displayed at its native resolution (not shrinked) and centered on the clicked point. Another click on the video will restore the previous visualization.


The focused area can also be zoomed (only on the first player for now).

#### **Note**

This function is not yet available in frame-by-frame mode.

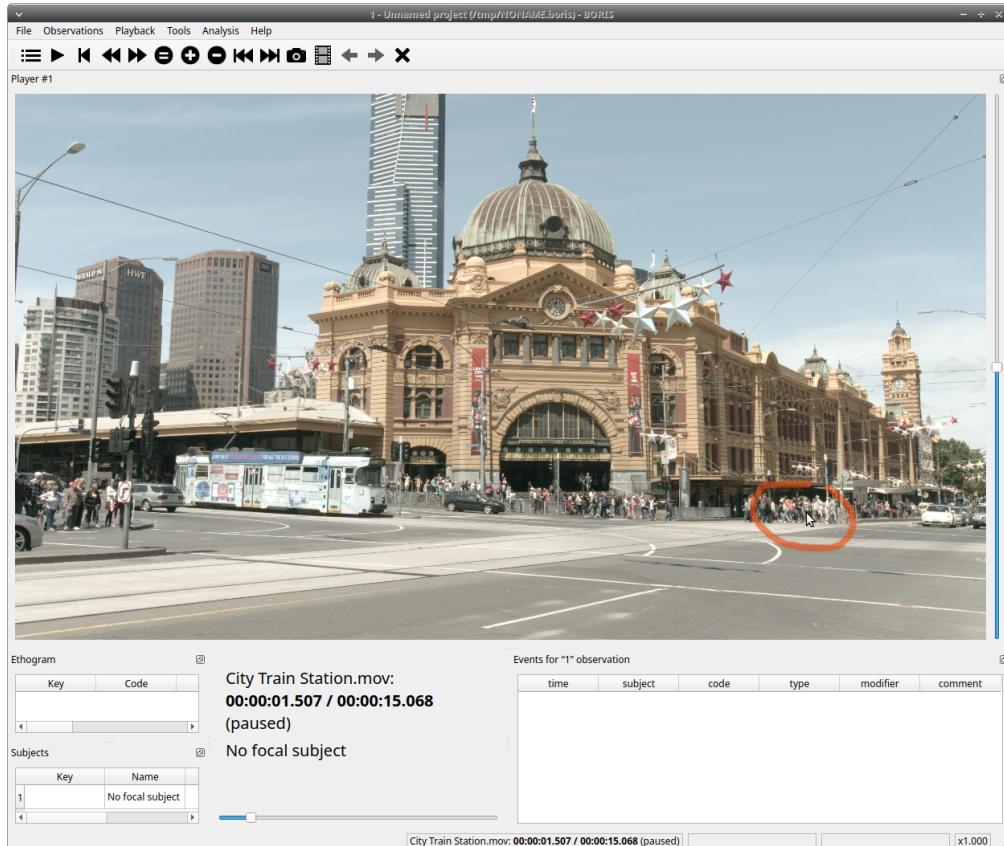


The visualization can be focused on a particular region of the video by clicking on the video (with the left button). See the mouse cursor on the screenshot above.

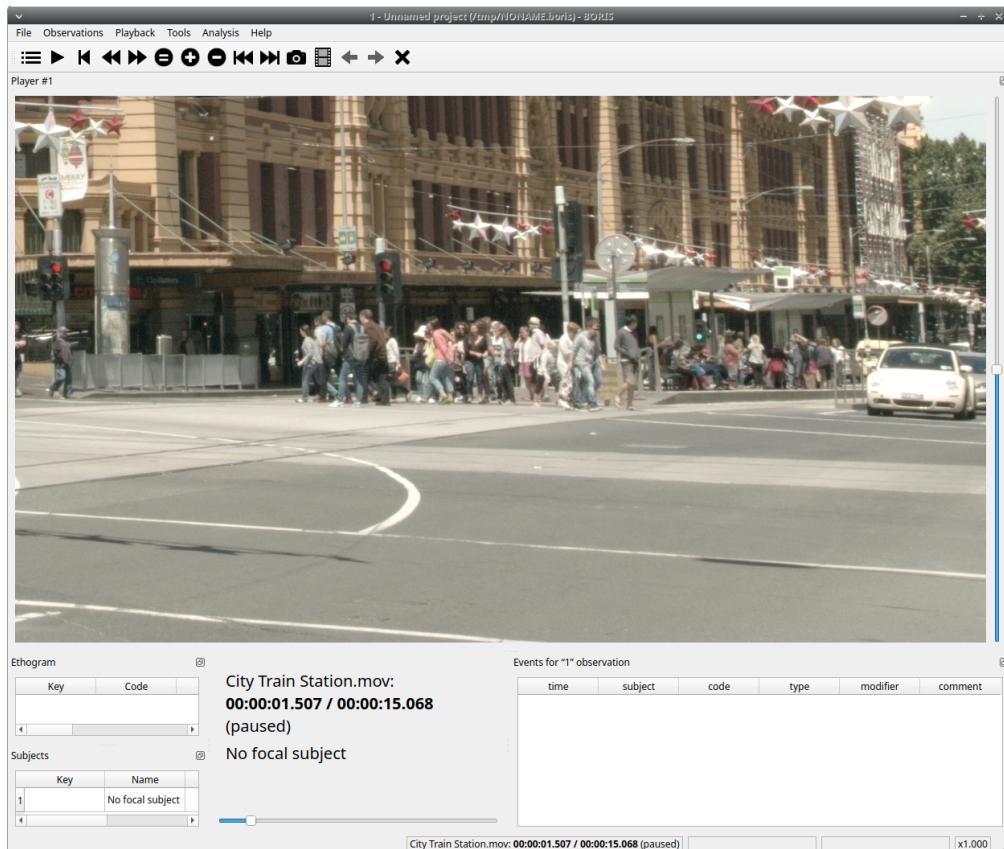

Now the video is displayed with its real resolution allowing to see details. Another click on the video will restore the previous visualization without focus.

## Geometric measurements

Geometric measurements can only be made in frame-by-frame mode. Distances, areas and angles can be measured. Click on **Tools > Geometric measurements** to activate the measurements. The **Measurements window** will be shown:

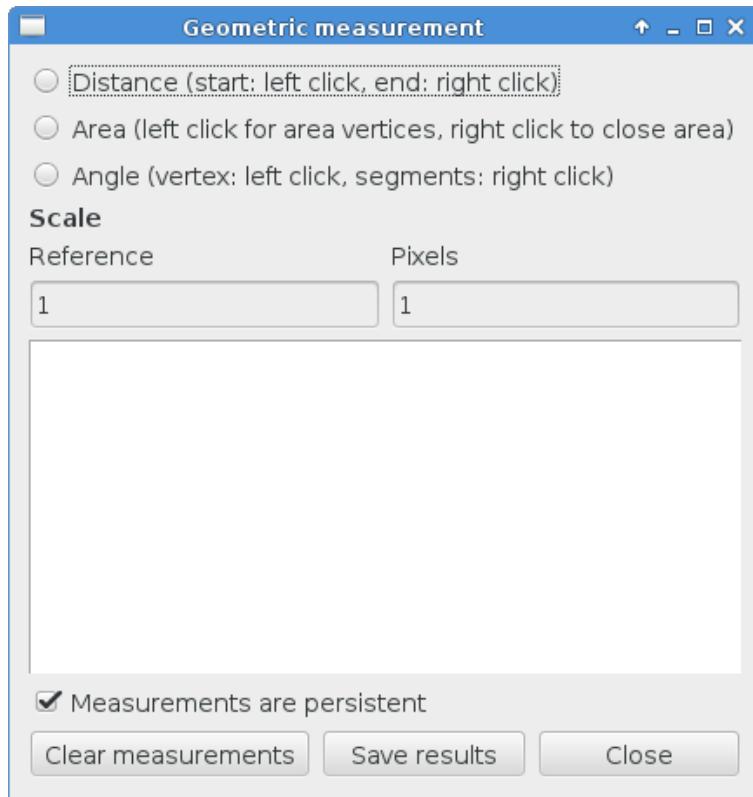


### Setting the scale

For distance and area measurements you can set a scale in order to have results of measurements in a real unit (like centimeters, meters etc).

1. measure a reference object (that have a known size) on the frame (with the distance tool. See next chapter for details) and set the pixel distance in the **Pixel** text box.
2. Set the real size of the reference object in the **Reference** text box (must be a number without unit).

### Distance measurements

Select the **Distance** radio button. Click the left mouse button on the frame bitmap to set the start of the segment that will be measured. A blue circle with a cross will be drawn. Click the right mouse button to set the end. A red circle with a cross will be drawn. The distance between the two selected points will be available in the text area of the **Measurements window**.

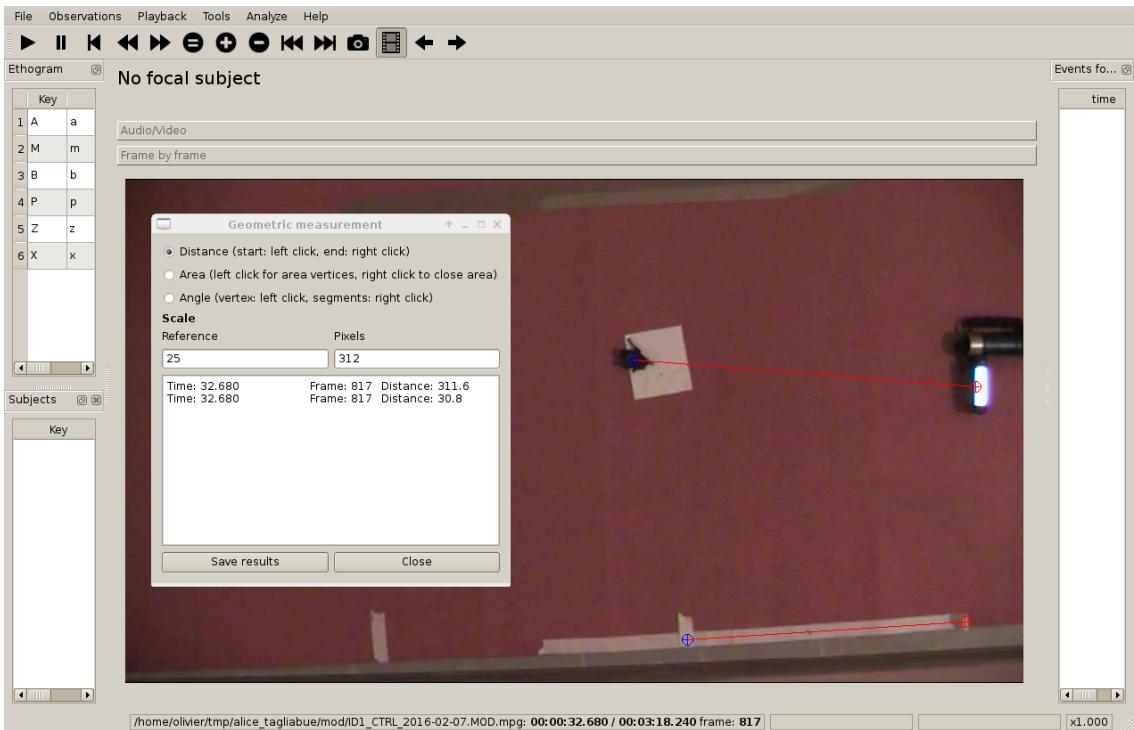

## Area measurements

Select the **Area** radio button. Click the left mouse button on the frame bitmap to set the area vertices. Blue circles with a cross will be drawn. Click the right mouse button to close the area. The area of the drawn polygon will be available in the text area of the **Measurements window**.

## Angle measurements


Select the **Angle** radio button. Click the left mouse button on the frame bitmap to set the angle vertex. A red circle with a cross will be drawn. Click the right mouse button to set the two segments. Blue circles with a cross will be drawn. The angle between the two drawn segments will be available in the text area of the **Measurements window**.

## Persistent measurements

If the **Measurements are persistent** checkbox is checked the measurement schemes will be available on all frames otherwise they will be deleted between frames.

## Coding pad

During observation a coding pad with the available behaviors can be displayed (**Tools > Coding pad**). This **Coding pad** allows the user to code using a touch-screen or by clicking on the buttons. When the **Coding pad** is displayed you can continue to code using the keyboard or the ethogram.

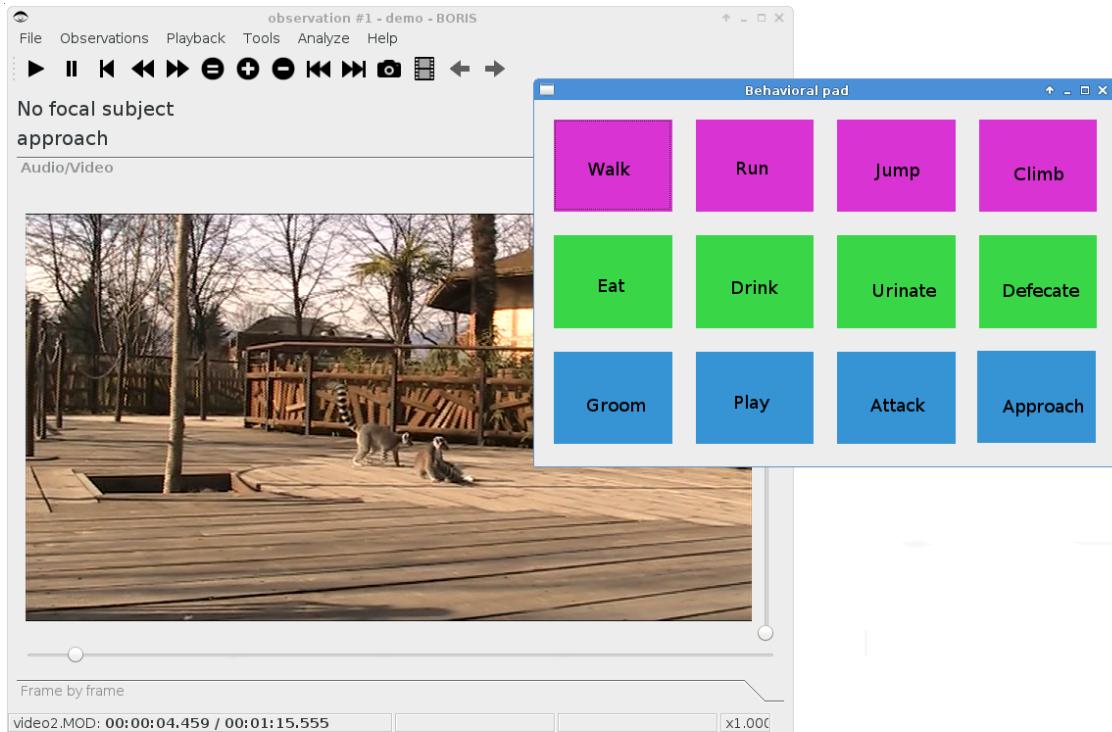


## Subjects pad

A pad with all defined subjects (or filtered subjects) can be displayed during the observation (**Tools > Subjects pad**). This **Subjects pad** allows the user to select the focal subject using a touch-screen or by clicking on the buttons. When the **Subjects pad** is displayed you can continue to select the focal subject using the keyboard or the subjects list.



## Image overlay on video

An image can be overlaid on the video (**Tools > Image overlay on video > Add an image overlay**).

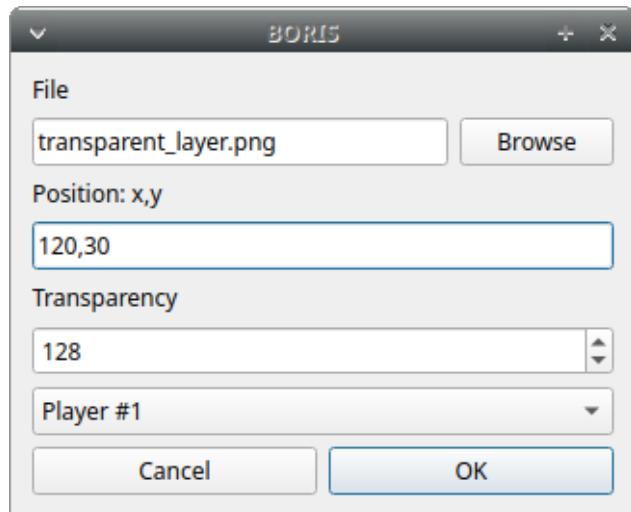

The image must be in PNG format, if the image is smaller than the video resolution the image position can be set from the top-left corner. The transparency of image can be set from 0 (full transparency) to 255 (no transparency).

The remove the image overlay click **Tools > Image overlay on video > Remove**

## Converters for external data values

Converters can be written using the Python 3 programming language.

The **INPUT** variable will be loaded with the original value of the external data file (for example 01:22:32).

The **OUTPUT** variable must contain the converted value in seconds (the dot must be used for decimal separator).

Example of a converter for converting HH.MM:SS format in seconds:

```
h, m, s = INPUT.split(':')
OUTPUT = int(h) * 3600 + int(m) * 60 + int(s)
```


The Python function **strptime()** from the **datetime** module can be useful for converting time values:  
<https://docs.python.org/3/library/datetime.html#strftime-strptime-behavior>

Example of a converter for converting ISO8601 format in seconds using the strptime() function:

```
import datetime
epoch = datetime.datetime.utcfromtimestamp(0)
datetime_format = "%Y-%m-%dT%H:%M:%S%z"

OUTPUT = (datetime.datetime.strptime(INPUT, datetime_format) - epoch).total_seconds()
```

**File > Edit project > Converters tab**

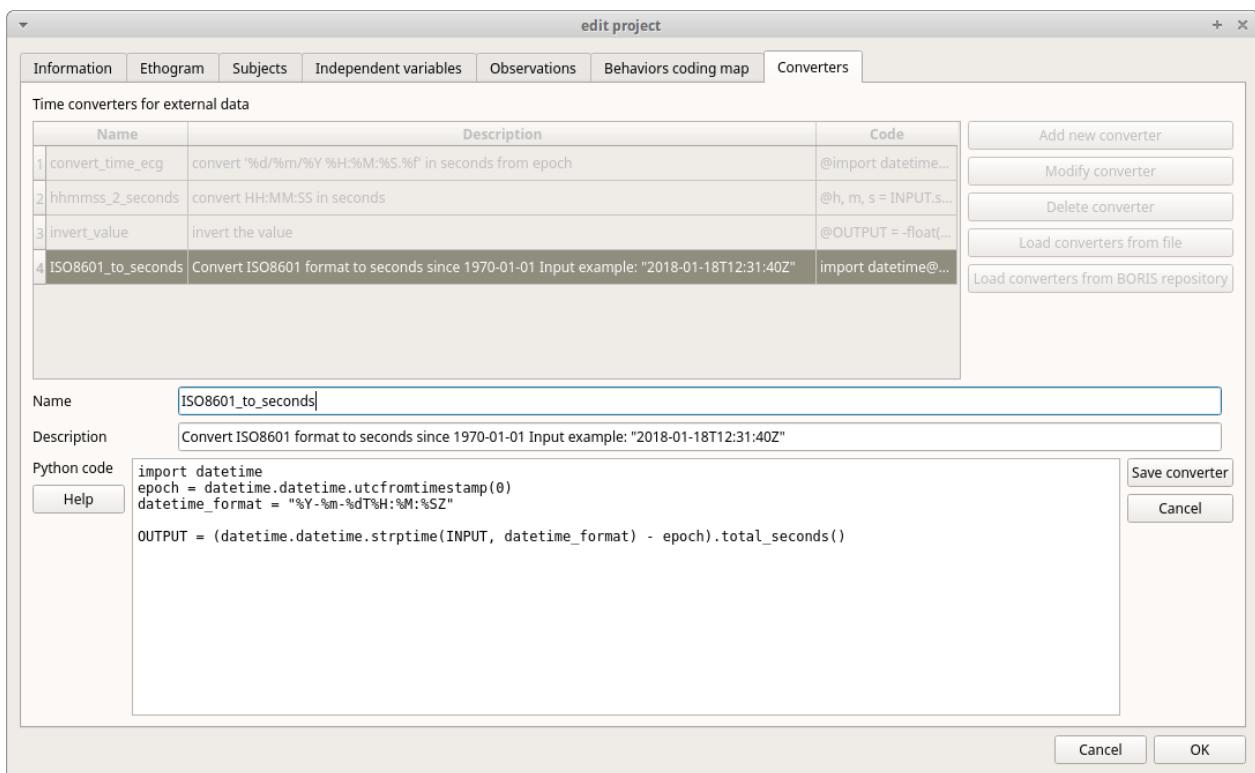

## Transitions flow diagram

BORIS can generate DOT scripts and flow diagrams from the transitions matrices (See Observations > Create transition matrix for obtaining the transitions matrices).

### DOT script (*Graphviz language*)

#### Tools > Transitions flow diagram > Create transitions DOT script

Choose one ore more transitions matrix files and BORIS will create the relative DOT script file(s).

The DOT script files can then be used with [Graphviz](#) (Graph Visualization Software) or [WebGraphviz](#) (Graphviz in the Browser) to generate flow diagram of transitions.

See [DOT \(graph description language\)](#) for details.


### Flow diagram

If [Graphviz](#) (Graph Visualization Software) is installed on your system (and the **dot** program available in the path) BORIS can generate flow diagram (PNG format) from a transitions matrix file.


#### Tools > Transitions flow diagram > Create transitions flow diagram

Choose one ore more transition matrix files and BORIS will create the relative flow diagram.


### **Flow diagram of frequencies of transitions**

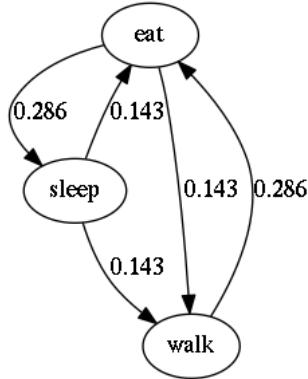

### **Flow diagram of frequencies of transitions after behavior**

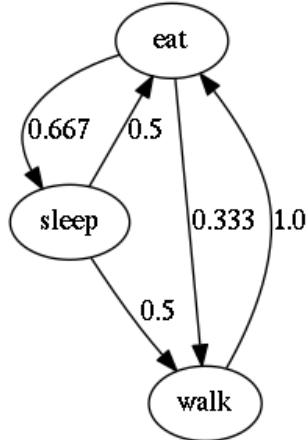

### **Flow diagram of number of transitions**

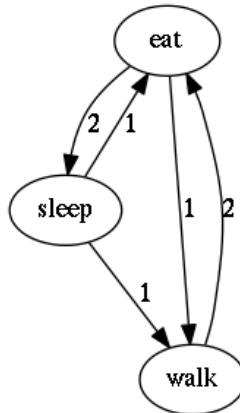

## **Re-encoding and resizing a video file**

BORIS can re-encode and resize your video files in order to reduce the size of the files and have a smooth coding (specially with two video files playing together). The re-encoding and resizing operations are done with the embedded ffmpeg program with high quality parameters (bitrate 2000k).

Select the files you want re-encode and resize and select the horizontal resolution in pixels (the default is 1024). The aspect ratio will be maintained.

You can continue to use BORIS during the re-encoding/resizing operation.

The re-encoded/resized video files are renamed by adding the re-encoded.avi extension to the original files.

### **Warning**

Please note that MTS video files should be re-encoded to be used in frame-by-frame mode. Otherwise the extracted frames are not reliable.

## **Rotating a video file**

BORIS can rotate your video files in order to code them using the right view. The rotating operation is done with the embedded ffmpeg program using the same quality parameters than the original video.

Select the files you want rotate and select the rotation between: **Rotate 90 clockwise**, **Rotate 90 counter clockwise** and **Rotate 180**.

The aspect ratio will be maintained.

You can continue to use BORIS during the rotation operation.

The rotated video files are renamed by adding the **rotated<ANGLE>** to the original file name.

## **Project server**

This function allow you to communicate with the [BORIS App](#) by sending project and receiving observations made with the [BORIS App](#). For details read the [BORIS App user guide](#)

### **Starting the project server**

- Open your project
- Select **File > Project server**

A window will open and show you the URL of the server that should be used on the [BORIS App](#). The project server will serve the project for 30 minutes.

### **Stopping the project server**

- Select **File > Stop serving project**

## **Coding map**

A coding map is a bitmap image with user-defined clickable areas that will help to code for behaviors or modifiers for a behavior.

2 types of coding maps are available:

- Behaviors coding map
- Modifiers coding map

## **The Behaviors coding map**

BORIS allows creating a **Behaviors coding map** using the **Map creator** tool (**Tools > Create a coding map > for behaviors**).


A **Behaviors coding map** can be created only if you have defined behaviors in your ethogram.

## **Creating a Behaviors coding map**

To create a new **Behaviors coding map** launch the **Behaviors coding map creator**


**Tools > Create a coding map > for behaviors).**

A new window will open

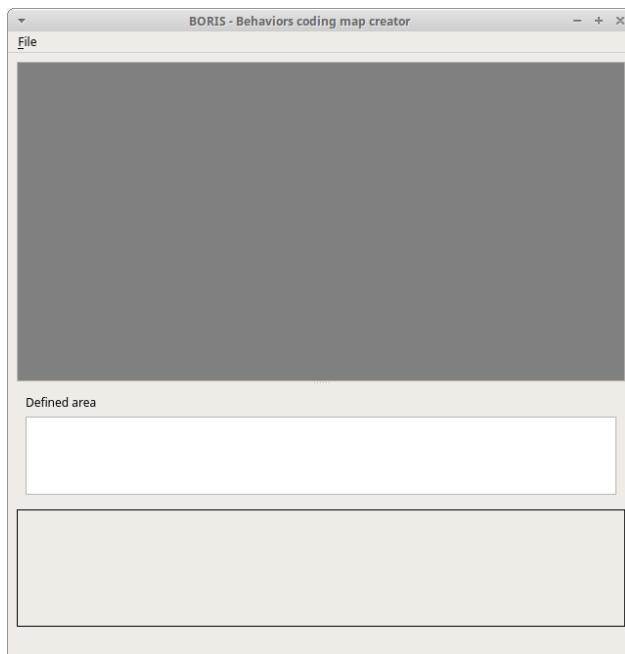

**File > New behaviors coding map**

Enter a name for the new **Behaviors coding map**

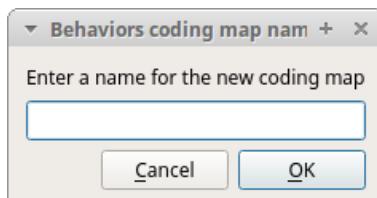


**Loading a bitmap for a behaviors coding map**

Click the **Load bitmap** button in the bottom of the window and select a bitmap image (PNG and JPEG formats are accepted).


If the size of your bitmap image is bigger than 640 x 640 pixels BORIS will resize it to 640 x 640 pixels keeping the aspect ratio and store the resized version in the coding map file.

The bitmap will be displayed

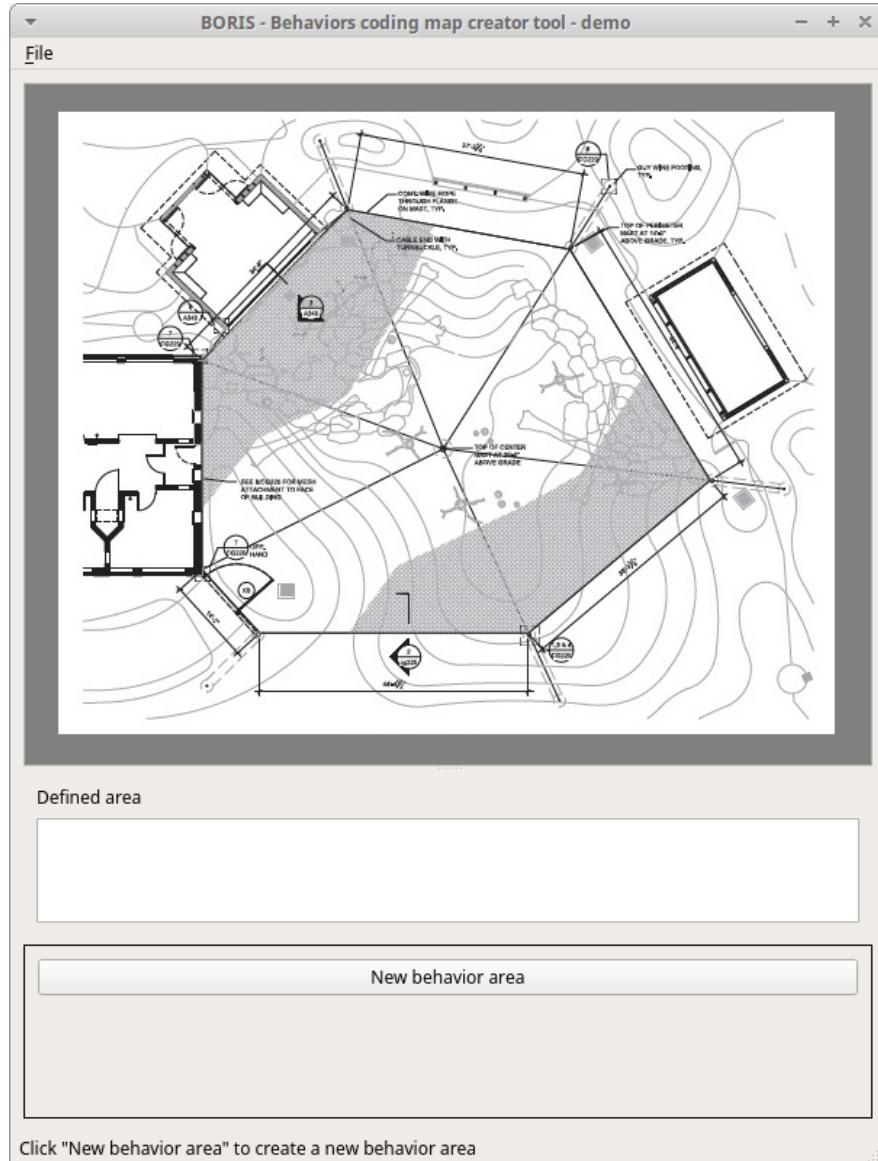

### **Adding areas corresponding the behaviors**

Click the **New behavior area** button in the bottom of the window and select a behavior by clicking on the **Select behavior** button.

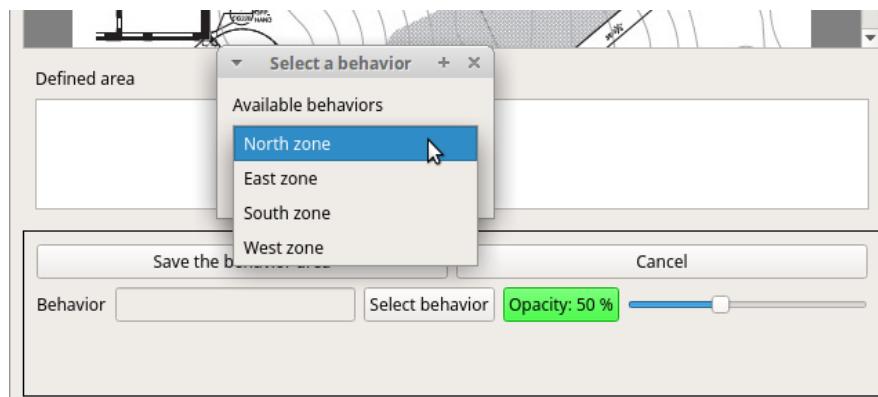

The available behaviors are taken from the ethogram of the current project.


Click on the bitmap to define the vertex on the area that will code the selected behavior. Close the area by clicking again on the first point.

The color of the new area can be changed using the **Opacity** button. The opacity can be changed (from 0 to 100%) using the slider.

Save the behavior area by clicking on the **Save the behavior area** button

The area will be added to the **Defined area** list

You can add more area and also add more than one area for a same behavior. Two or more areas can overlap. In this case all corresponding behaviors will be triggered.

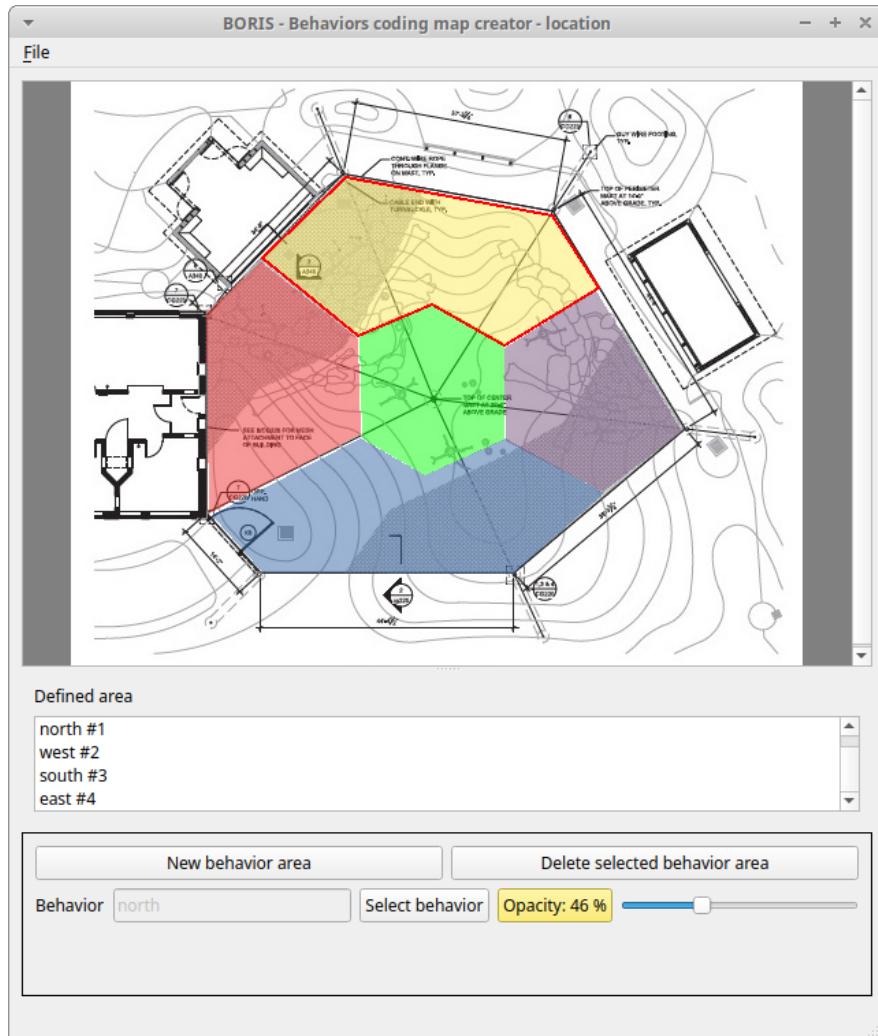

#### **Add the Behaviors coding map to the current project**

**File > Add coding map to project**

The coding map will be added to the current project

You can add a **Behaviors coding map** to the current project from a file containing the coding map:

(**File > Edit project > Behaviors coding map > Add a behaviors coding map** )

#### **Saving the Behaviors coding map**

Saving the **Behaviors coding map** will create a file containing the **Behaviors coding map** including the bitmap image.

**File > Save the current Behaviors coding map**

The file containing the **Behaviors coding map** can be then reloaded in the **Behaviors coding map creator** or added to a BORIS project (**File > Edit project > Behaviors coding map > Add a behaviors coding map** )


## The Modifiers coding map

BORIS allows creating a modifiers coding map using the **Modifiers Map creator** tool (**Tools > Create a coding map > for modifiers**.) Clickable areas may correspond to specific modifiers that can be meaningful for the behavioral coding. Facial expression is the case we thought to when developing this function.

### *Creating a modifiers coding map*

#### *Loading a bitmap for a modifiers coding map*

To create a new **Modifiers coding map**, launch the **Modifiers Map creator** tool (**Tools > Create a coding map > for modifiers**). The BORIS main window will be replaced by the **Modifiers Map creator** window. Click on **Modifiers Map creator > New Modifiers map** and enter a name for the new map in the edit box. You have to load a bitmap image (JPEG or PNG) using the **Load bitmap** button. The loaded image will be displayed.

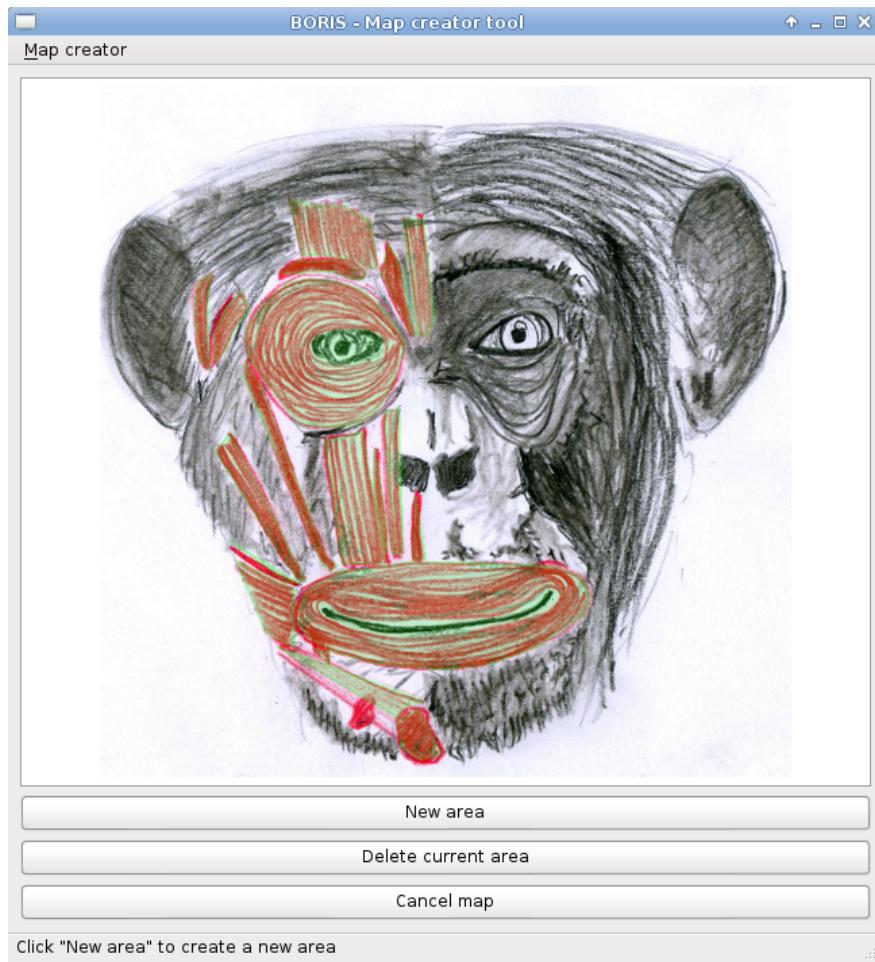

If the size of your bitmap image is bigger than 640 x 640 pixels BORIS will resize it to 640 x 640 pixels keeping the aspect ratio and store the resized version in the coding map file.

### **Adding areas corresponding to the modifiers**

To create clickable areas on a coding map, you have to click on the **New area** button and enter an **Area code** in the edit box. The new area can now be defined by clicking on the image. The drawing tool allows defining a irregular polygon (a plane shape with straight sides, which does not have all sides equal and all angles equal) by clicking to determine subsequent vertices. It can be convex or concave. Straight sides must not cross each other. Once selected an area can be deleted using the **Delete area** button. When an area is closed and its name has been defined in the **Area code** field, it can be saved by using the **Save area** button. The areas can partially overlap each other. See the **Using a Coding map** section for more details. Once all areas are added the entire map can be saved using the **Save map** option menu (**Map creator > Save map**). The map is now saved in its own file (.boris\_map) which is NOT part of the BORIS project. A map can be edited at anytime by opening the map file from the **Open map** menu option (**Map creator > Open map**).

### **Adding a modifiers coding map to your project**

Creating a Coding map is not automatically adding the map to your project. The Coding map have to be added to your project by selecting the corresponding **Behavior type (Point event with coding map, State event with coding map)**. BORIS will ask to select the file name containing the coding map (.boris\_map) and load the coding map in the project. The coding map name will appear in the **Coding map** column and will be saved in the BORIS project file.

NOTE: If you later modify your coding map you must reload the new version in your BORIS project.

## **Analysis**

### **Time budget**


The **Time budget** analysis can be done by behavior (including or not the modifiers) or by category of behaviors. Choose the option from the **Analysis** menu.

The **Time budget** analysis can be done on one or more observations. If you select more than one observation you must then choose for a global time budget analysis that will contain all selected observations or a time budget analysis for every single observation.



Choose **Yes** to group all observations in one time budget analysis

The **Analysis > Time budget** menu option shows the time budget for the events of the selected observations. Select the subjects and behaviors you want to include in the time budget analysis:

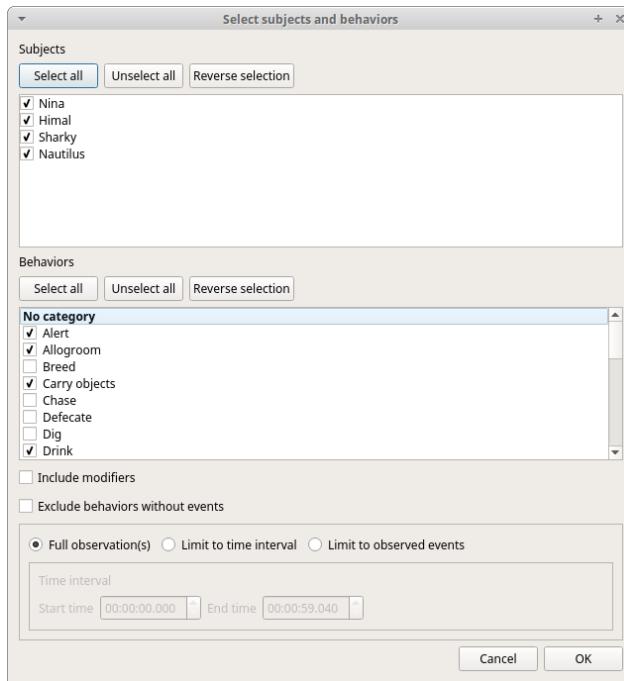


All behaviors can be selected or unselected by clicking on the Category (bold).

You can choose to include or not the behavior modifiers in the Time budget analysis and to exclude behaviors without coded events.

The Time budget analysis can be restricted to a part of the observation:

- Full observation(s): the analysis will be made on the full observation length.
- Limit to time interval: use the **Start time** and **End time** boxes to select starting time and ending time.
- Limit to observed events: the analysis will be made from the first observed event to the last observed event.

The last dialog window will allow you to subtract the duration of one or more behaviors from the total duration of the observation. This can help you if you have defined a "out-of-sight" code in your ethogram for example:

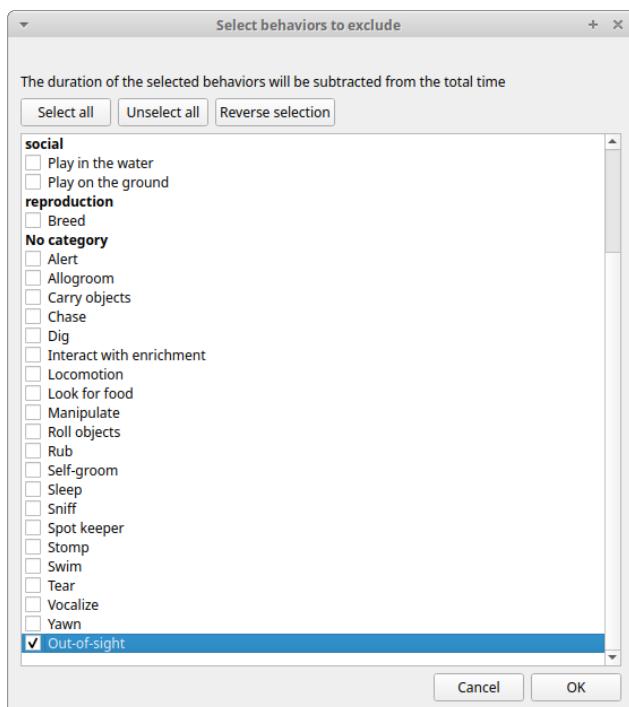

The results contain for each subject and behavior the **total number of occurrences**, the **total duration** (for the behaviors defined as state events), The **duration mean** (for the behaviors defined as state events), the **standard deviation of duration**, the **inter-events intervals duration mean**, th **standard deviation of the inter-events intervals duration** and the **percent of total duration of observation(s)**. All duration times are expressed in seconds (s).

| Time budget           |         |                    |           |              |                    |                   |                  |                                |                               |                   |
|-----------------------|---------|--------------------|-----------|--------------|--------------------|-------------------|------------------|--------------------------------|-------------------------------|-------------------|
| Selected observations |         |                    |           |              |                    |                   |                  |                                |                               |                   |
|                       | Subject | Behavior           | Modifiers | Total number | Total duration (s) | Duration mean (s) | Duration std dev | inter-event intervals mean (s) | inter-event intervals std dev | % of total length |
| 1                     | Nina    | Alert              | -         | 1            | 1.882              | 1.882             | NA               | NA                             | NA                            | 0.4               |
| 2                     | Nina    | Locomotion         | -         | 2            | 2.913              | 1.457             | 0.824            | NA                             | NA                            | 0.7               |
| 3                     | Nina    | Manipulate         | -         | 2            | 65.048             | 32.524            | 15.414           | NA                             | NA                            | 14.6              |
| 4                     | Nina    | Rest               | -         | 1            | 25.549             | 25.549            | NA               | NA                             | NA                            | 5.7               |
| 5                     | Himal   | Alert              | -         | 4            | 18.131             | 4.533             | 2.542            | 10.95                          | NA                            | 4.1               |
| 6                     | Himal   | Allogroom          | -         | 1            | 3.603              | 3.603             | NA               | NA                             | NA                            | 0.8               |
| 7                     | Himal   | Drink              | -         | 1            | 1.133              | 1.133             | NA               | NA                             | NA                            | 0.3               |
| 8                     | Himal   | Locomotion         | -         | 10           | 36.633             | 3.663             | 5.777            | 6.124                          | 5.88                          | 8.2               |
| 9                     | Himal   | Manipulate         | -         | 3            | 66.702             | 22.234            | 17.18            | 4.723                          | NA                            | 15.0              |
| 10                    | Himal   | Play on the ground | -         | 1            | 16.46              | 16.46             | NA               | NA                             | NA                            | 3.7               |
| 11                    | Himal   | Rub                | -         | 1            | 11.031             | 11.031            | NA               | NA                             | NA                            | 2.5               |
| 12                    | Himal   | Self-groom         | -         | 1            | 3.571              | 3.571             | NA               | NA                             | NA                            | 0.8               |

Save results Close

The time budget results can be saved in various formats for further analysis: \* Plain text files: TSV, CSV or HTML \* Spreadsheets: Open Document (ODS), Microsoft Excel (XLSX) and Legacy Microsoft (XLS)

## Note

If a STATE behavior has an odd number of coded events, BORIS will report "UNPAIRED" instead of results"

## Synthetic time budget

The synthetic time budget is similar to time budget but with fewer parameters and a different organization of results. Results of all selected observations are organized in columns on a single page. Two parameters are provided for now: **number of occurrences** and **total duration** (for the behaviors defined as state events)

|    | A                | B       | C              | D                     | E              | F                    | G              | H                    | I              | J                    |
|----|------------------|---------|----------------|-----------------------|----------------|----------------------|----------------|----------------------|----------------|----------------------|
| 1  |                  |         | Nina           | Nina                  | Nina           | Nina                 | Nina           | Nina                 | Nina           | Nina                 |
| 2  |                  |         | Alert          | Alert                 | Allogroom      | Allogroom            | Breed          | Breed                | Carry objects  | Carry objects        |
| 3  | Total length (s) |         | Total duration | Number of occurrences | Total duration | Number of occurrence | Total duration | Number of occurrence | Total duration | Number of occurrence |
| 40 | 0037             | 32.160  |                |                       |                |                      |                |                      |                |                      |
| 41 | 0038             | 86.880  | 1.675          | 1                     |                |                      |                |                      |                |                      |
| 42 | 0039             | 96.960  | 13.471         | 2                     |                |                      |                |                      |                |                      |
| 43 | 0040             | 335.520 | 86.615         | 12                    | 15.84          | 2                    |                |                      |                |                      |
| 44 | 0041             | 130.560 | 9.069          | 2                     | 5.85           | 1                    |                |                      |                |                      |
| 45 | 0042             | 36.960  | 24.298         | 3                     |                |                      |                |                      |                |                      |
| 46 | 0043             | 107.040 | 42.928         | 5                     |                |                      |                |                      |                |                      |
| 47 | 0044             | 109.920 | 11.769         | 4                     | 12.471         | 2                    |                |                      |                |                      |
| 48 | 0045             | 102.240 | 72.857         | 6                     |                |                      |                |                      |                |                      |
| 49 | 0046             | 84.960  | 44.459         | 6                     |                |                      |                |                      |                |                      |
| 50 | 0047             | 77.280  | 68.929         | 2                     |                |                      |                |                      |                |                      |
| 51 | 0048             | 23.520  | 3.1            | 1                     | 10.797         | 1                    |                |                      |                |                      |
| 52 | 0049             | 34.080  | 9.897          | 2                     | 8.431          | 1                    |                |                      |                |                      |
| 53 | 0050             | 23.040  |                |                       |                |                      |                |                      |                |                      |
| 54 | 0051             | 44.160  | 1.079          | 1                     |                |                      |                |                      |                |                      |
| 55 | 0052             | 18.240  |                |                       |                |                      |                |                      |                |                      |
| 56 | 0053             | 39.840  |                |                       |                |                      |                |                      |                |                      |
| 57 | 0054             | 35.040  |                |                       |                |                      |                |                      |                |                      |
| 58 | 0055             | 256.800 | 61.321         | 15                    | 19.73          | 1                    |                |                      |                |                      |
| 59 | 0056             | 26.400  | 15.241         | 5                     |                |                      |                |                      |                |                      |
| 60 | 0057             | 45.120  | 6.887          | 3                     | 10.984         | 2                    |                |                      |                |                      |
| 61 | 0058             | 45.120  | 2.647          | 1                     |                |                      |                |                      |                |                      |
| 62 | 0059             | 41.760  |                |                       |                |                      |                |                      |                |                      |
| 63 | 0060             | 292.800 | 1.543          | 2                     |                |                      |                |                      |                |                      |
| 64 | 0061             | 25.920  |                |                       |                |                      |                |                      |                |                      |

All duration times are expressed in seconds (s).


The time budget results can be saved in various formats for further analysis:  
\* Plain text files: TSV, CSV or HTML  
\* Spreadsheets: Open Document (ODS), Microsoft Excel (XLSX) and Legacy Microsoft (XLS)

## Plot events

The coded events from many observations can be plotted along a time axis.

### Analysis > Plot events

The subjects and behaviors you want to include in the plot can be selected in the following window:

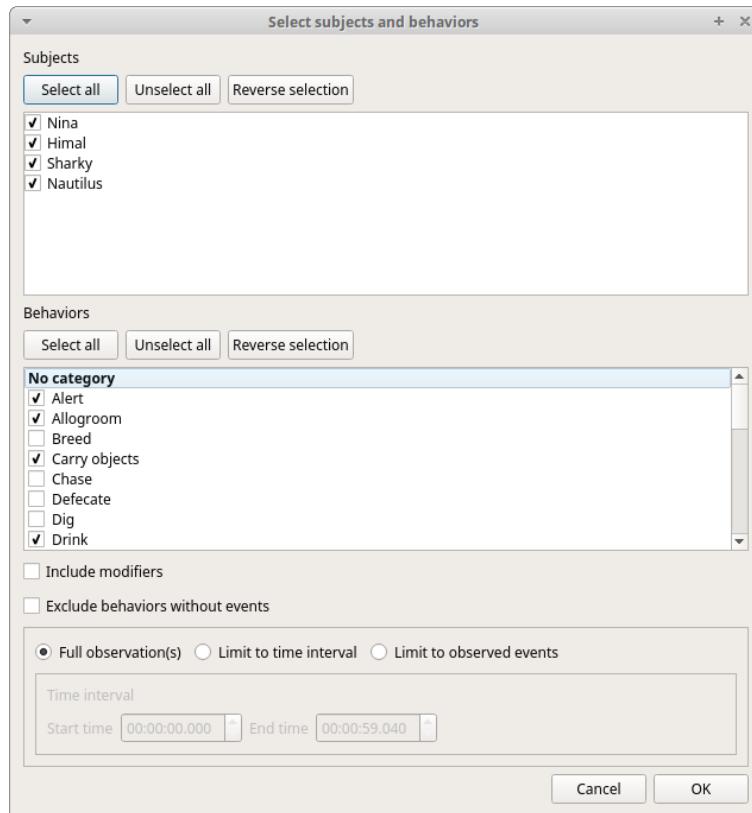

You can choose to include or not the behavior modifiers (if any) and to exclude behaviors without coded events.

The time interval can be selected (See time budget)


The plot can be exported in various formats like bitmap (PNG, JPG, TIFF) or vectorial graphic (SVG, PDF, EPS, PS). The SVG format can be further edited with the [Inkscape vector graphics editor](#).

### Note

If a STATE behavior has an odd number of coded events, you will see this error message: "The STATE behavior XXX is not paired"

This function creates one plot by subject on one figure.

The color of behaviors can be customized. See [plot colors](#)

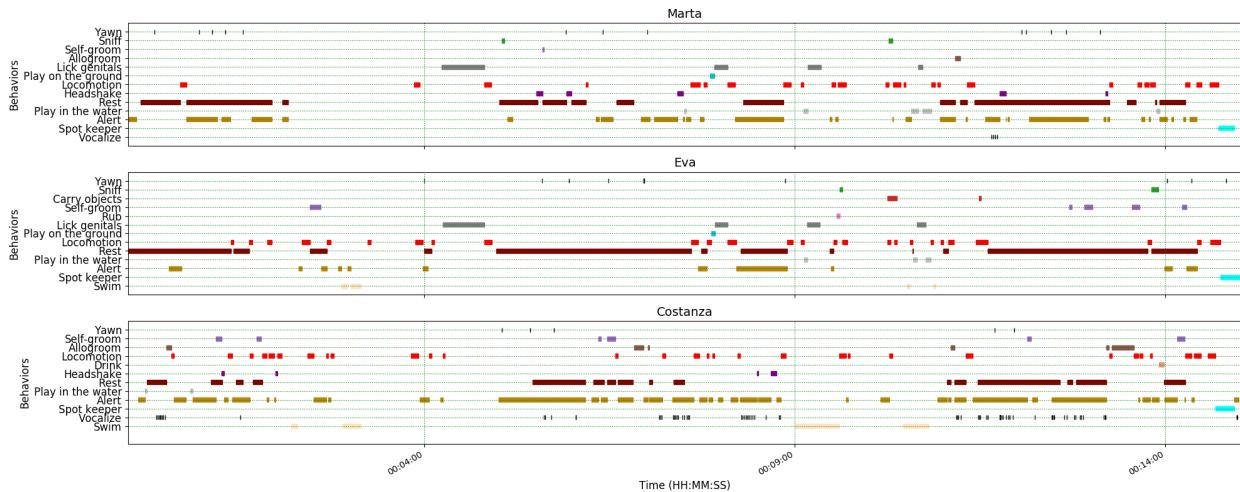


## Inter-rater reliability

The Cohen's kappa coefficient can be calculated (**Analysis > Inter-rater reliability > Cohen's kappa**).

[Cohen's kappa on Wikipedia](#)

After selecting 2 observations and a time window (in seconds) for the analysis (the default value is 10 seconds) the Cohen's kappa will be displayed in the results window.

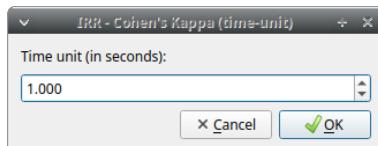

### Implementation of IRR Cohen's Kappa in BORIS

If a time window of n seconds is set the 2 selected observations will be checked every n seconds for agreement/disagreement from the first event to the last event of the 2 observations . In case of a point event the presence of a corresponding event in the other observation will be verified using a time window of n seconds centered on the point event.

A IRR Cohen's Kappa analysis is available in the GSEQ program (<http://www2.gsu.edu/~psyrab/gseq>). For this the coded events can be exported as aggregated events in SDIS format. See [export aggregated events](#)

## Advanced event filtering

### Warning

This function is still experimental in v. 7.8.5

This function allows to filter events from one or more observations by selecting subjects, behaviors and logical operators.

To use the filter, select a subject, select a behavior and click on the button with the green arrow on the side of the behaviors list. The tuple subject/behavior will be added in the **filter text edit**. A complex filter query can be constructed by adding parenthesis and logical operator (**AND &** / **OR |**) for combining various subjects and behaviors.

Example of a summarized output showing the occurrences of Himal resting while Nautilus in alert:

Advanced event filtering

Filter: "Himal|Rest" & "Nautilus|Alert"

Subjects: Himal, Nautilus, Nina, Sharky

Behaviors: Alert, Allogroom, Breed, Carry objects, Chase, Defecate, Dig, Drink, Eat, Interact with enrichment, Locomotion

Logical operators: AND

Results (197 observations)

|   | Observation id | Number of occurrences | Total duration | Mean   | Std Dev |
|---|----------------|-----------------------|----------------|--------|---------|
| 1 | 0015           | 1                     | 3.351          | 3.351  | NA      |
| 2 | 0019           | 2                     | 8.581          | 4.29   | 2.461   |
| 3 | 0028           | 1                     | 3.905          | 3.905  | NA      |
| 4 | 0030           | 1                     | 8.231          | 8.231  | NA      |
| 5 | 0032           | 1                     | 14.245         | 14.245 | NA      |
| 6 | 0043           | 1                     | 9.219          | 9.219  | NA      |
| 7 | 0048           | 1                     | 2.016          | 2.016  | NA      |

Save results Close

Example of a detailed output showing the overlapping intervals while Himal rests and Nautilus is in alert:

Advanced event filtering

Filter: "Himal|Rest" & "Nautilus|Alert"

Subjects: Himal, Nautilus, Nina, Sharky

Behaviors: Alert, Allogroom, Breed, Carry objects, Chase, Defecate, Dig, Drink, Eat, Interact with enrichment, Locomotion

Logical operators: AND

Results (146 events)

|   | Observation id | Comment | Start time | Stop time | Duration |
|---|----------------|---------|------------|-----------|----------|
| 1 | 0019           |         | 3.559      | 6.088     | 2.529    |
| 2 | 0019           |         | 14.907     | 17.457    | 2.550    |
| 3 | 0028           |         | 0.0        | 3.905     | 3.905    |
| 4 | 0048           |         | 18.333     | 18.627    | 0.294    |
| 5 | 0054           |         | 28.77      | 30.019    | 1.249    |
| 6 | 0138           |         | 174.796    | 178.699   | 3.903    |
| 7 | 0138           |         | 180.77     | 181.407   | 0.727    |

Save results Close

The same subject can be used many times in the query with OR or AND (in case of non exclusive behaviors):

**Advanced event filtering**

**Filter**  
"Himal|Drink" | "Himal|Eat"

Summary       Details

**Subjects**  
Himal  
Nautilus  
Nina  
Sharky

**Behaviors**  
Carry objects  
Chase  
Defecate  
Dig  
Drink  
**Eat**        
Interact with enrichment  
Locomotion  
Look for food  
Manipulate  
Play in the water  
Sleep on the ground

**Logical operators**  
AND  
**OR**      

**Results (188 events)**

| Observation id | Comment | Start time | Stop time | Duration |
|----------------|---------|------------|-----------|----------|
| 1 0004         |         | 49.275     | 50.408    | 1.133    |
| 2 0009         |         | 5.579      | 10.48     | 4.901    |
| 3 0009         |         | 18.683     | 40.19     | 21.507   |
| 4 0013         |         | 24.324     | 30.869    | 6.545    |
| 5 0055         |         | 3.65       | 5.326     | 1.676    |
| 6 0055         |         | 21.362     | 29.117    | 7.755    |
| 7 0056         |         | 22.885     | 26.37     | 3.485    |

**Save results** **Close**

An unlimited number of conditions can be used:

**Advanced event filtering**

**Filter**  
"Himal|Rest" & "Nautilus|Rest" & "Nina|Rest" & "Sharky|Rest"

Summary       Details

**Subjects**  
Himal  
Nautilus  
Nina  
Sharky

**Behaviors**  
Rest        
Roll objects  
Rub  
Self-groom  
Sleep  
Sniff  
Spot keeper  
Stomp  
Swim  
Tear  
Urinate

**Logical operators**  
AND  
OR      

**Results (2 observations)**

| Observation id | Number of occurrences | Total duration | Mean  | Std Dev |
|----------------|-----------------------|----------------|-------|---------|
| 1 0999         | 1                     | 0.403          | 0.403 | NA      |
| 2 1000         | 2                     | 2.769          | 1.384 | 1.592   |

**Save results** **Close**

Parenthesis can be used to group logical conditions into block(s):

Advanced event filtering

**Filter**

Summary  Details

| Subjects                                   | Behaviors   | Logical operators |
|--|---|-------------------|
| Himal<br>Nautilus<br>Nina<br><b>Sharky</b> | Carry objects<br>Chase<br>Defecate<br>Dig<br>Drink<br><b>Eat</b><br>Interact with enrichment<br>Locomotion<br>Look for food<br>Manipulate<br>Play in the water<br>Sleep on the ground | AND<br><b>OR</b>  |

**Results (30 events)**


| Observation id | Comment | Start time | Stop time | Duration |
|----------------|---------|------------|-----------|----------|
| 1 0055         |         | 8.689      | 12.435    | 3.746    |
| 2 0055         |         | 19.778     | 21.361    | 1.583    |
| 3 0055         |         | 29.118     | 30.692    | 1.574    |
| 4 0055         |         | 66.079     | 69.591    | 3.512    |
| 5 0055         |         | 199.744    | 200.81    | 1.066    |
| 6 0239         |         | 149.63     | 152.719   | 3.089    |
| 7 0230         |         | 162.807    | 164.850   | 1.052    |

The results can be saved in a Tab Separated Values (TSV) file using the **Save results** button. Other formats will be added in future.

## Preferences

You can customize BORIS using the Preferences window (**File > Preferences**)

## General preferences

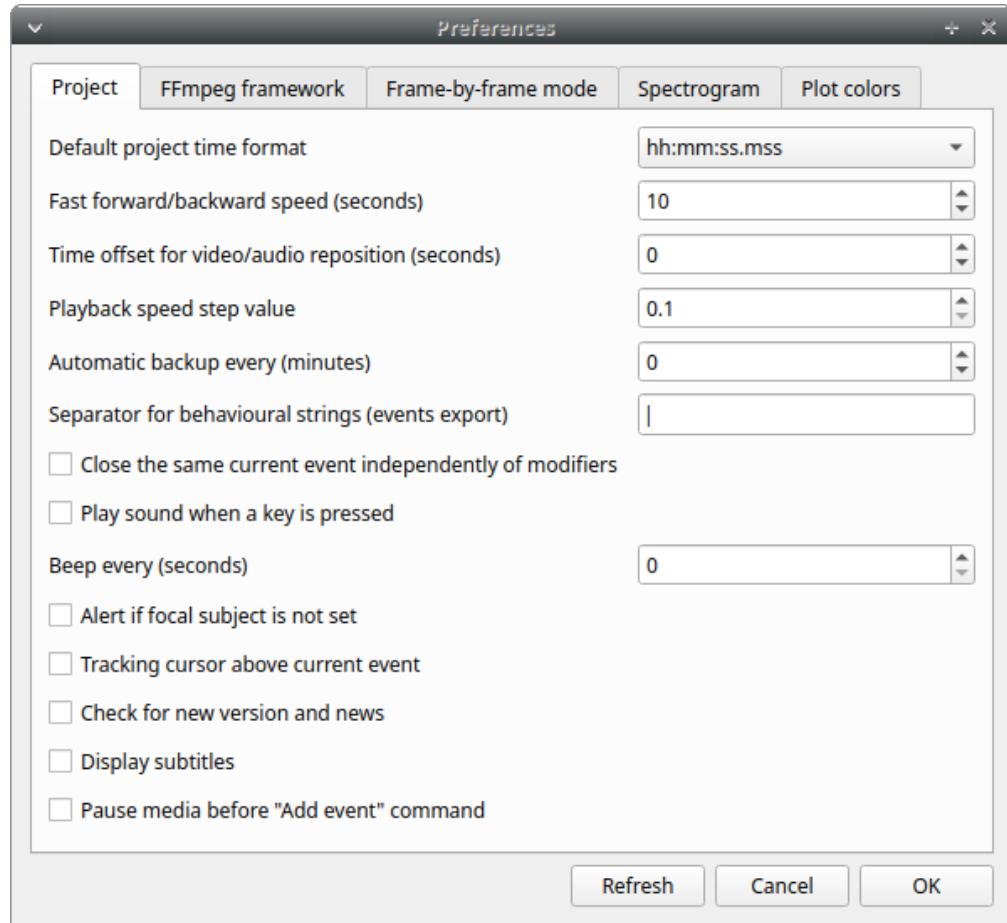

### Refresh button

Option to reinitialize the configuration to default. BORIS will be closed.

### Default project time format

This option allows the user to choose the format for displaying time in the project. Please note that time is internally always saved in seconds with a precision of 3 decimal digits

### Fast forward/backward speed (seconds)

This option allows the user to customize the amount of time for "jumping" forward or backward in media.

### Time offset for media reposition (seconds)

This value indicates the time offset for repositioning the media after double-click on a row event of the *Events* table. 'for example -4 seconds indicates that after a double-click the media will be repositioned 4 seconds before the recorded event.'

### Playback speed step value

This value indicate how much the speed will be increased or decreased after pressing the *change playback speed* buttons.

### Auto-save project every (minutes)

If set BORIS will save your project automatically every n minutes. 0 indicate no automatic backup. The project will be saved if the project is already saved and an observation is open.

### Separator for behavioural strings

Character (or string) used to separate behaviors when exporting events as behavioural strings. See also Behatrix

### Close the same current event independently of modifiers

Option used to STOP the current behavior without regarding the modifiers

#### **Play sound when a key is pressed**

Activate a sound signal after every keypress event

#### **Alert if focal subject is not set**

If this option is activated BORIS will show an alert box if no focal subject is selected

#### **Tracking cursor above current event**

Check this box to position the tracking cursor above the current event in events list table.

#### **Check for new version**

Check for new version on BORIS web site every 15 days (internet access required)


#### **Display subtitles**

Option to display or hide the visualization of subtitles. In case of separate file, the file containing subtitles must have the same base name than the video files with a .srt extension.

#### **Pause media before "Add event" command**

Option to pause the media before manually adding an event.

## **FFmpeg framework**

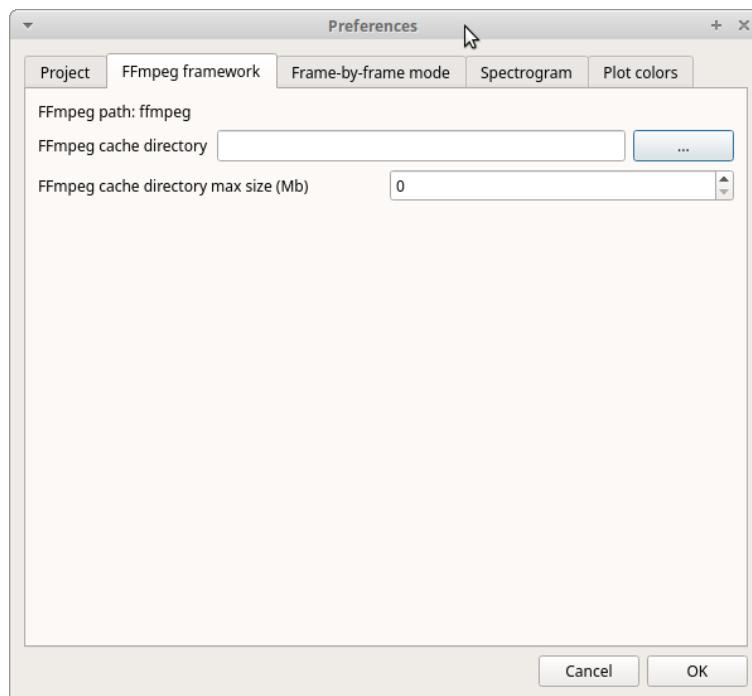

The path for the ffmpeg executable program is displayed. The FFmpeg executable is now included with BORIS in Windows and Mac OS versions.


#### **FFmpeg cache directory**

This indicates the directory that will be used as image cache for frame-by-frame mode and spectrogram visualization. If you do not specify a path, BORIS will use the default temporary directory of your system.

#### **FFmpeg cache directory max size**

Indicate a size limit (in Mb) for the image caching. 0 indicates no limit

## Frame-by-frame mode

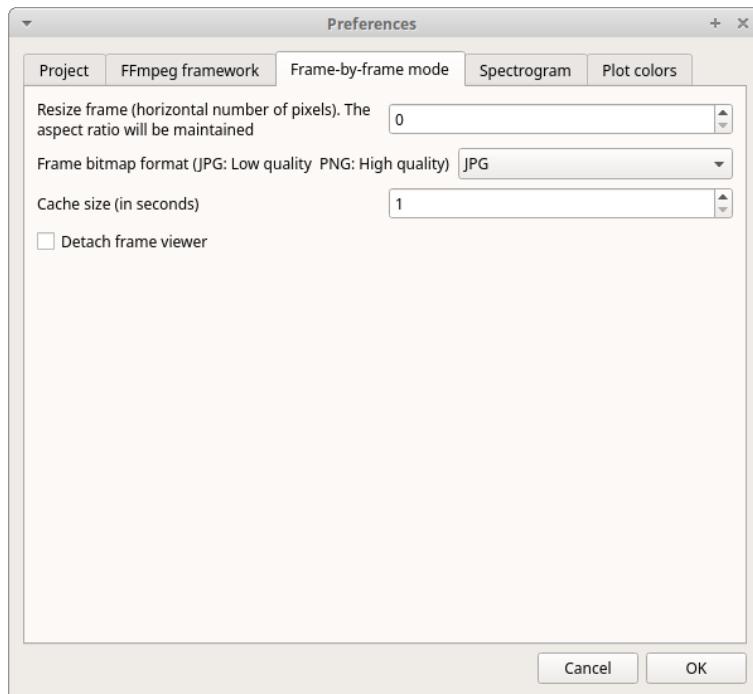

### **Resize frame**

In case of high resolution video (for example 4K 1920x1080) the extracted frames can be resized before visualization to improve the speed. Select the width in pixel for the extracted frames. The aspect ratio will be maintained. 0 will indicate no frame resizing.

### **Frame bitmap format**

A bitmap format can be selected between JPG (JPEG low quality image) and PNG (Portable Network Graphic - high quality image). The extracted frames will be saved in the directory defined in previous tab (FFmpeg cache directory).


### **Cache size**

In frame-by-frame mode the frames are extracted from the video with the **ffmpeg** program every second by default.

This value can be increased to extract frames for a longer sequence and avoiding waiting for frames extraction every second.

Note for high resolution video this operation can be long.

# Spectrogram

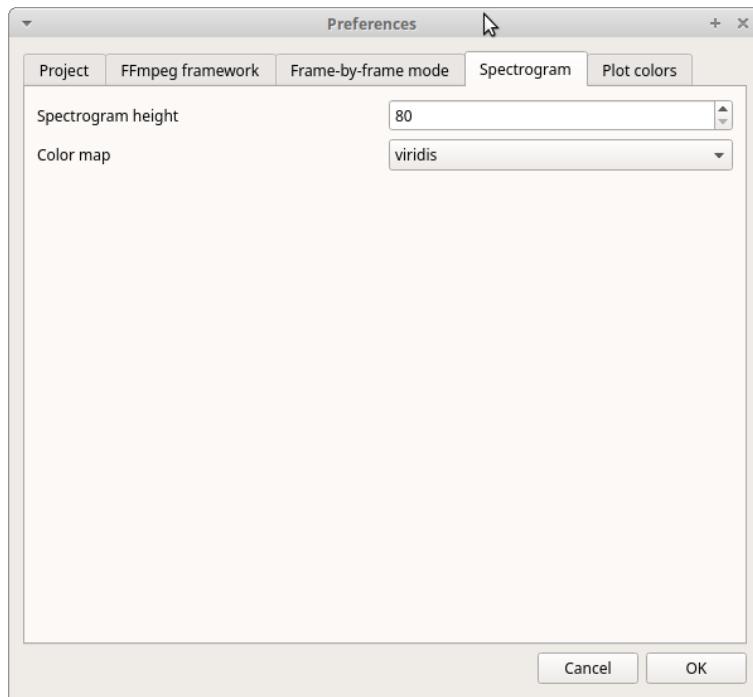

## **Spectrogram height**

Select the height of generated spectrogram (in pixels). You will need to restart the current observation to apply changes.


## **Color map**

Select the color map for displaying the generated spectrogram. See [Matplotlib colormaps](#) for details.

## **Plot colors**

The color of behaviors in the plot events functions can be customized. The first color will be associated to the first behavior in your ethogram, the second color to the second behavior and so on. Various color formats can be used to specify a color: **named color** or **hex RGB** (like #0F0F0F). See [https://matplotlib.org/api/colors\\_api.html](https://matplotlib.org/api/colors_api.html) and [https://matplotlib.org/examples/color/named\\_colors.html](https://matplotlib.org/examples/color/named_colors.html) for details

The **reset colors to default** button will reload the default colors.

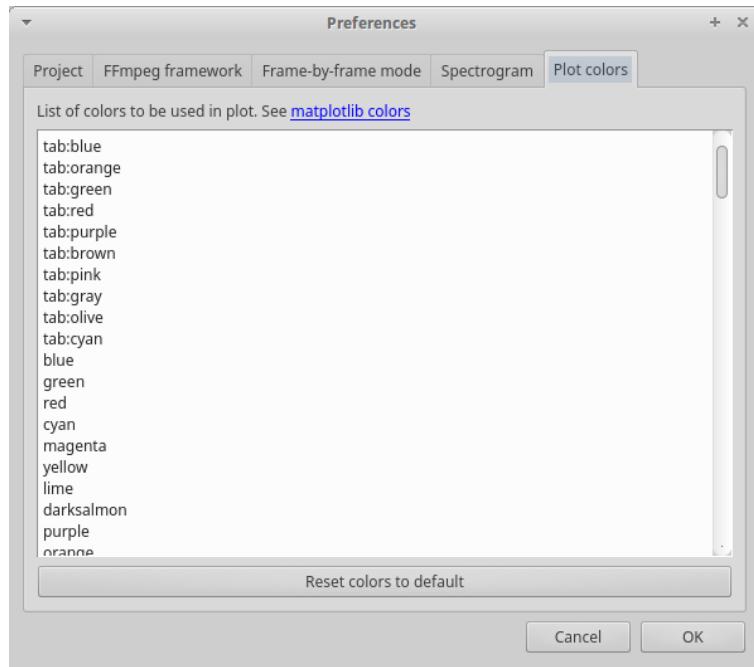

## Various

### Removing path of media files

Using BORIS you can choose to store the full path of the media files into the file project (for example: /home/user/Video/video\_n1.mp4 or c:\Users\user\Documents\video1.avi).

If you want to move your project on a different computer or if you want to move your media files you may want to do not store the full path. For this you can choose to add media files without the full path (See **Add media files** section). You can also remove the full path of your media files from the all observations of the current project (**File > Remove path from media files**). Please note that this operation is irreversible. After remotion the full path of your media will be lost and will not be recoverable.

**If you choose to do not store the full path of media files the media files must located in the same directory than your BORIS project file.**

### Citing BORIS


If you have used BORIS for publications, please cite:

Olivier Friard and Marco Gamba. (2016) BORIS: a free, versatile open-source event-logging software for video/audio coding and live observations. *Methods in Ecology and Evolution*, 7(11), 1324–1330  
DOI: 10.1111/2041-210X.12584

### Docking


All elements, including all the media players can be undocked from the main window and positioned where you prefer (e.g. they can be on the same desktop over one or two screens).

The position of the various widgets is saved in the [configuration file](#) at the end of the work session.

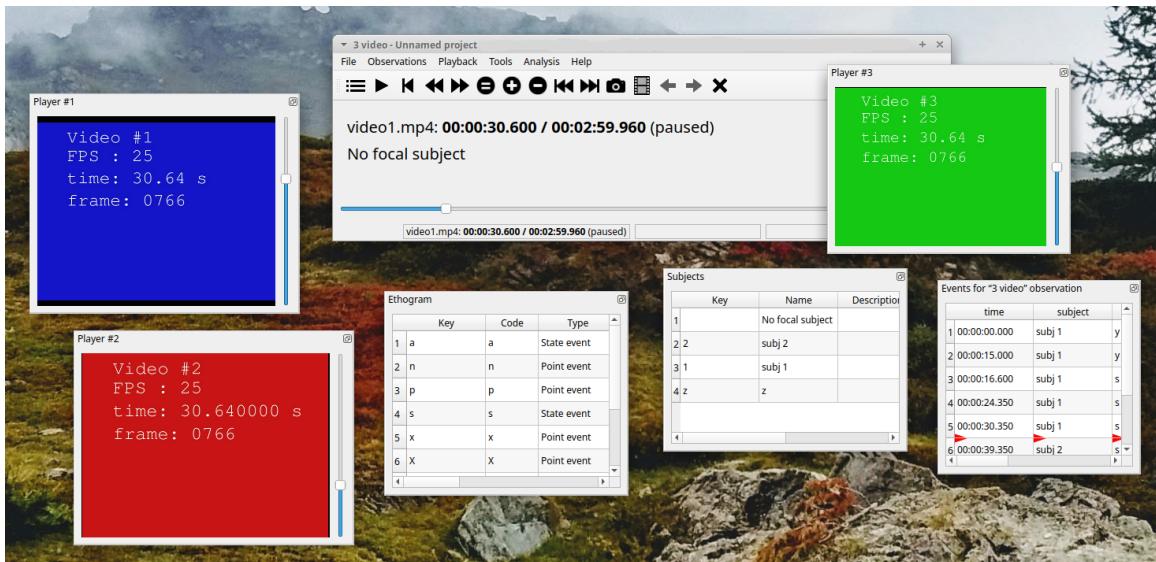


Click the icon present at the top-right corner of the widget (for MacOS the icon is located at the left-top corner) will undock the widgets that can be repositioned on another docking area or moved out of the main window. A double-click on the top bar of the widget will reposition it on the main window.


For Linux:

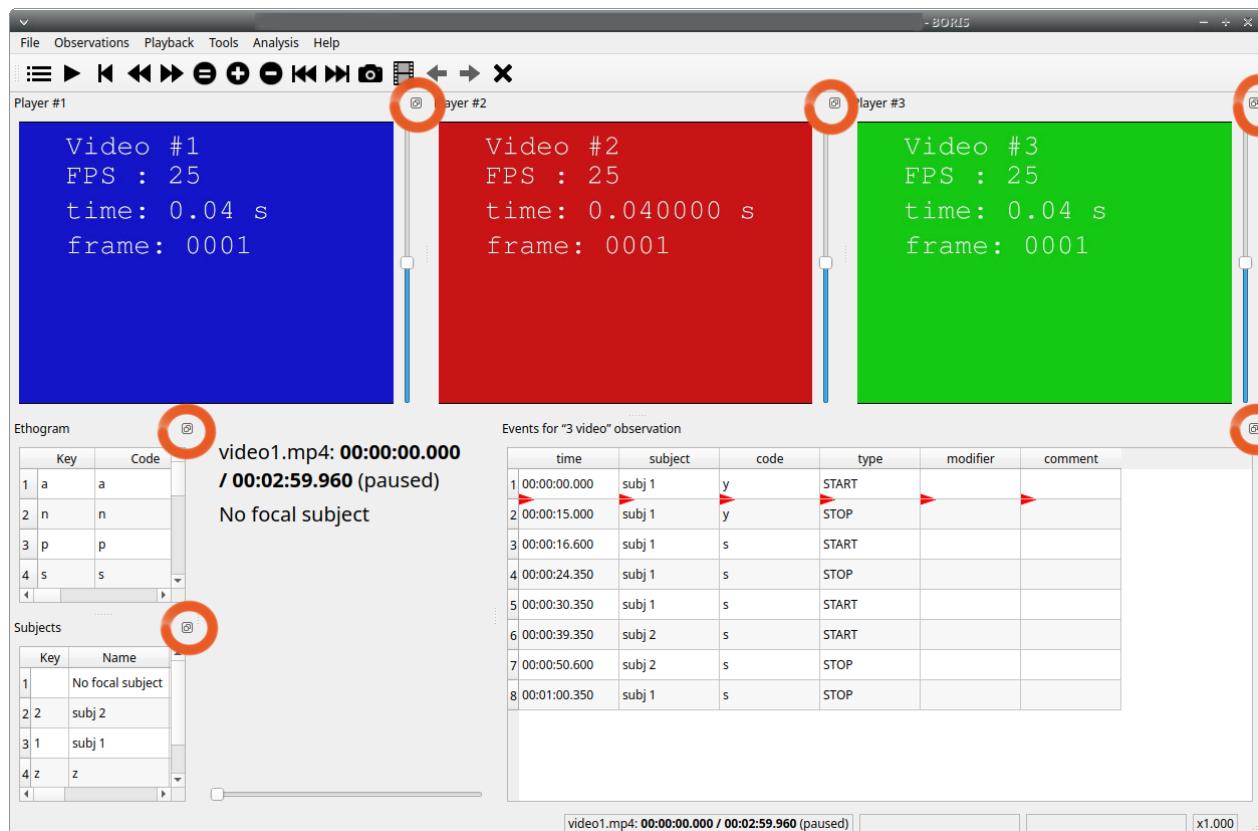

For Microsoft-Windows:

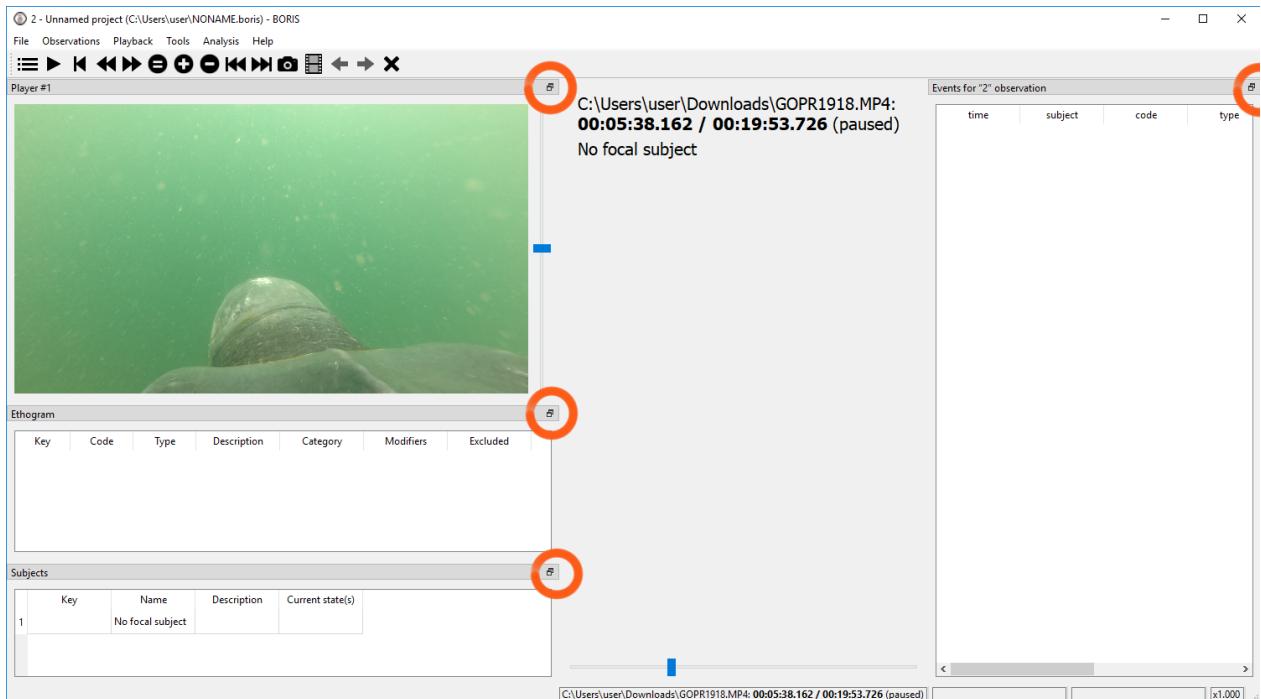

For MacOS:

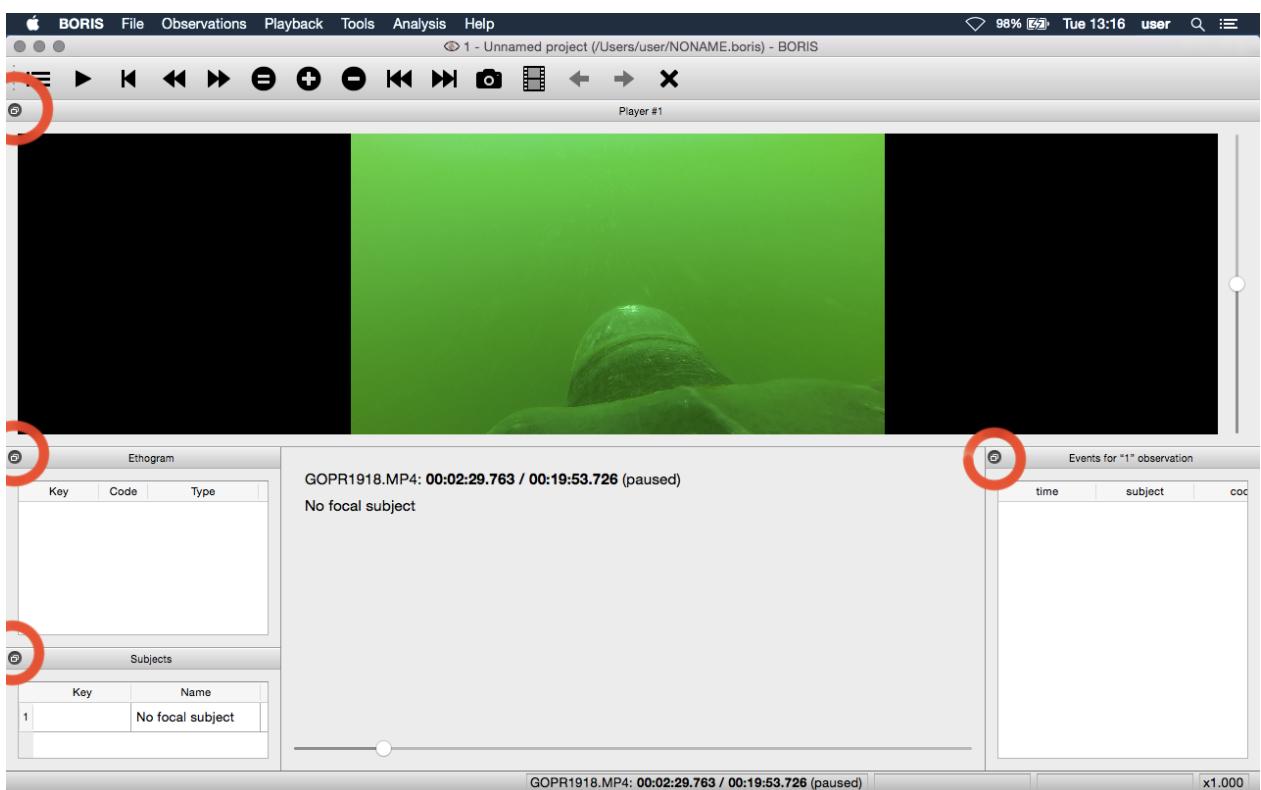

## Configuration files

BORIS saves the configuration (user preferences, windows position, widgets position) in a configuration file. This file is named **.boris** and it is saved on the home directory of the current user:

```
for Linux:  
/home/USERNAME/.boris
```

```
for Microsoft-Windows:  
C:\Users\USERNAME\.boris  
  
for MacOS:  
/Users/USERNAME/.boris
```

If you have some trouble using BORIS try to close the program, delete this file and relaunch BORIS.

The **recent projects list** is saved on the **.boris\_recent\_projects** file in the home directory of the current user.

## Bug reports and features request

Please send bug reports and features request using the BORIS GitHub repository (<https://github.com/olivierfriard/BORIS>) or by e-mail in english, french or italian (see web site <http://www.boris.unito.it/pages/about.html>).

In case of bug report please verify that you are using the **lastest version of BORIS** and indicate your operating system, its version and the CPU architecture (32/64 bits). You may also include the BORIS project that gave you an error. Any information you will provide will not be disclosed to any third party.

A debugging version (for Windows 64-bit) is available if you want to collaborate to BORIS.

## Acknowledgement

The authors would like to acknowledge Valentina Matteucci, Laura Ozella, Giovanni Di Panfilo and all others users for their precious help.

## Valid keys for triggering behavior

BORIS makes difference between lower case and upper case characters

- keys from a to z
- keys from A to Z
- keys from 0 to 9
- function keys from F1 to F12
- à é è ù ì ç
- ! " £ \$ % & / ( ) = ? ^ [ ] { } @ | § ° #